

Neuronal Input Strategies for Functional MRI

Peter A. Bandettini, Ph.D

Unit on Functional Imaging Methods
&
3T Neuroimaging Core Facility

Laboratory of Brain and Cognition
National Institute of Mental Health

Contrast in Functional MRI

- **Blood Volume**

- Contrast agent injection and time series collection of T2* or T2 - weighted images

- **BOLD**

- Time series collection of T2* or T2 - weighted images

- **Perfusion**




- T1 weighting
- Arterial spin labeling

- **CMRO₂**

- BOLD and Perfusion w/
Normalization to Global Perfusion Change

Resting

Active

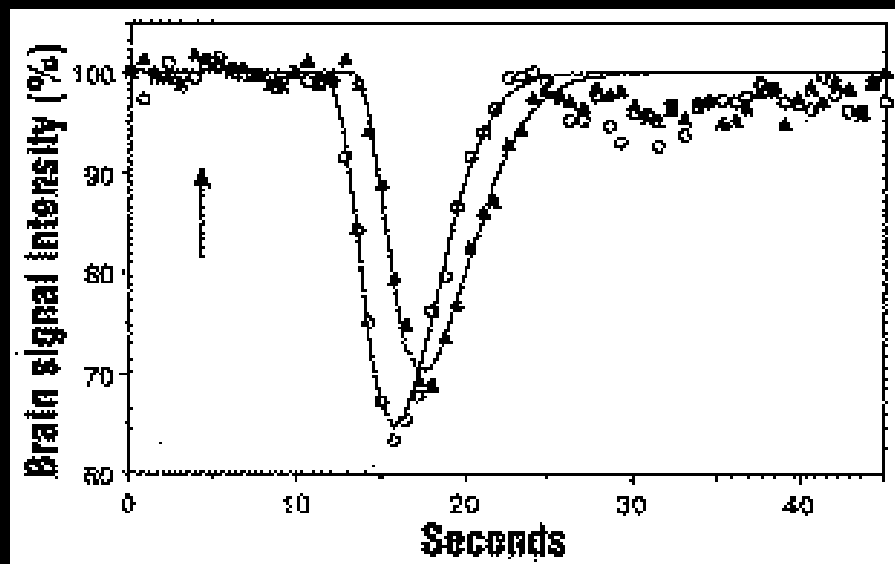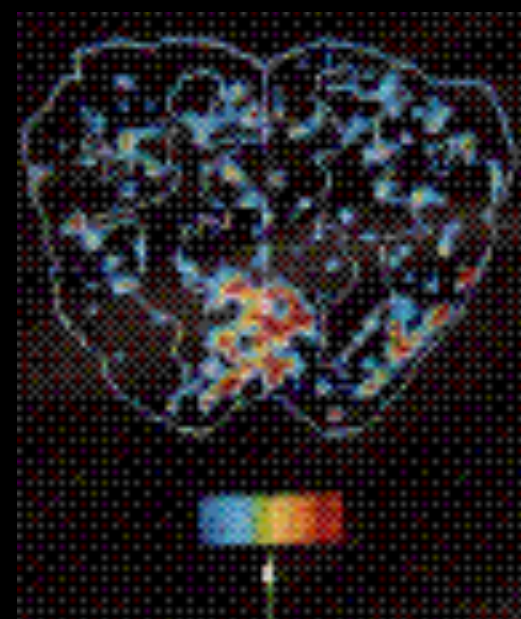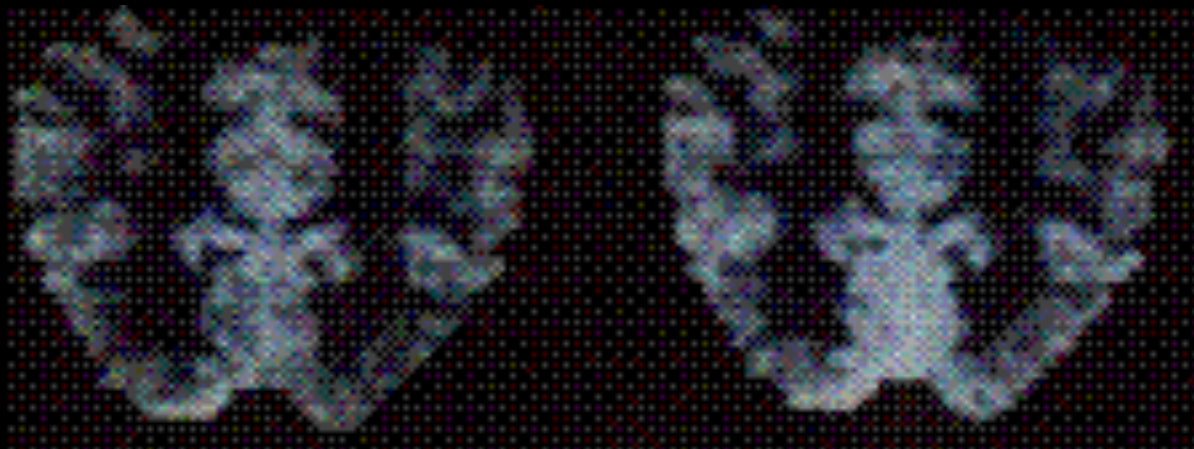


Photic Stimulation


MRI Image showing
activation of the
Visual Cortex

From Belliveau, et al.
Science Nov 1991


MSC - perfusion

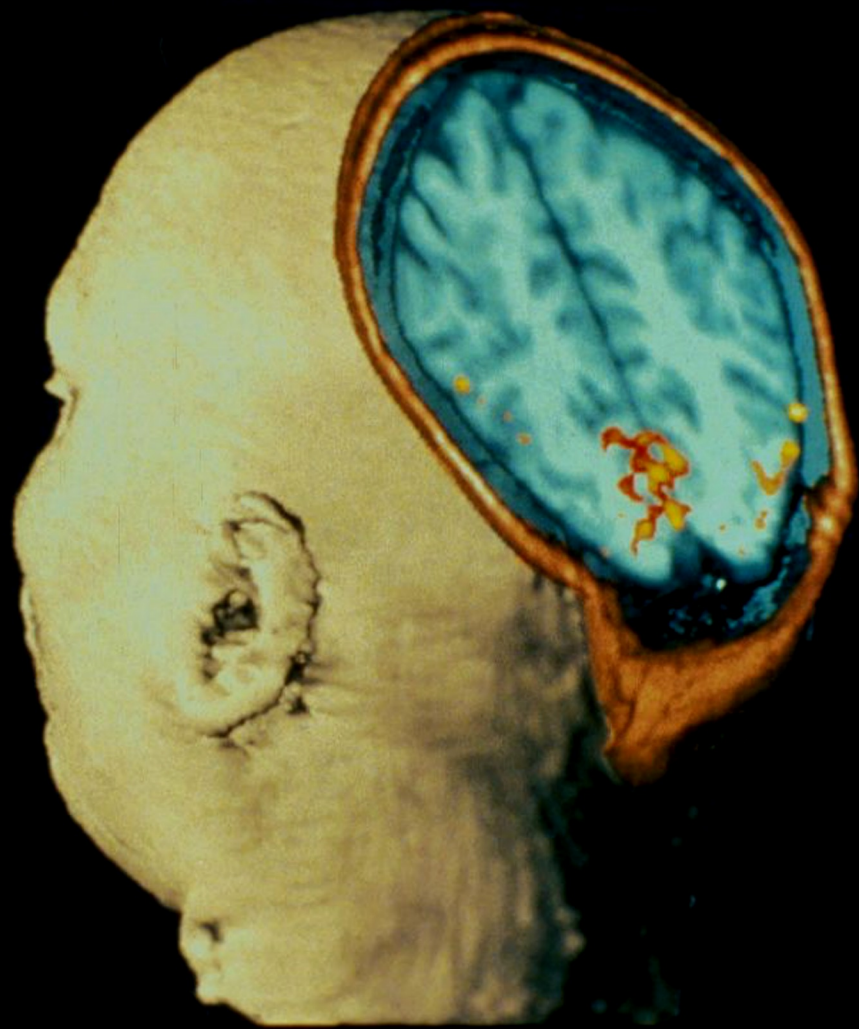

Susceptibility-Induced Field Distortion in the
Vicinity of a Microvessel \perp to B_0 .

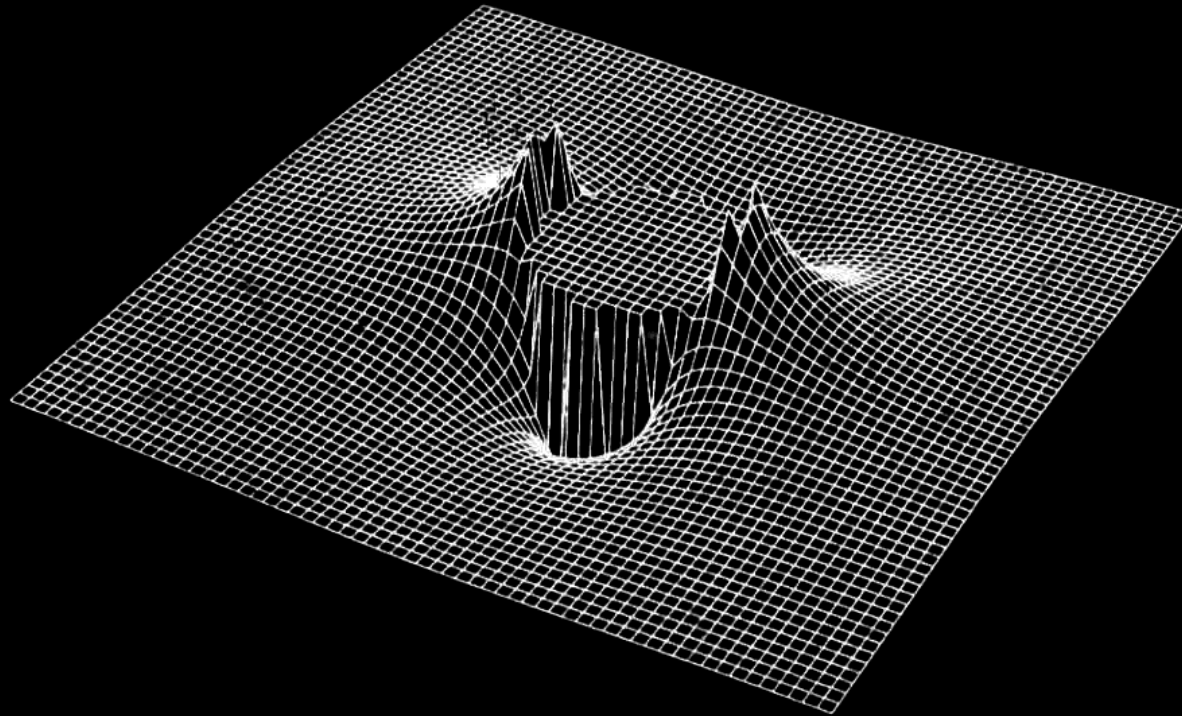


BOLD Contrast in the Detection of Neuronal Activity

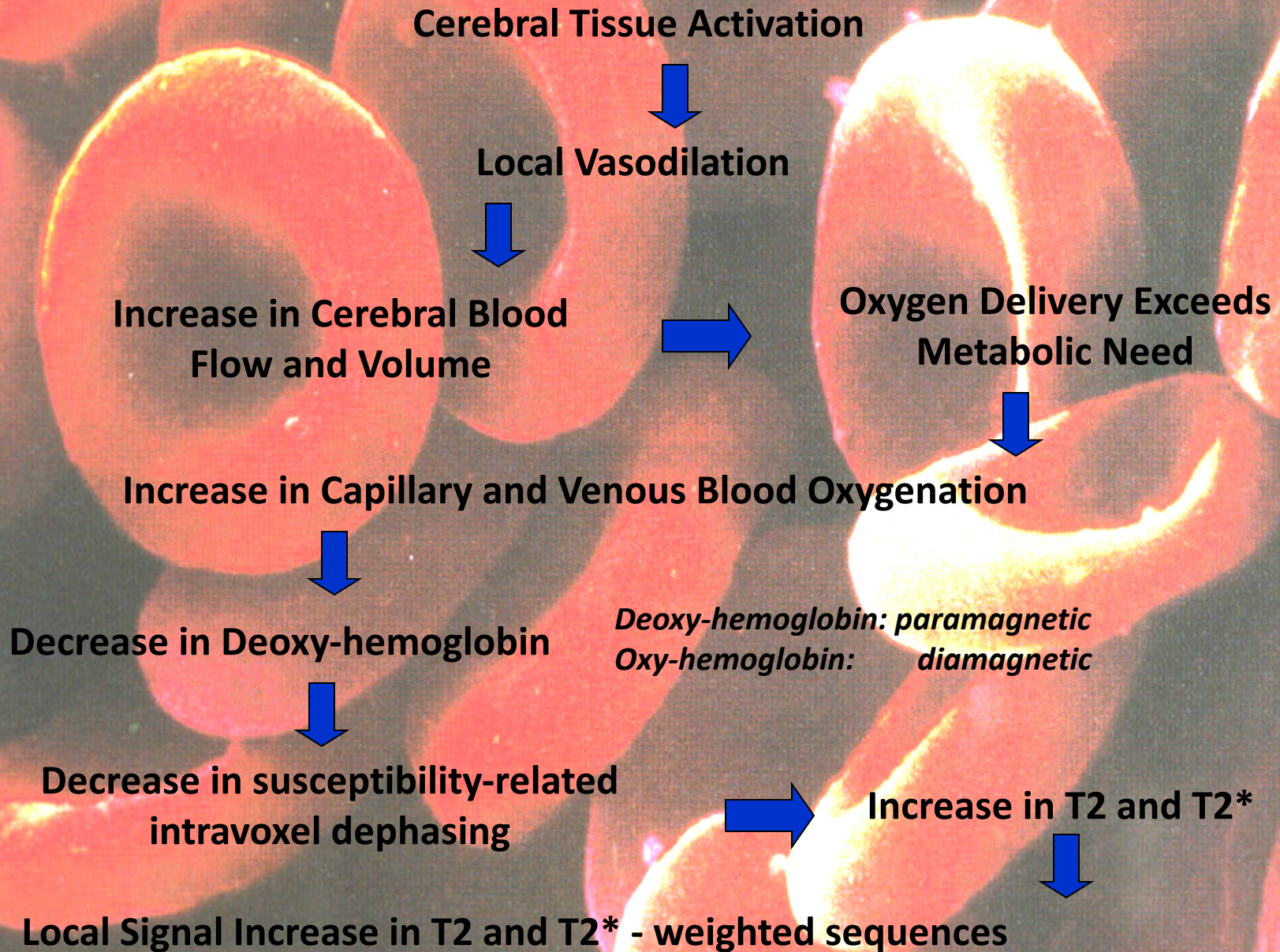


The BOLD Signal

Blood Oxygenation Level Dependent (BOLD) signal changes

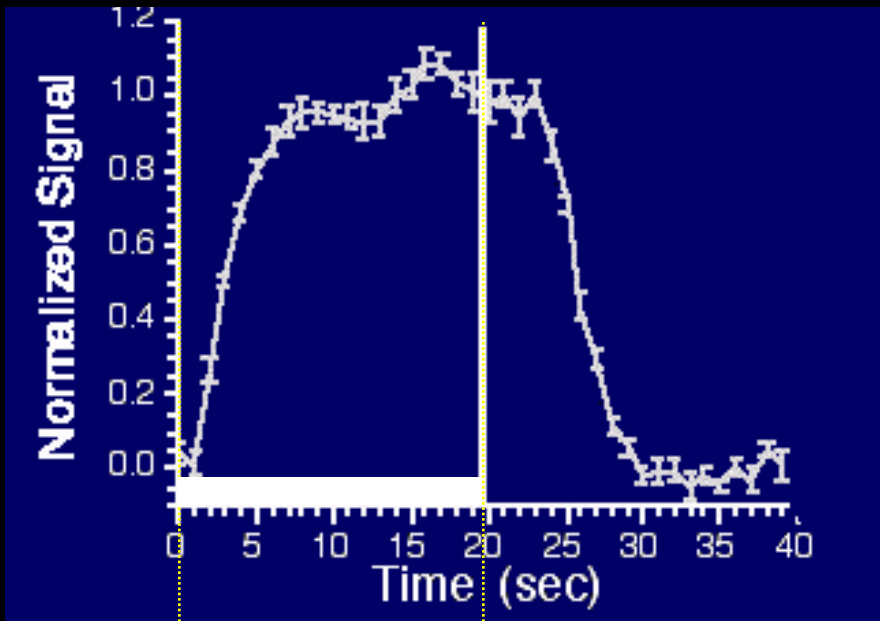

task

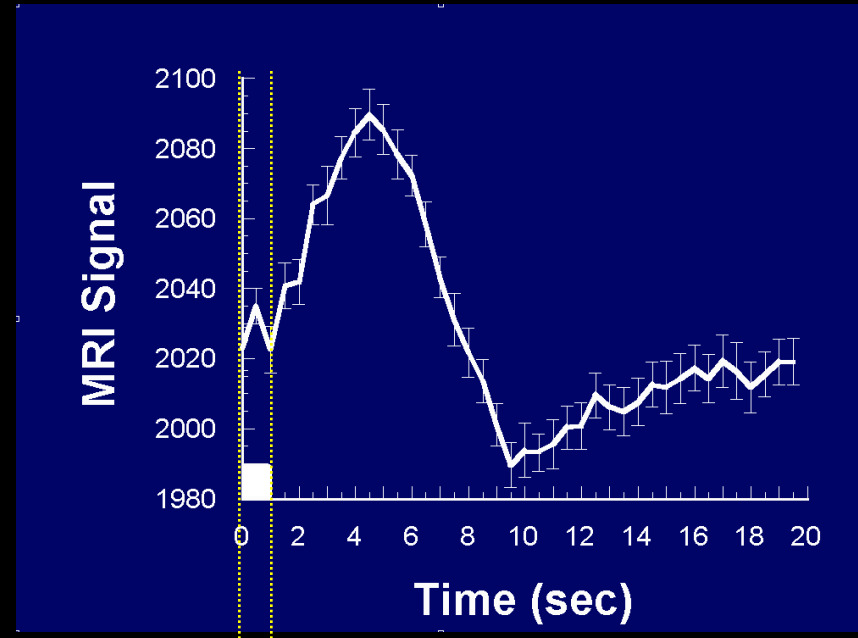

task




Alternating Left and Right Finger Tapping






~ 1992

Perfusion / Flow Imaging

EPISTAR

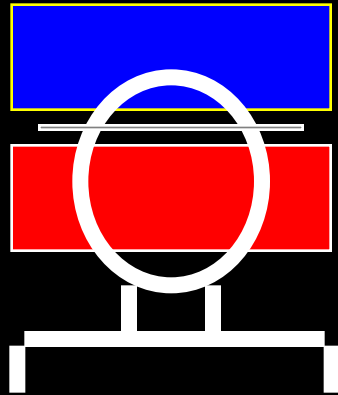

FAIR

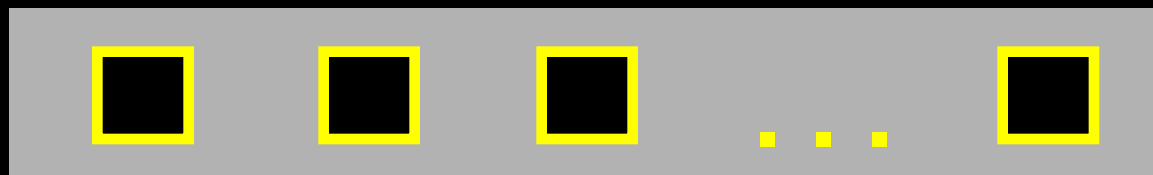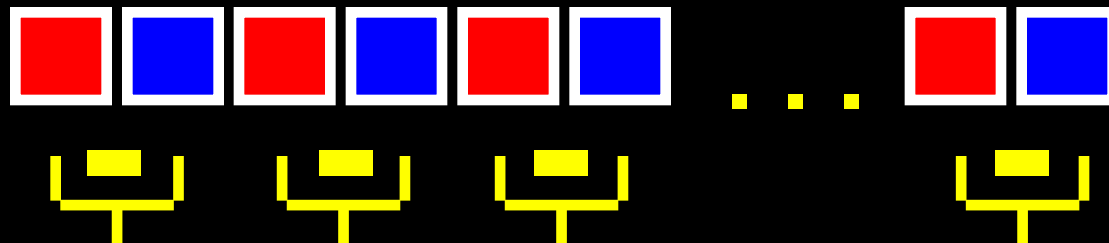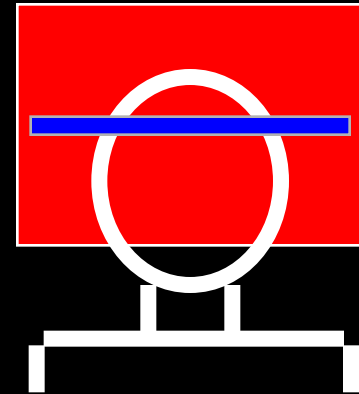

Perfusion
Time Series

TI (ms)

FAIR

EPISTAR

200


400

600


800

1000


1200

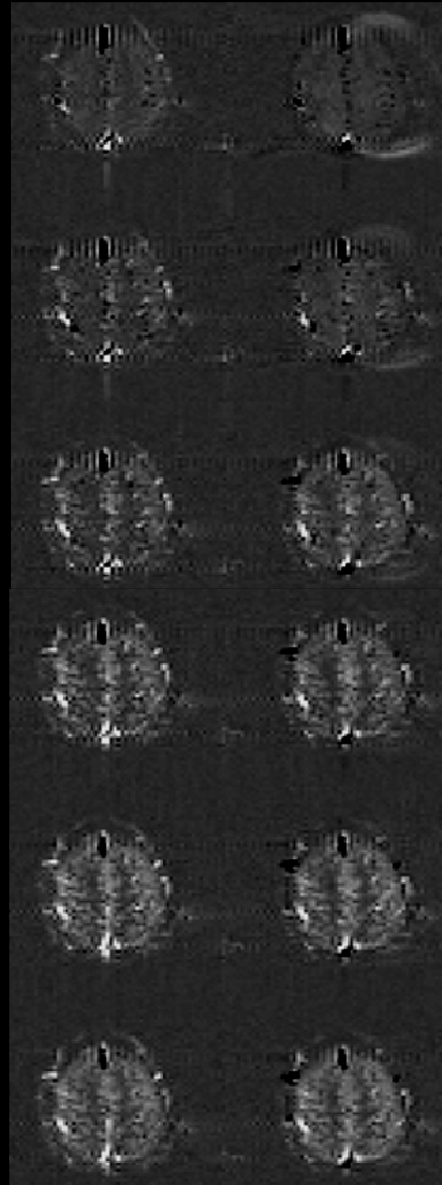


Resting ASL Signal

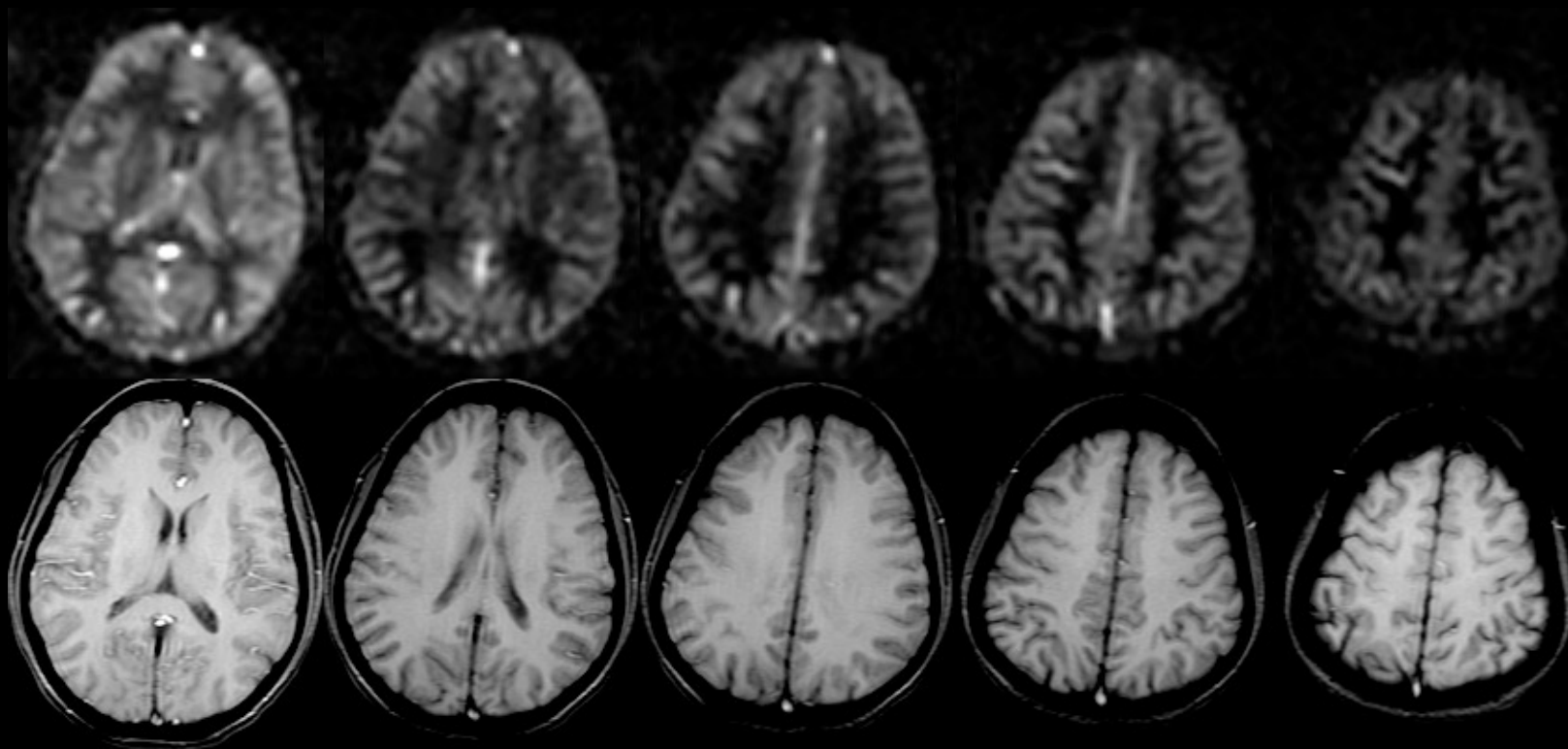

Comparison with Positron Emission Tomography

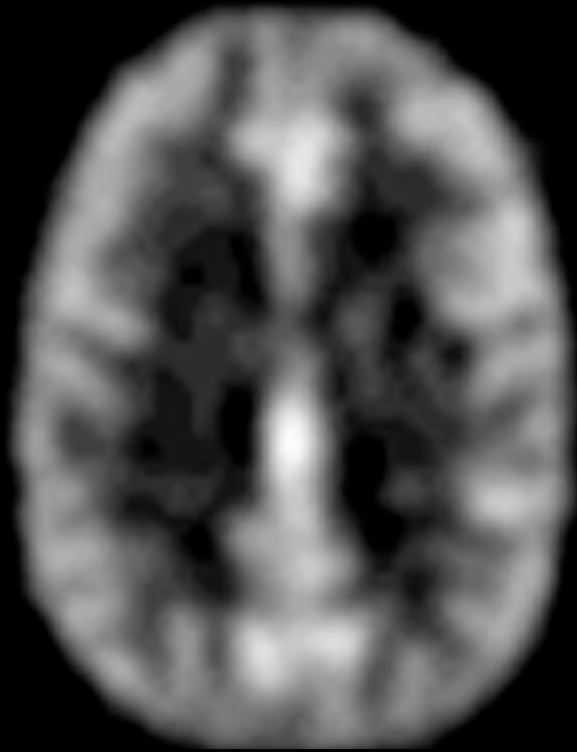

PET: $H_2^{15}O$

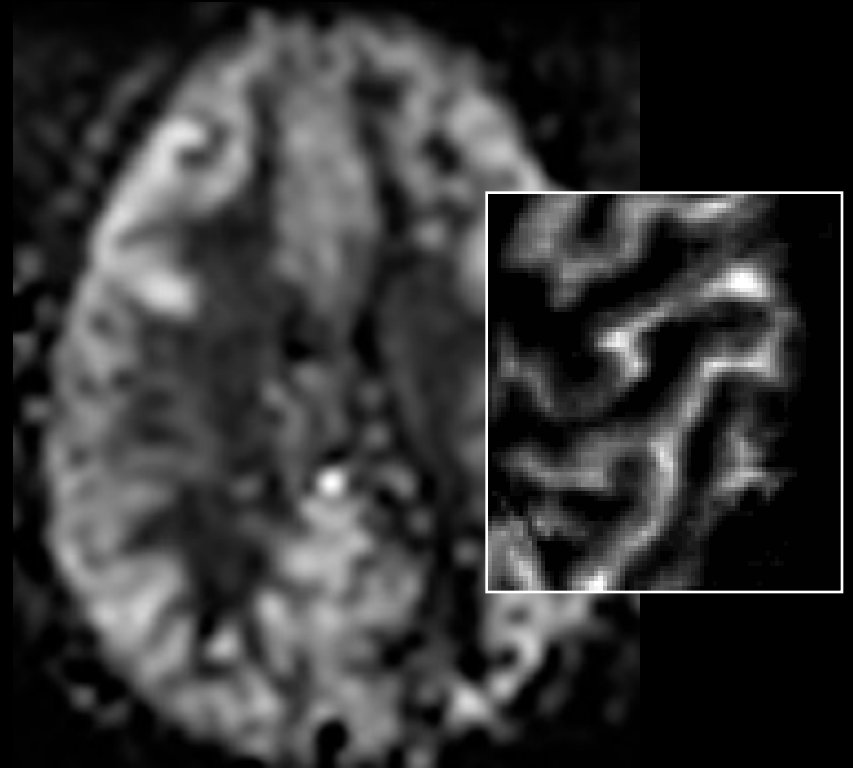





MRI: ASL

Perfusion

BOLD

Rest

Activation

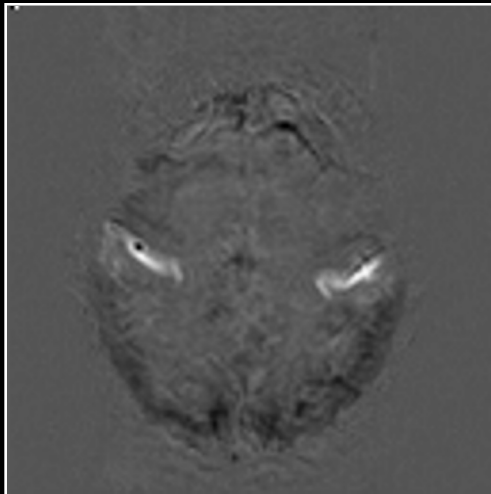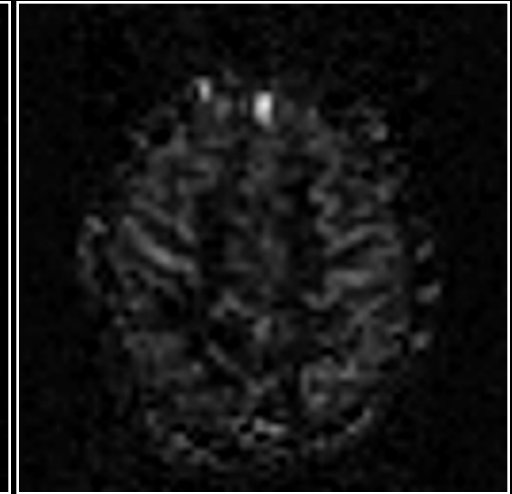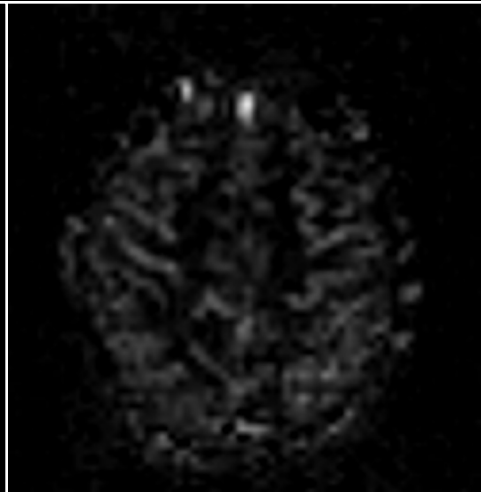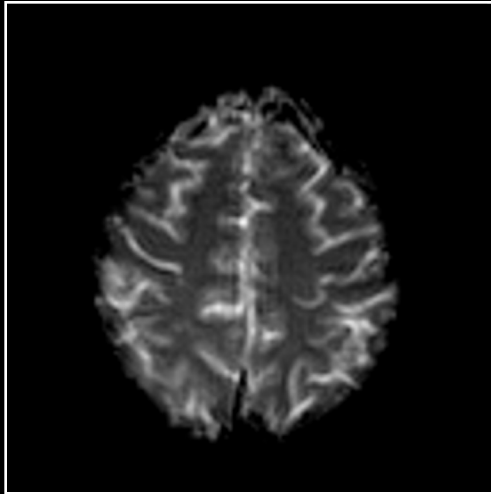

Anatomy




BOLD



Perfusion





CMRO₂-related BOLD signal deficit:

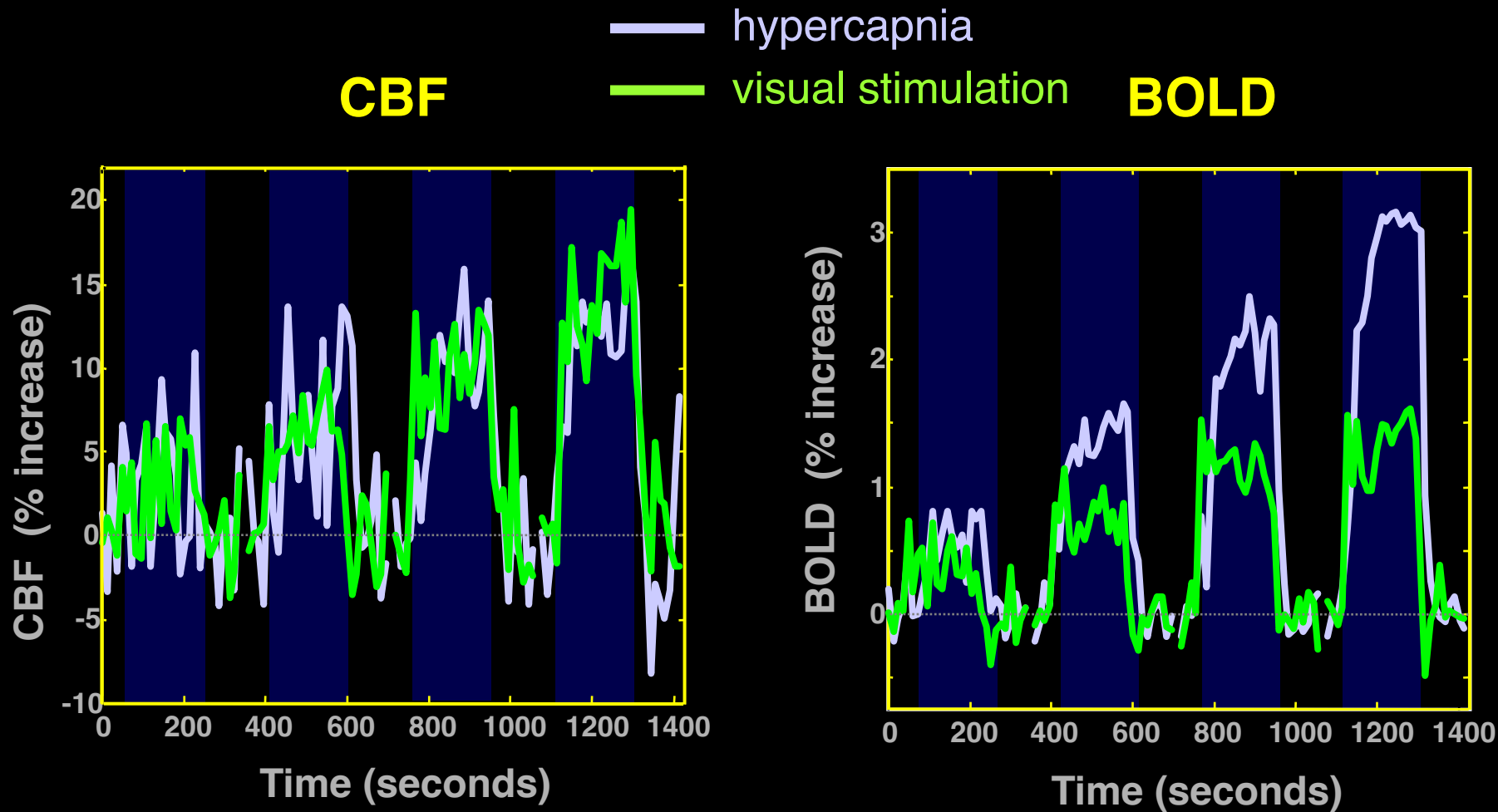

Simultaneous Perfusion and BOLD imaging during graded visual activation and hypercapnia


CBF-CMRO₂ coupling

Hoge, et al.

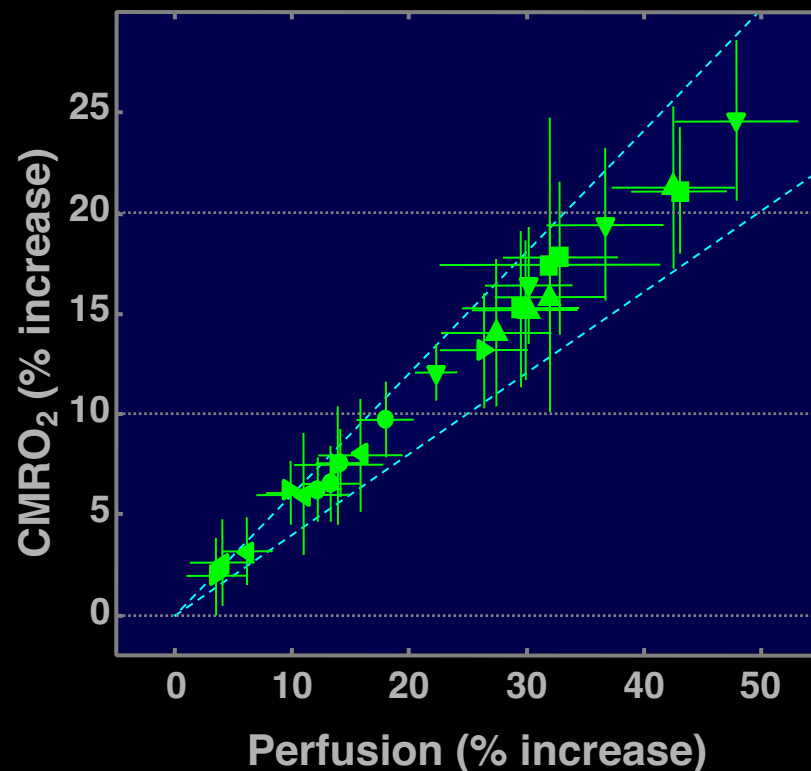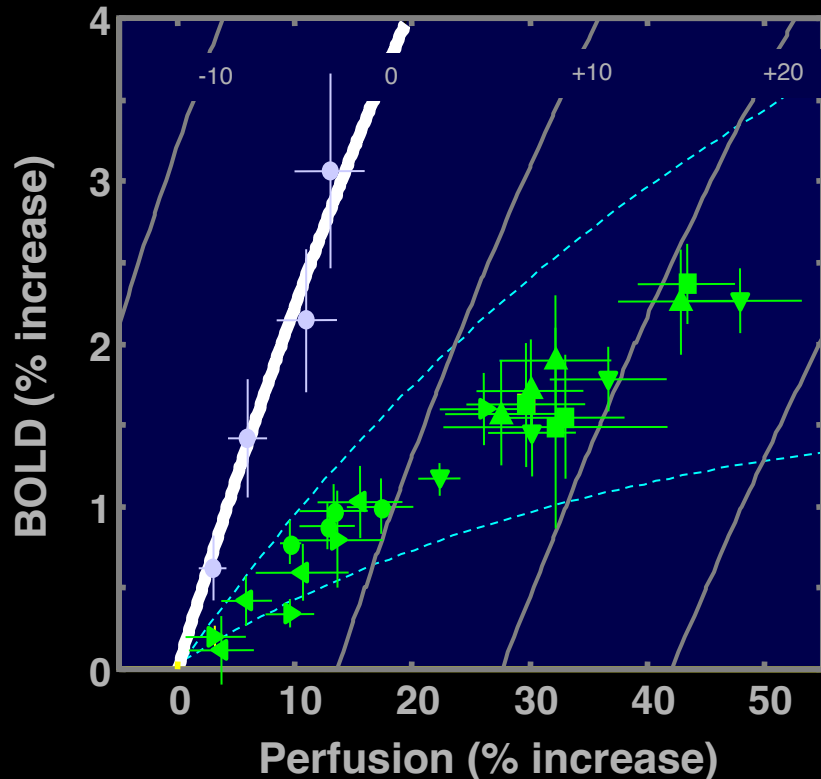


Characterizing Activation-induced CMRO₂ changes using calibration with hypercapnia


Computed CMRO₂ changes

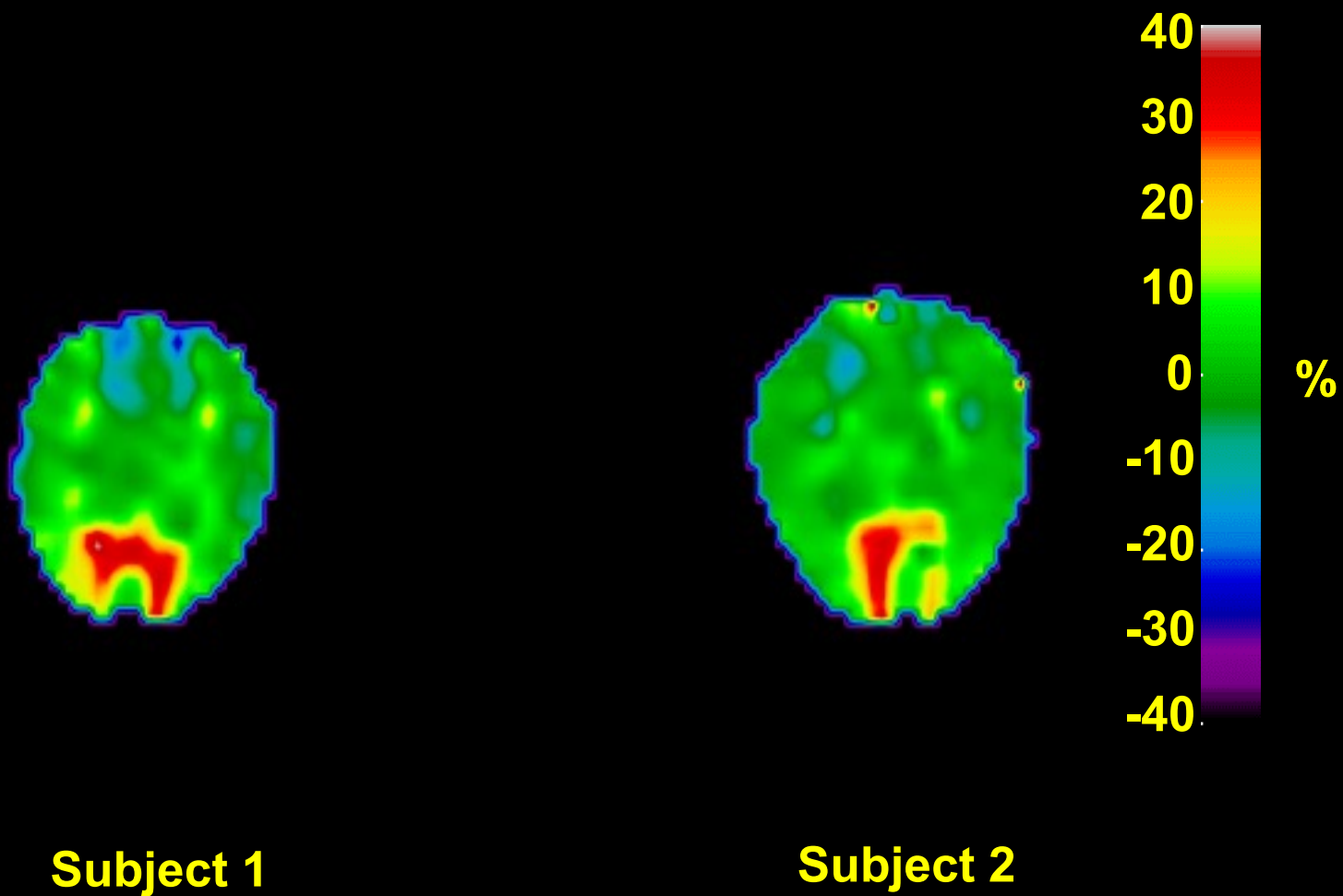



Hemodynamic Specificity

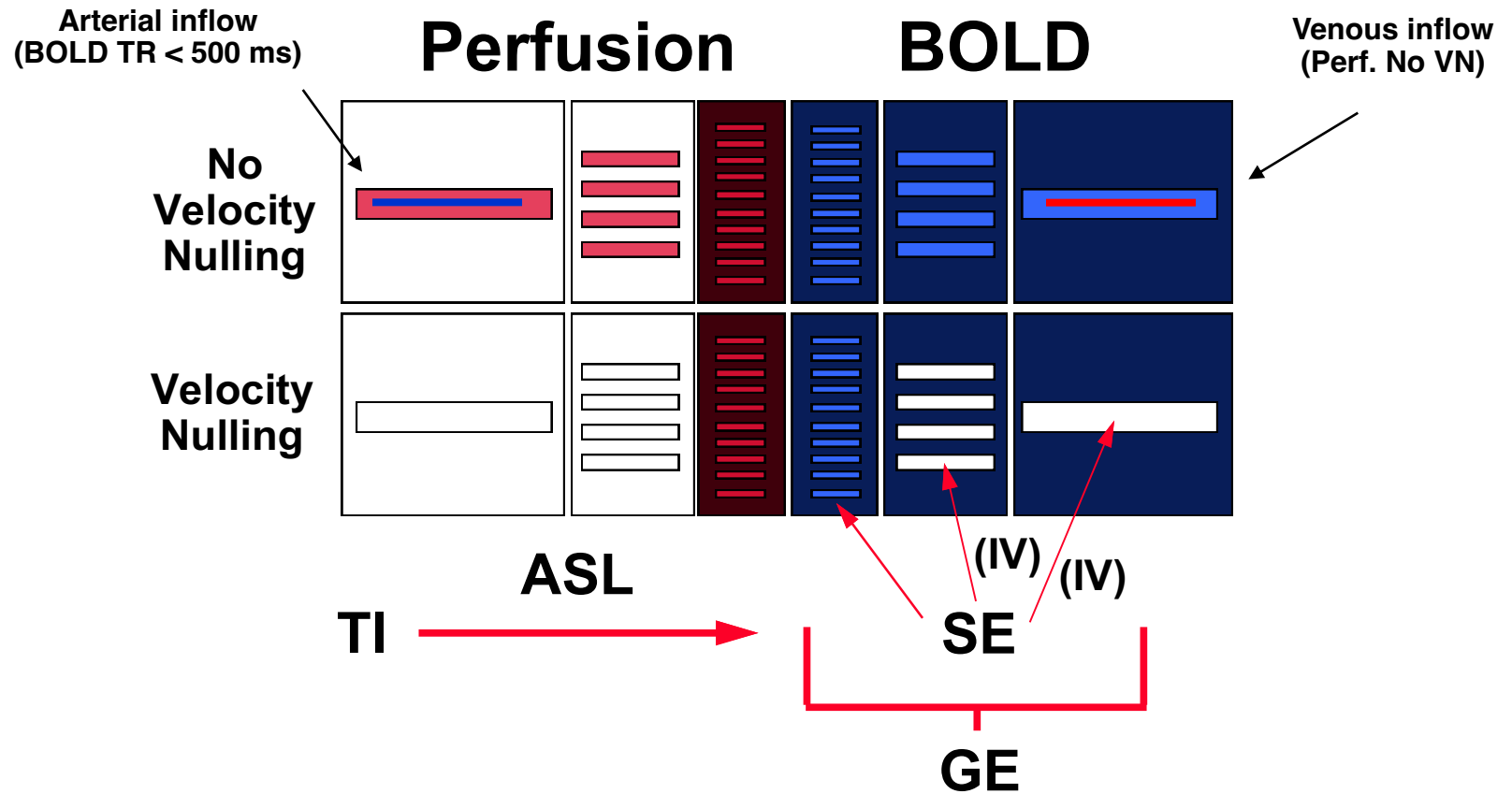

Hemodynamic Stress Calibration

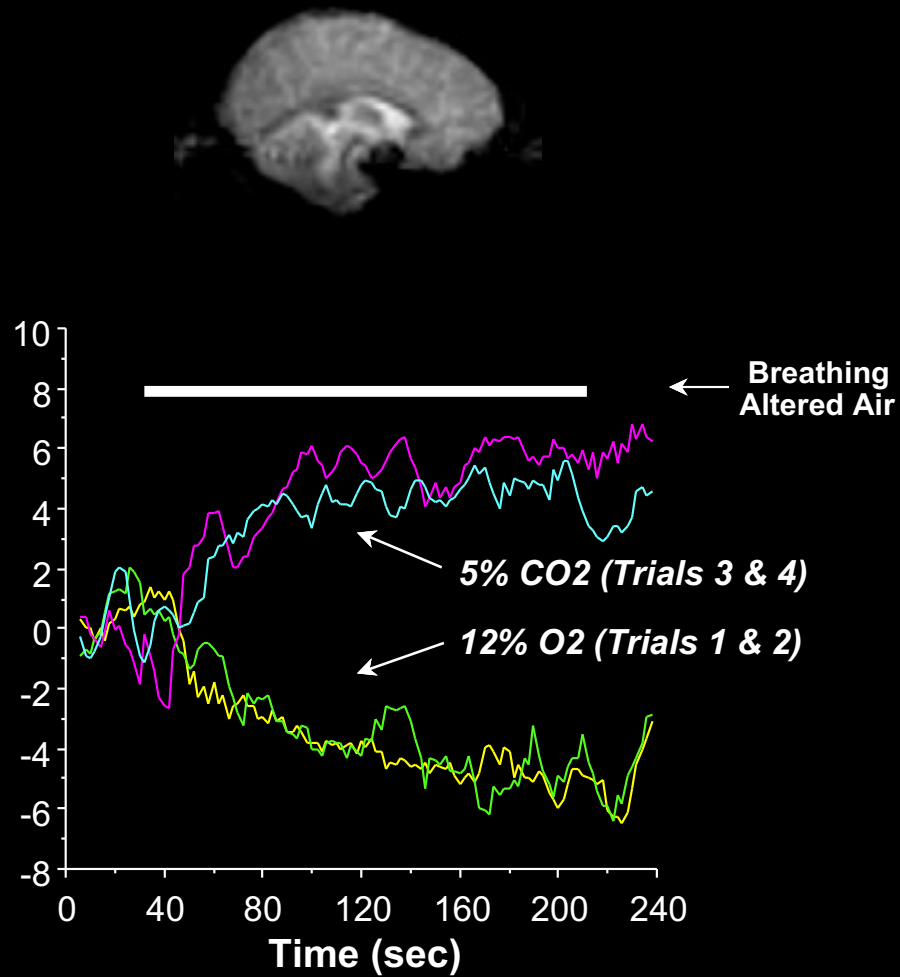


5% CO2

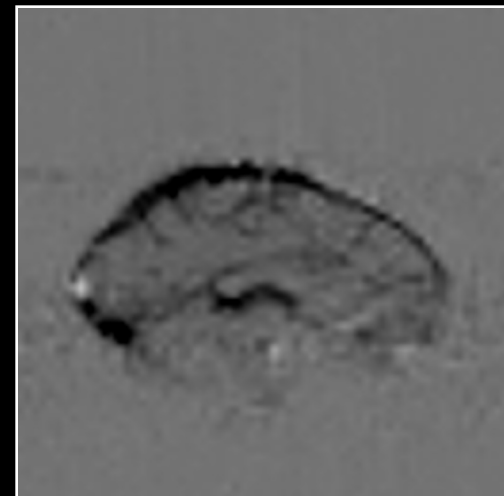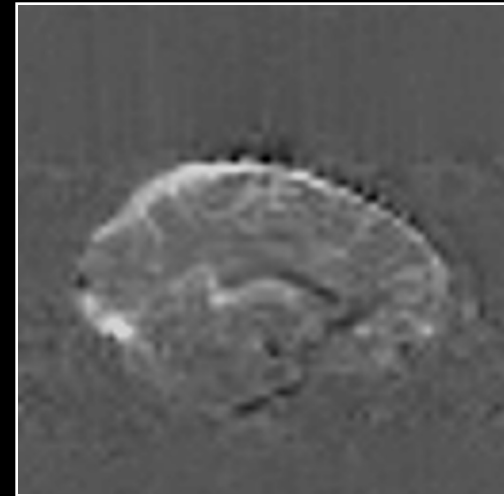

12% O2



CMRO₂-related BOLD signal deficit:

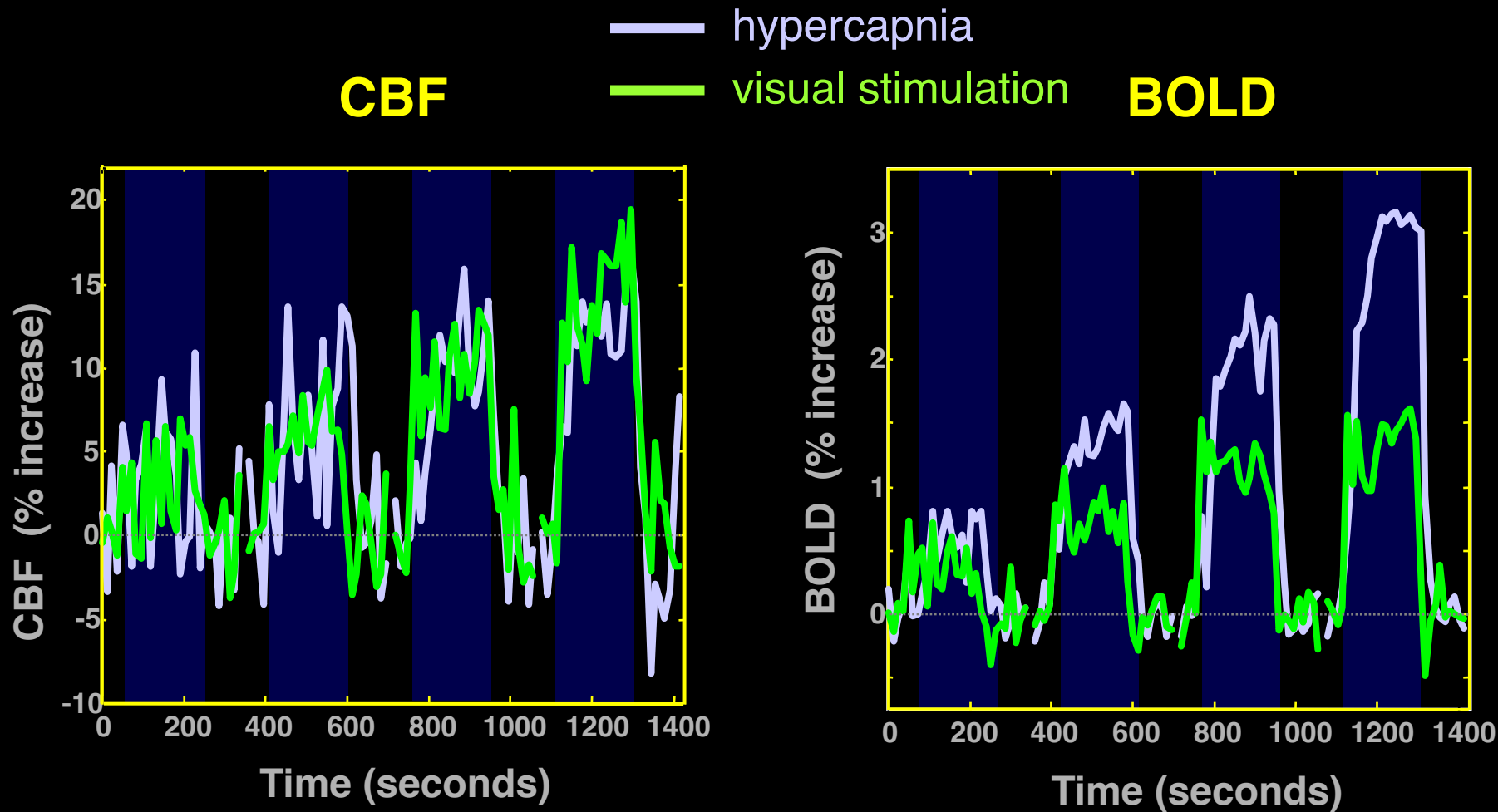

Simultaneous Perfusion and BOLD imaging during graded visual activation and hypercapnia


CBF-CMRO₂ coupling

Hoge, et al.

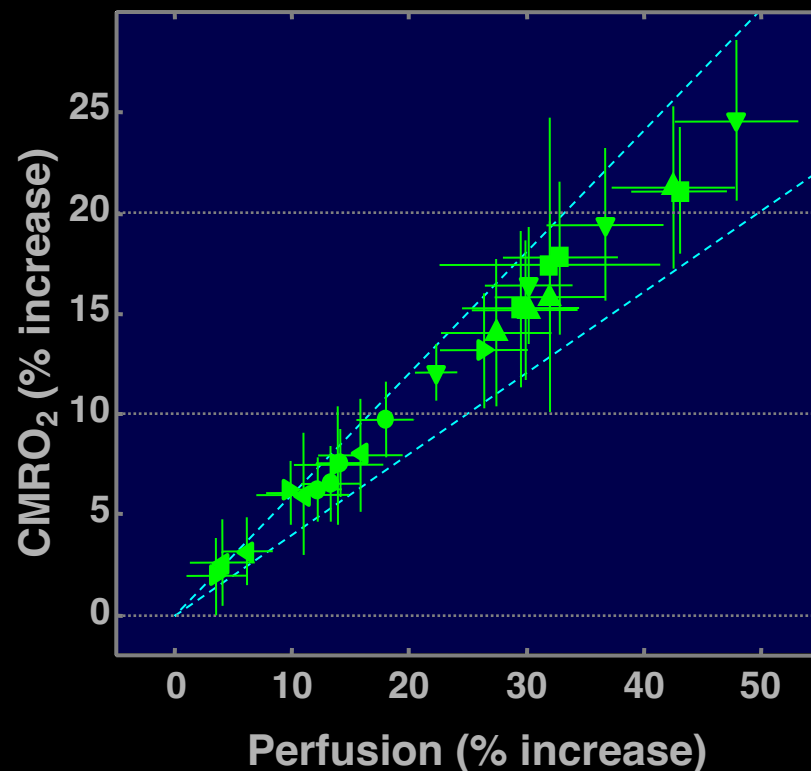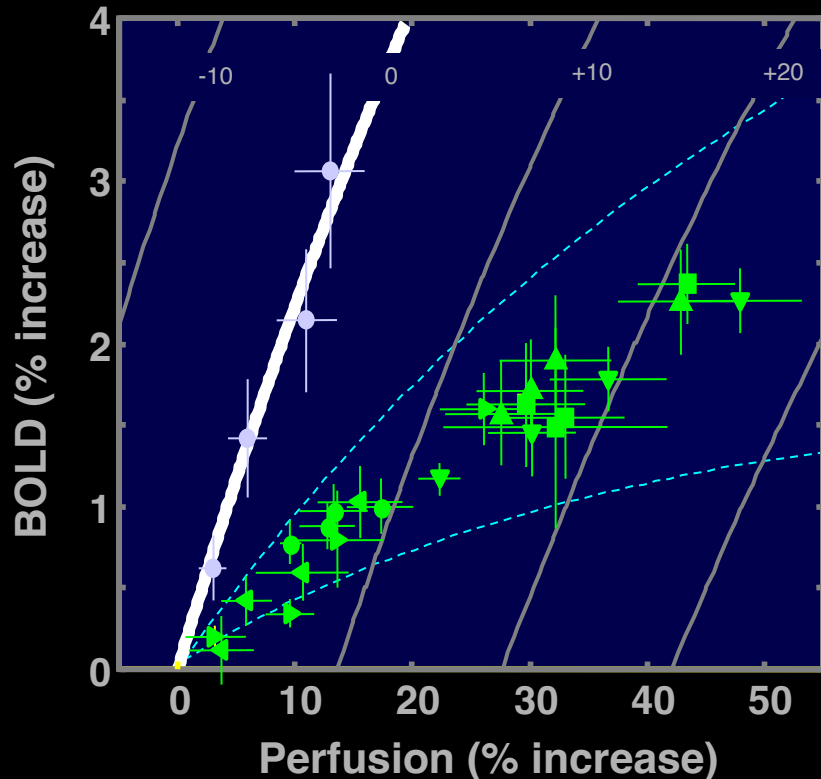

Characterizing Activation-induced CMRO₂ changes using calibration with hypercapnia


Computed CMRO₂ changes

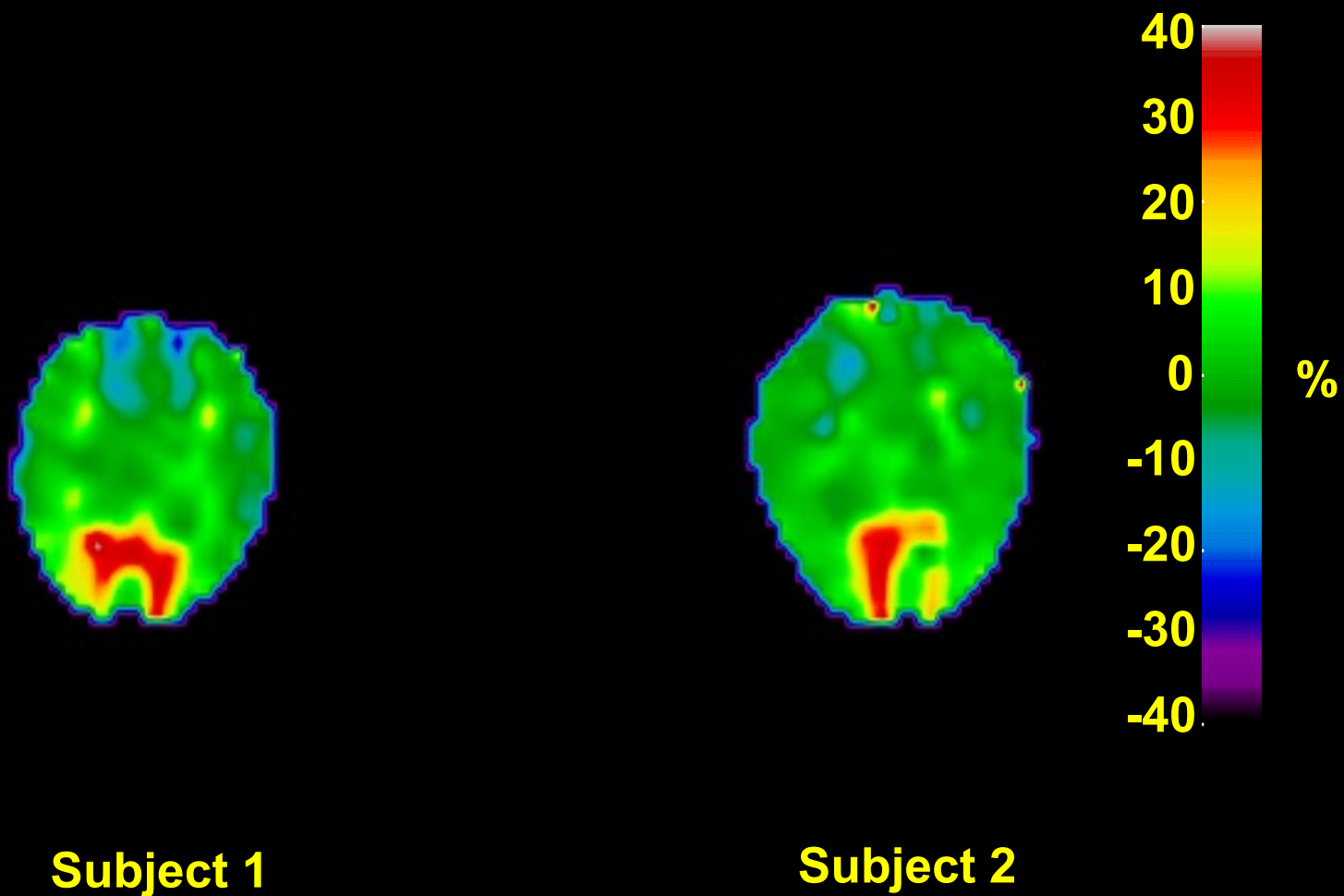


Neuronal Activation Input Strategies

1. Block Design

2. Frequency Encoding

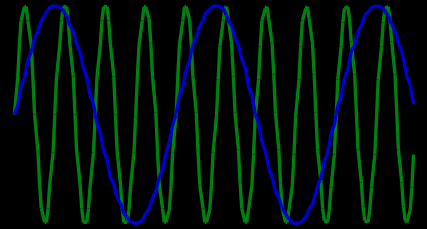


3. Phase Encoding

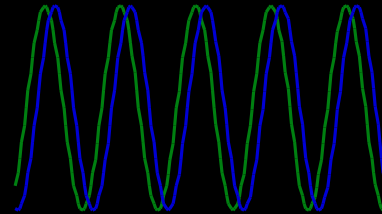

4. Single Event

5. Orthogonal Block Design

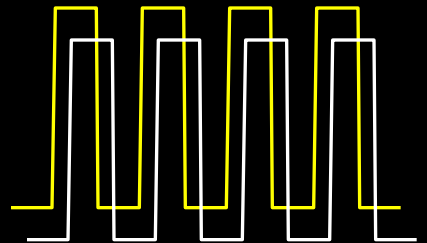


6. Free Behavior Design.


Neuronal Activation Input Strategies

1. Block Design

2. Frequency Encoding

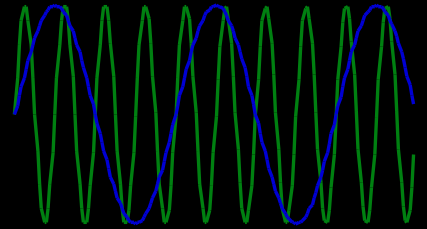


3. Phase Encoding

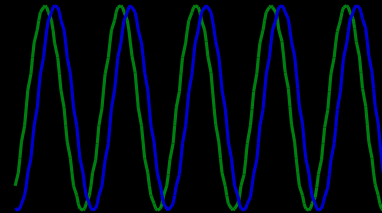



4. Single Event

5. Orthogonal Block Design

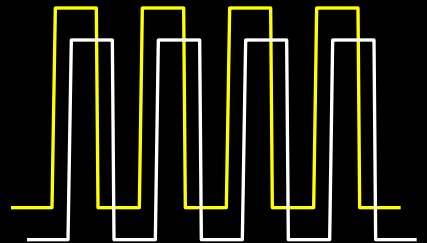


6. Free Behavior Design.


MRI Signal

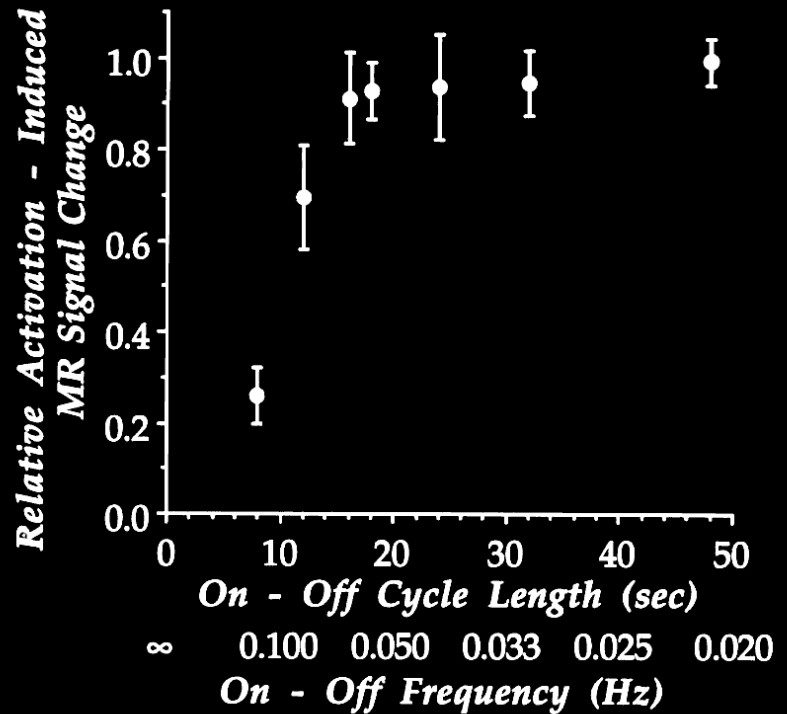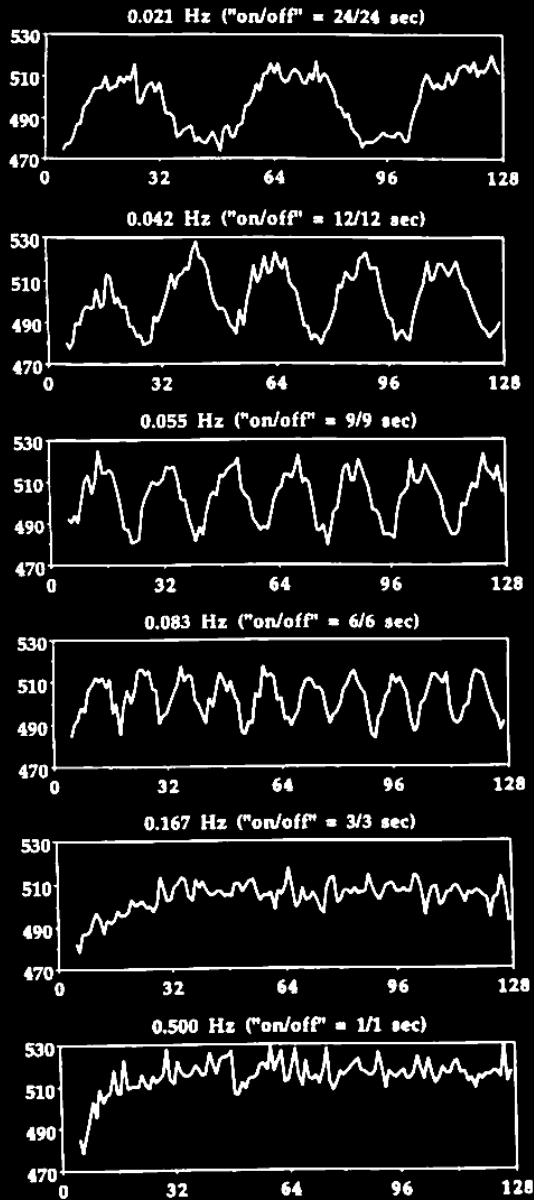


Motor Cortex

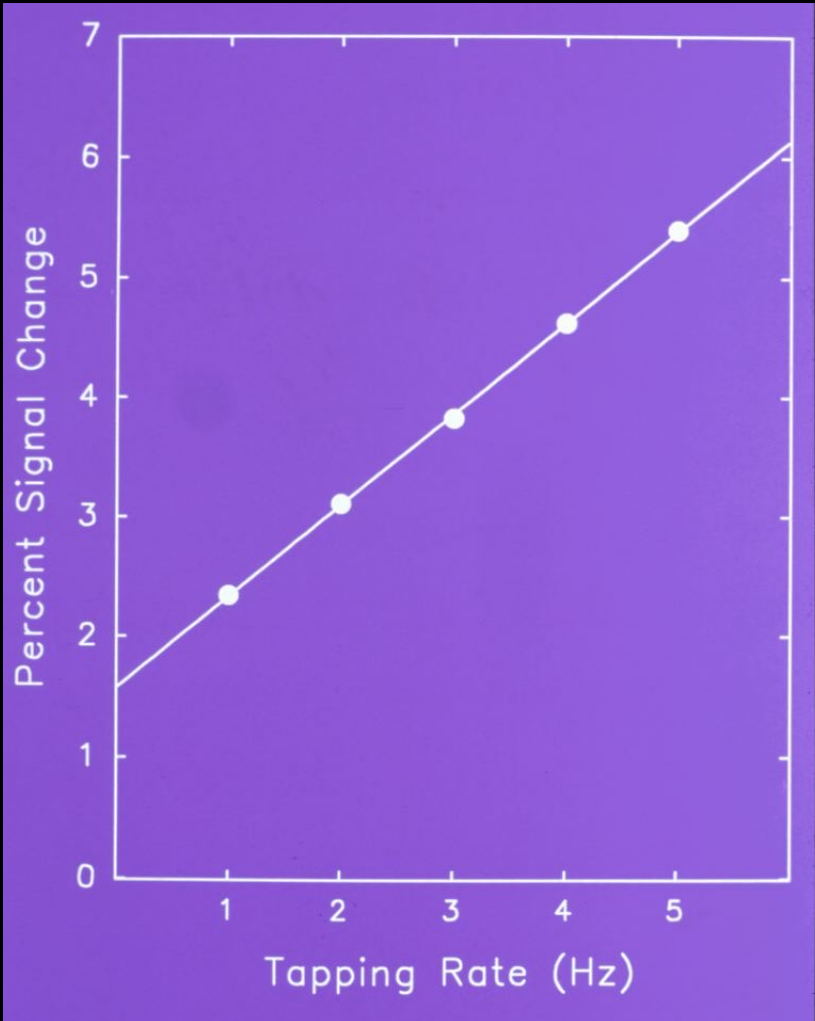



Auditory Cortex

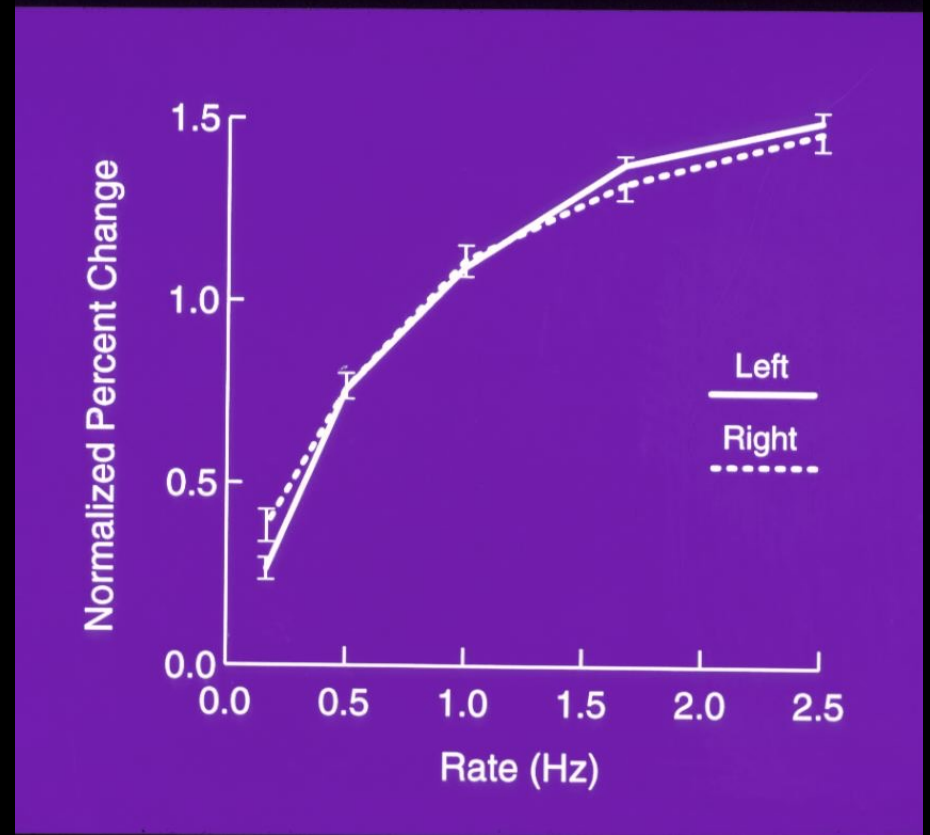

Both the Task and Presentation Rate Affect the fMRI Response

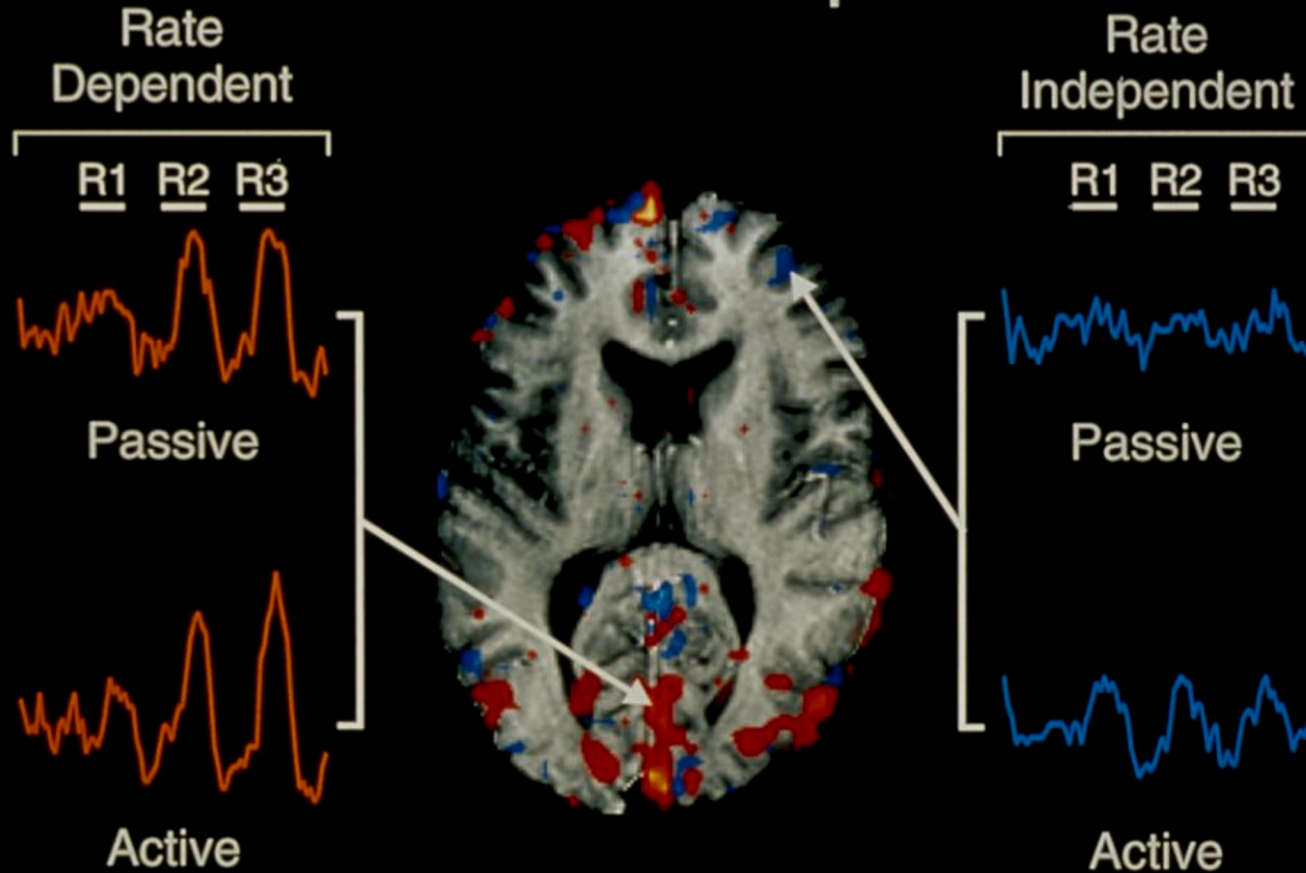

DeYoe et al.

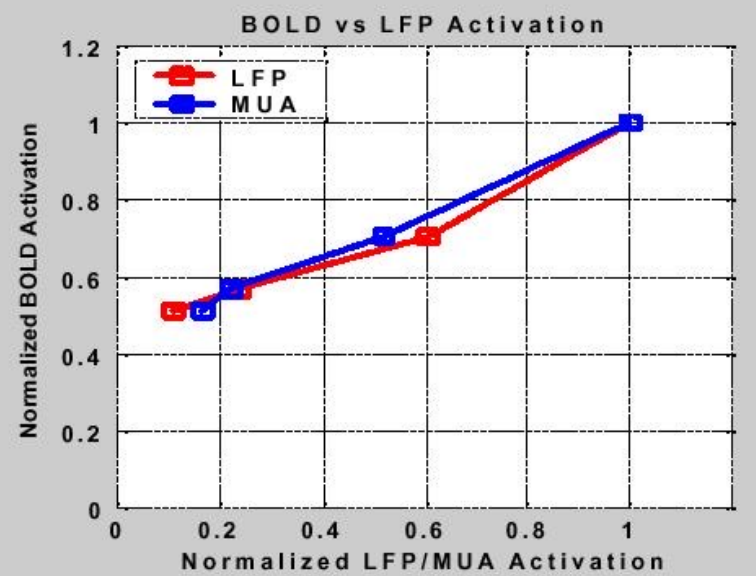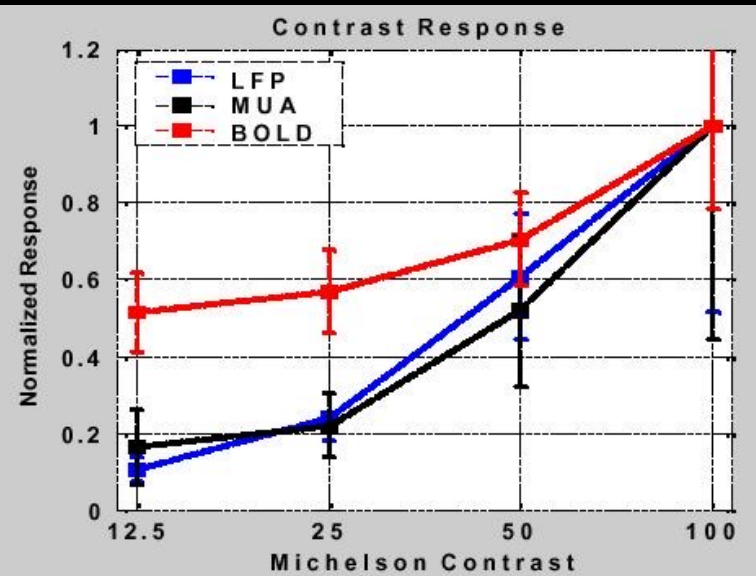


Logothetis et al. Nature, 412, 150-157


Neuronal Activation Input Strategies

1. Block Design

2. Frequency Encoding

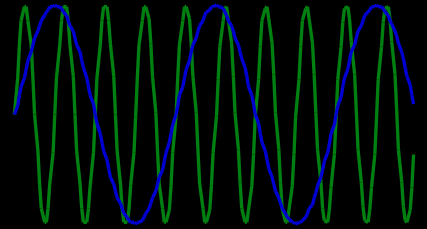


3. Phase Encoding



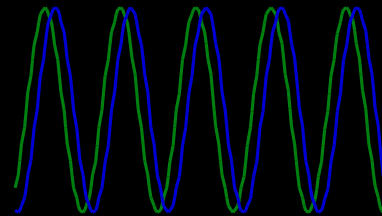


4. Single Event

5. Orthogonal Block Design

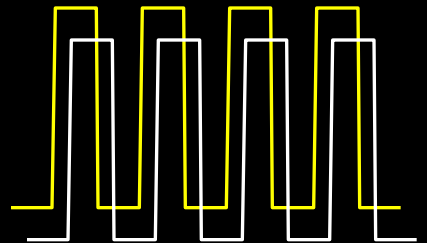


6. Free Behavior Design.

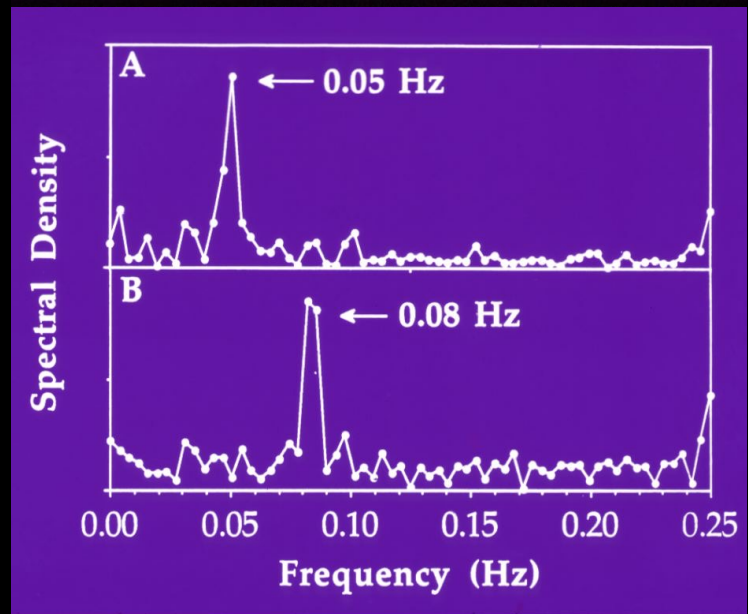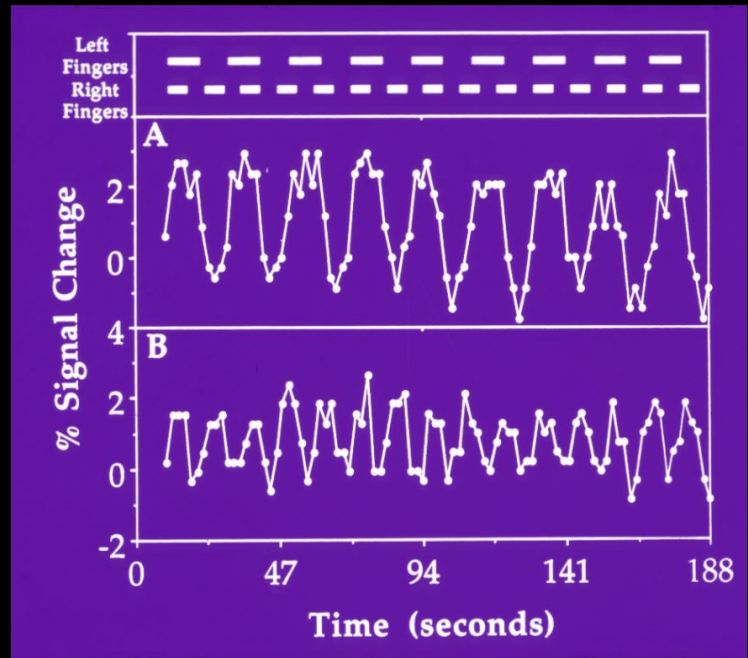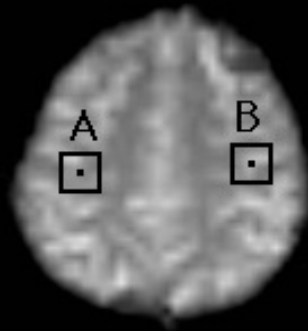



0.08 Hz

0.05 Hz

**spectral
density**

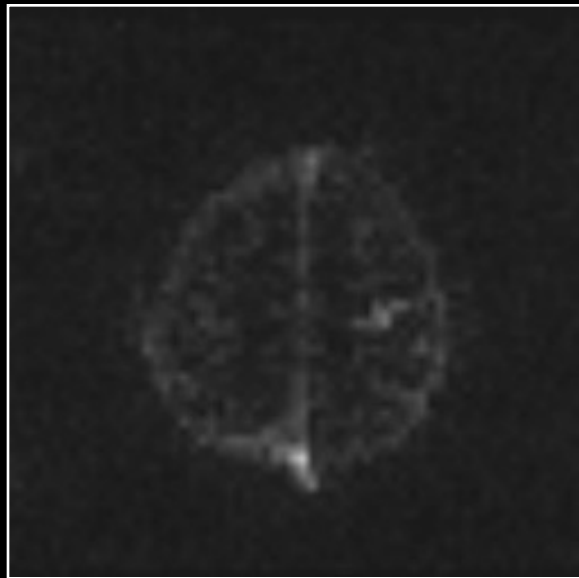

**c.c. > 0.5
with spectra**

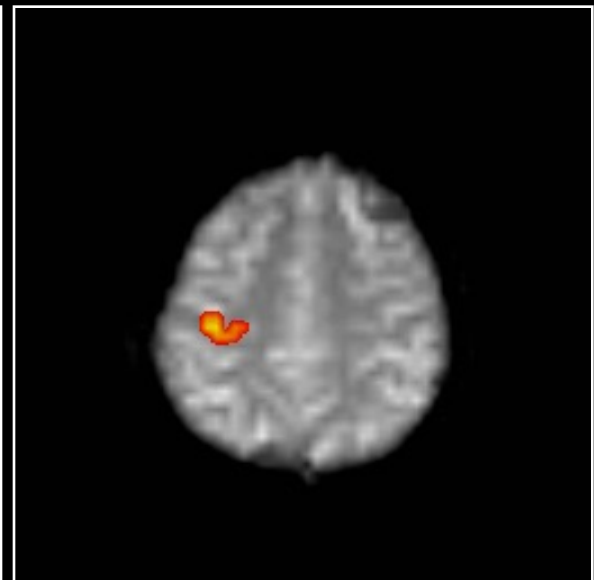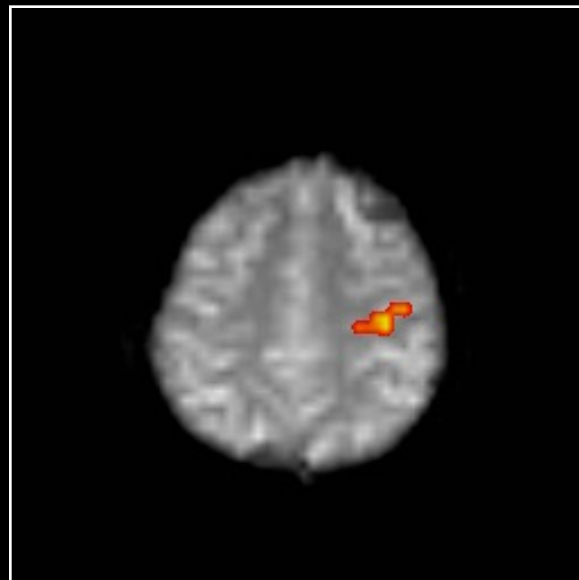

Neuronal Activation Input Strategies

1. Block Design



2. Frequency Encoding







3. Phase Encoding







4. Single Event

5. Orthogonal Block Design

6. Free Behavior Design.

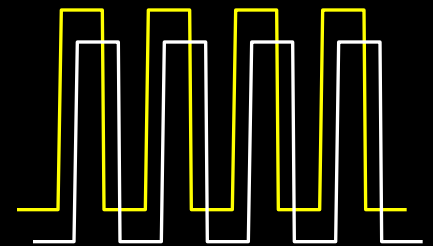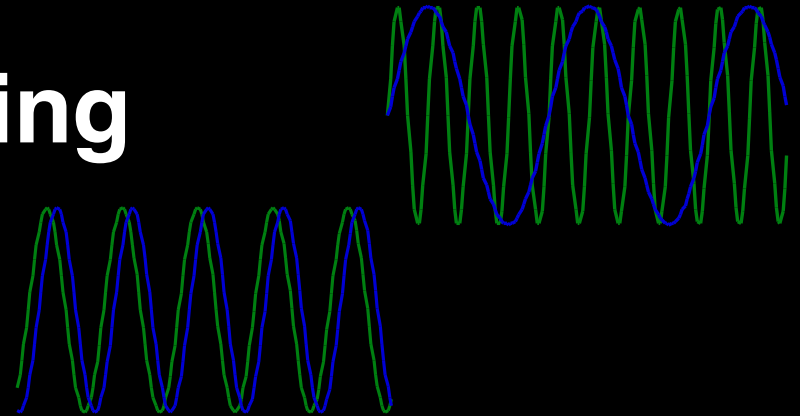

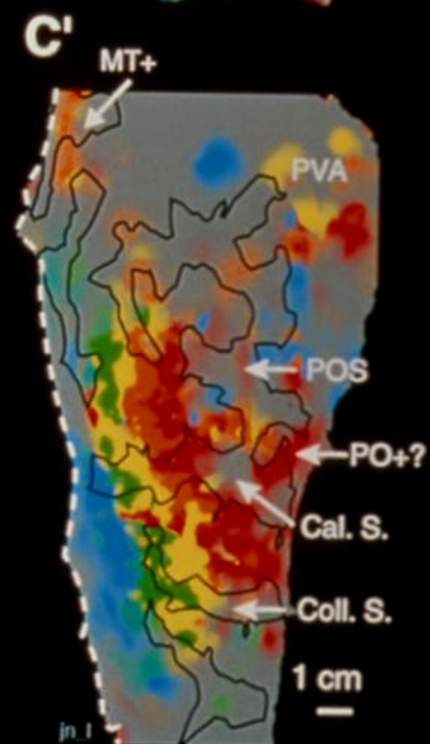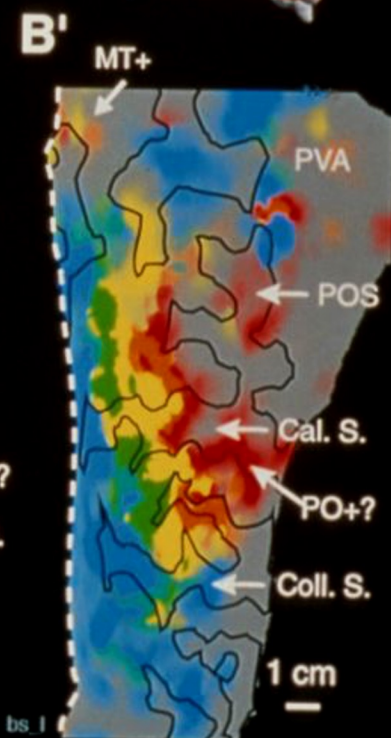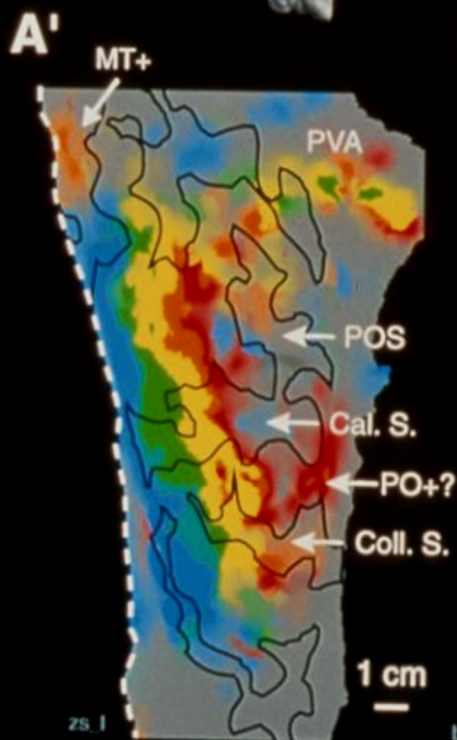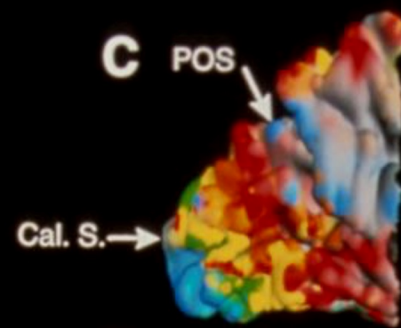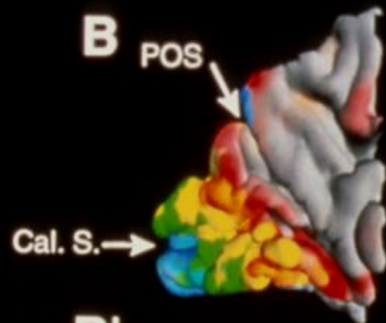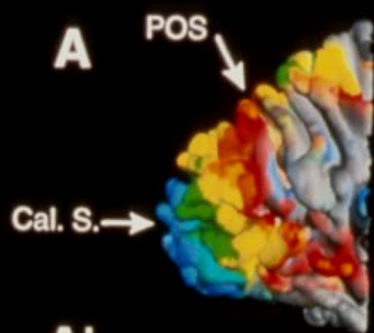


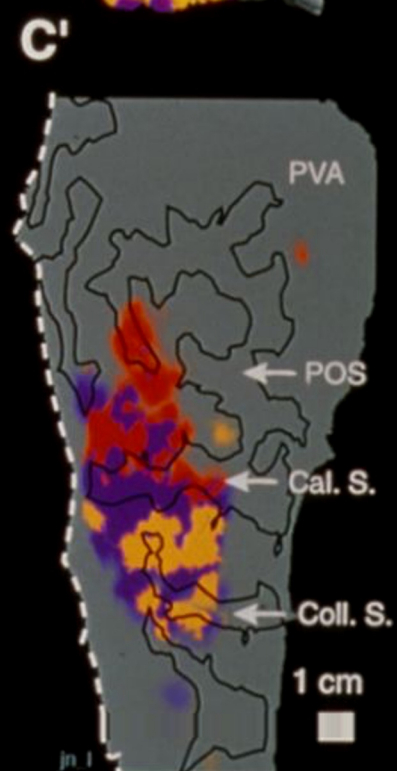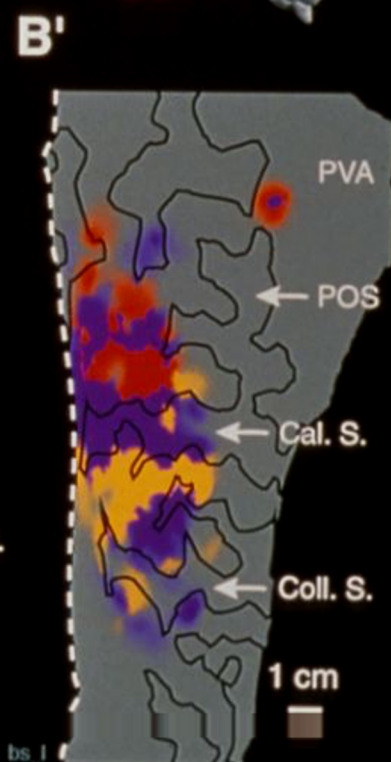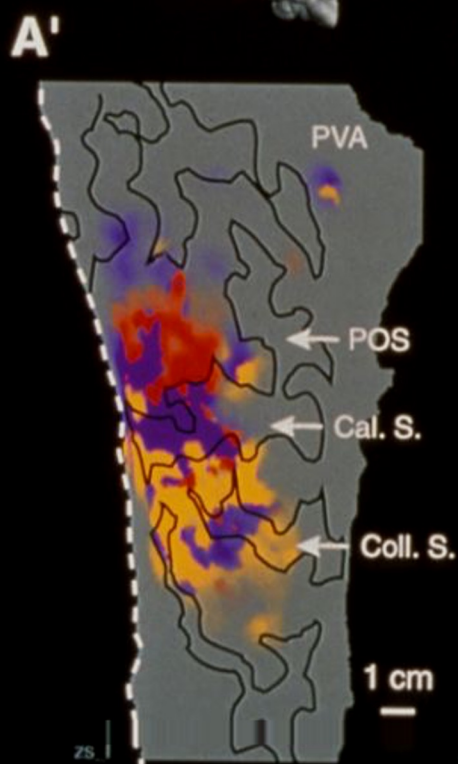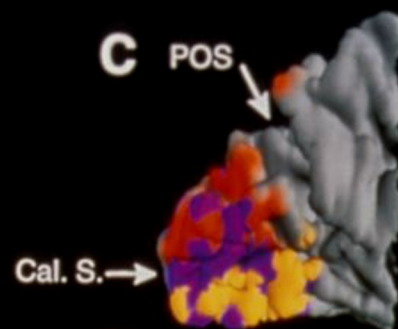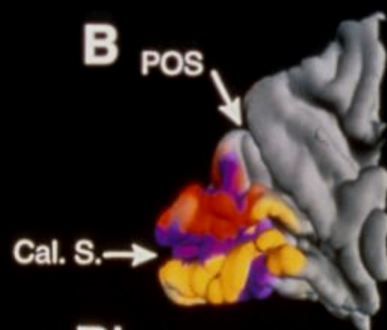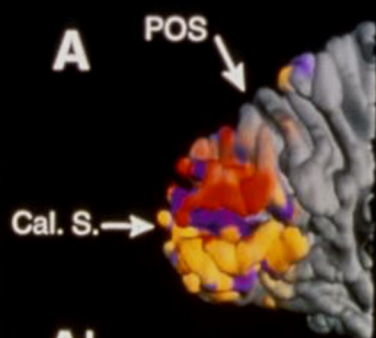


Neuronal Activation Input Strategies

1. Block Design

2. Frequency Encoding

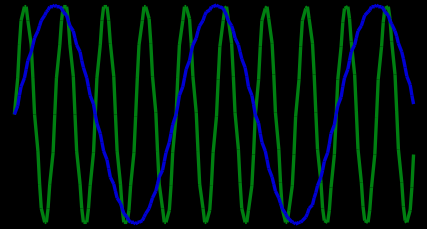


3. Phase Encoding

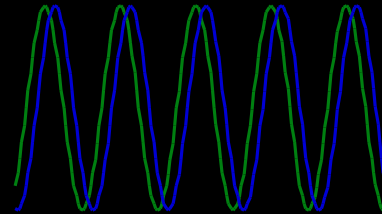


4. Single Event


5. Orthogonal Block Design

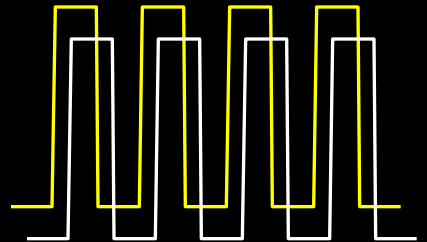


6. Free Behavior Design.

Motor Cortex

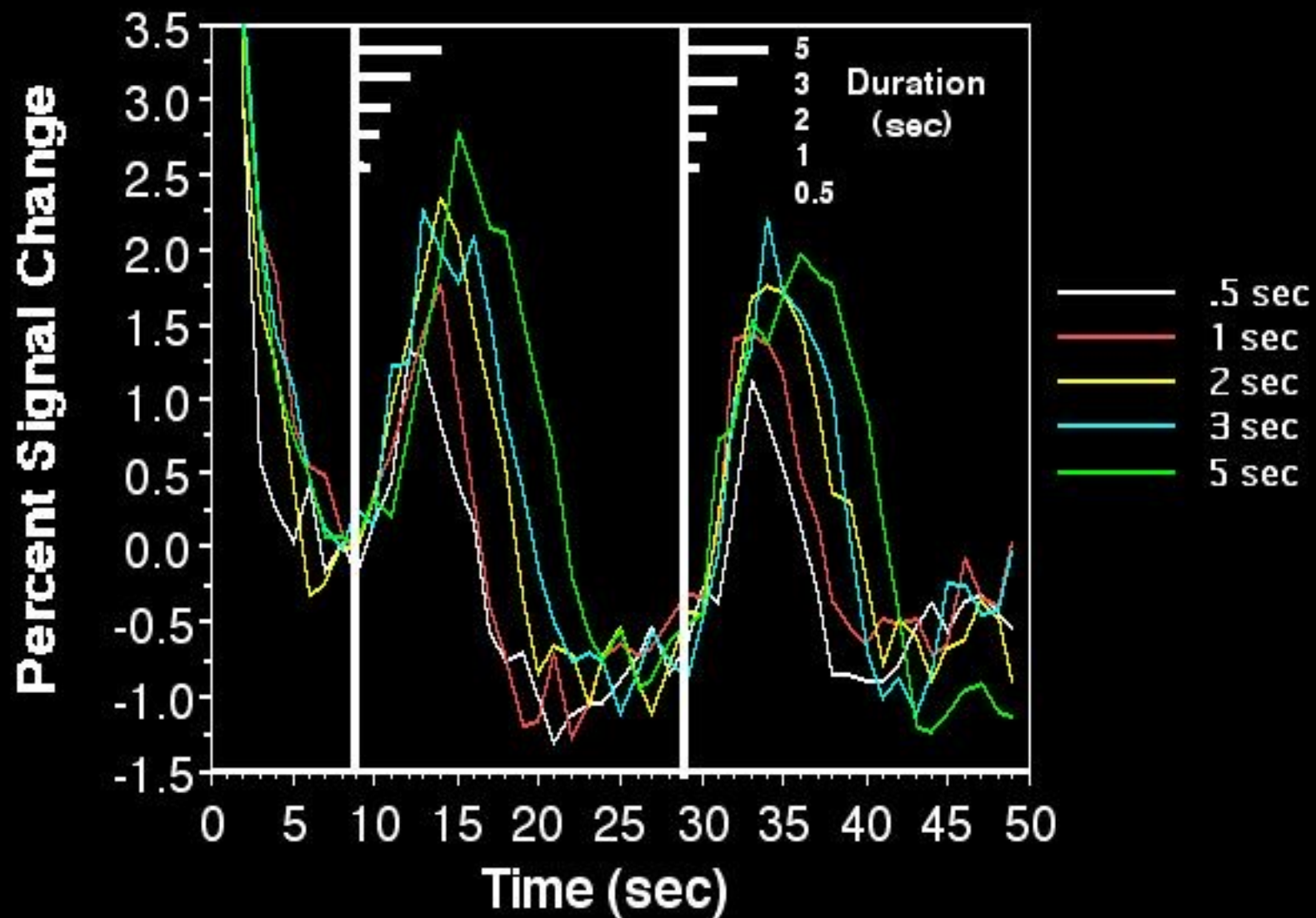



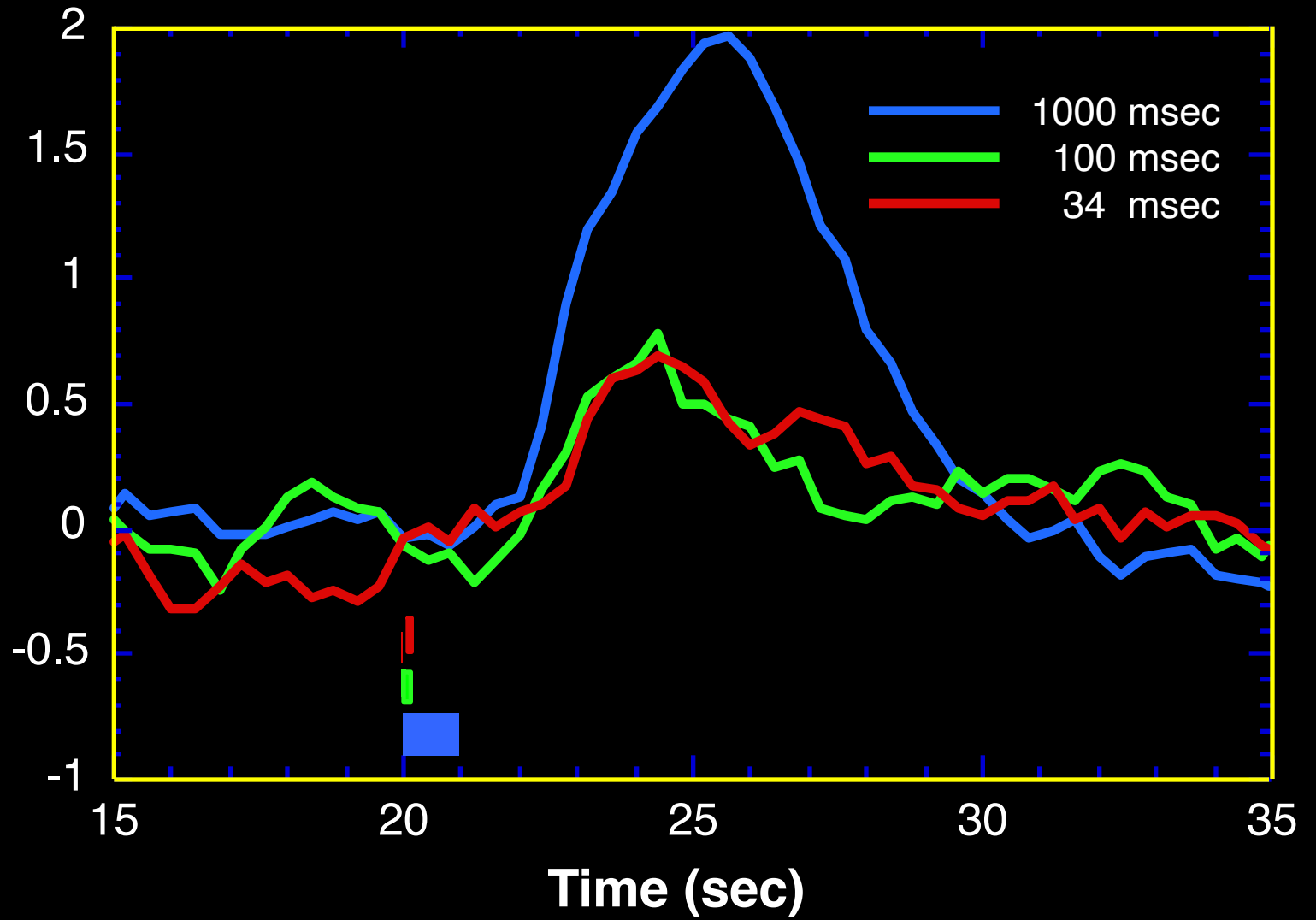


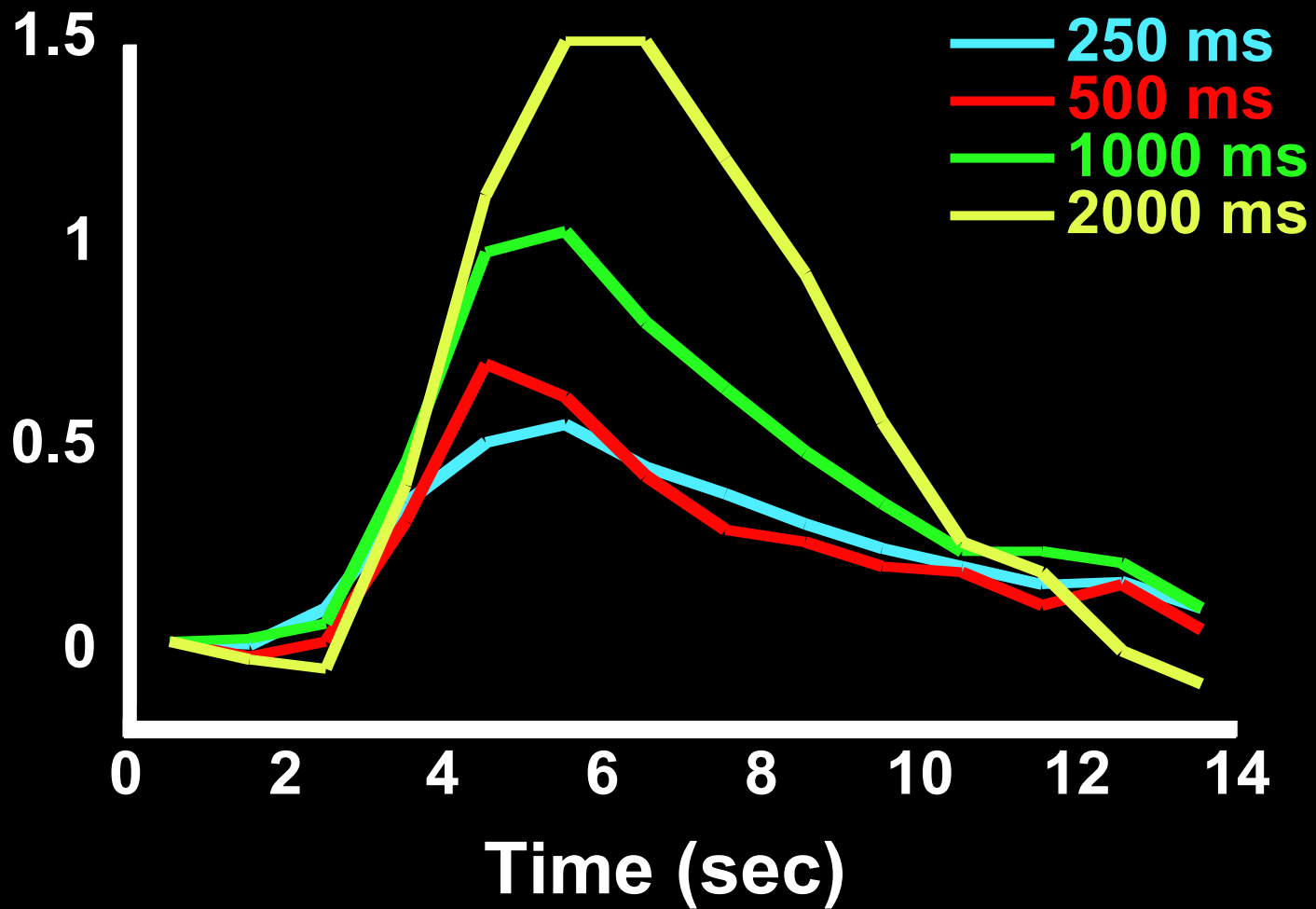


BLOCKED:



SINGLE TRIAL:


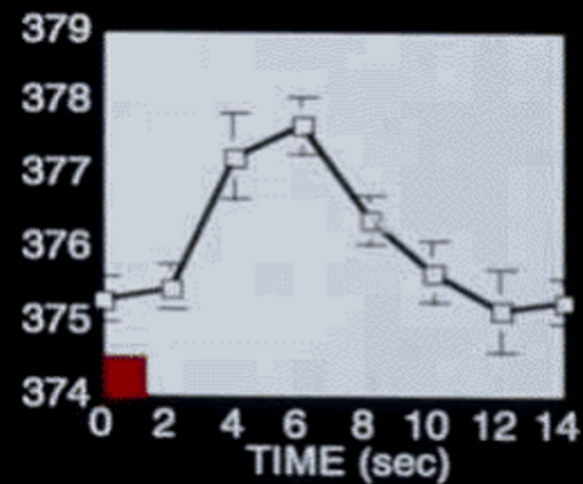
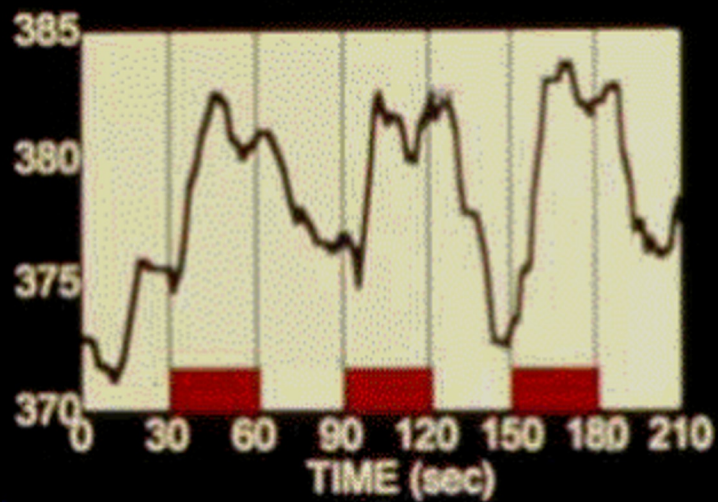
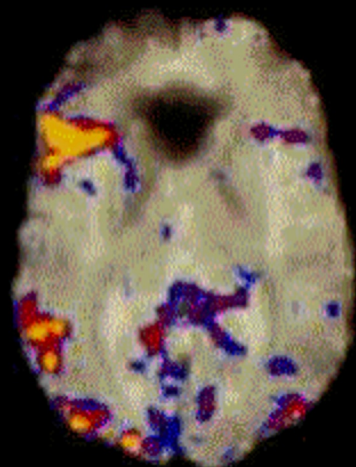


Time Course Comparison Across Brain Regions

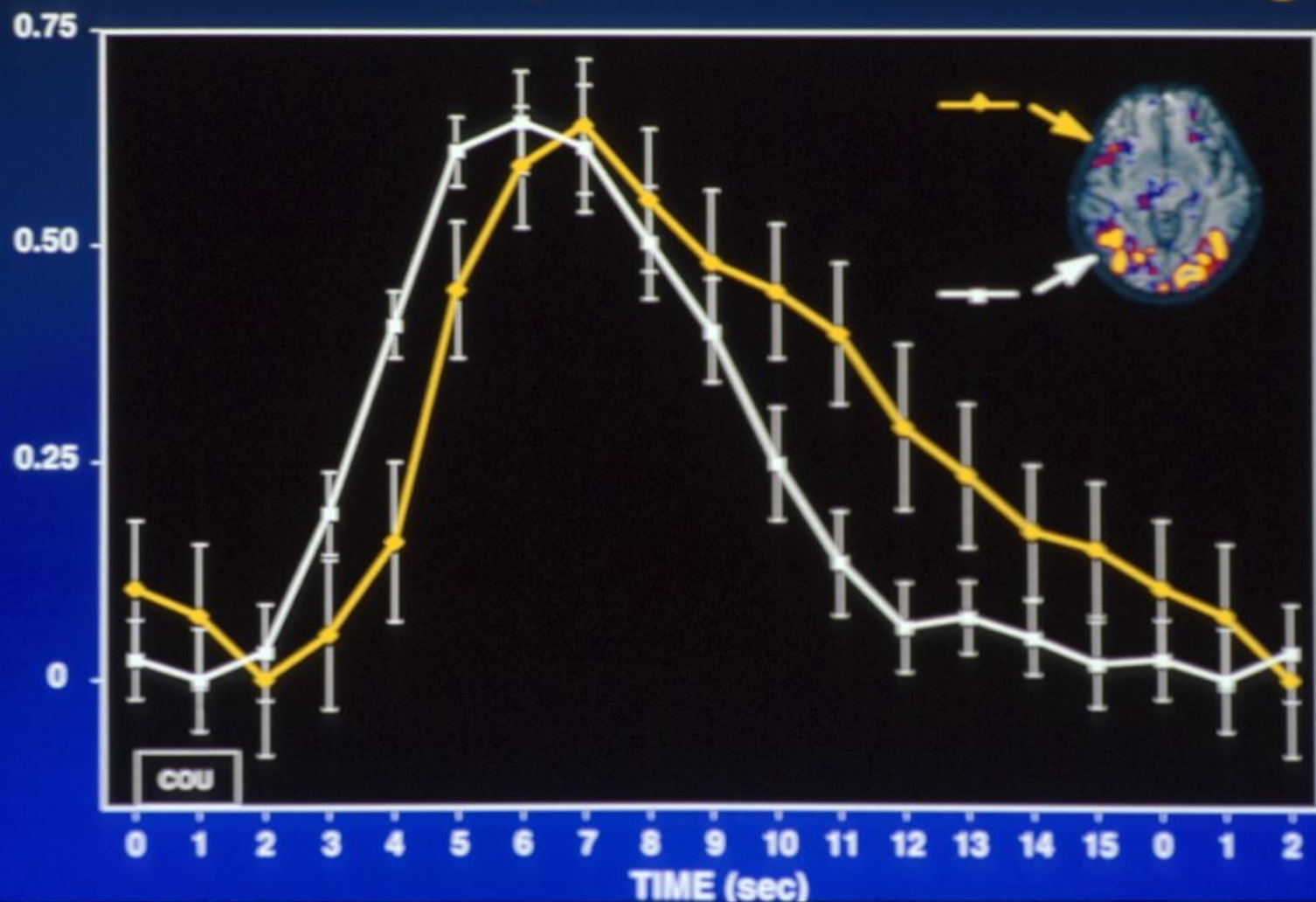


Latency

Magnitude

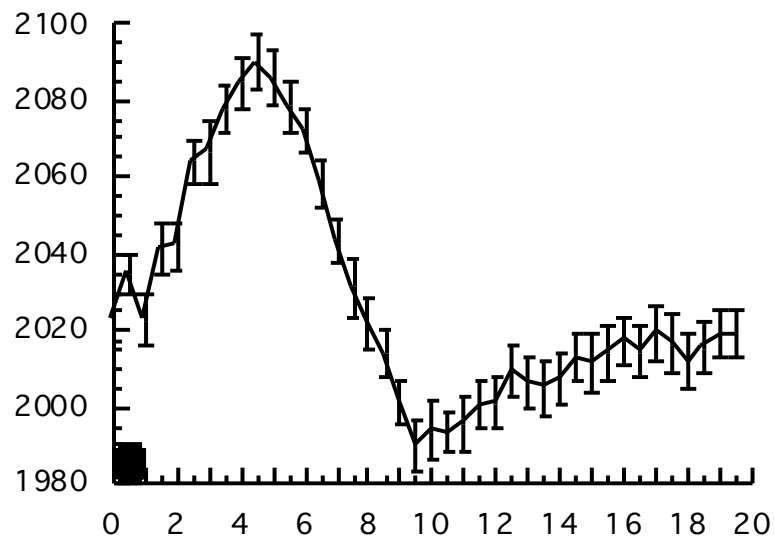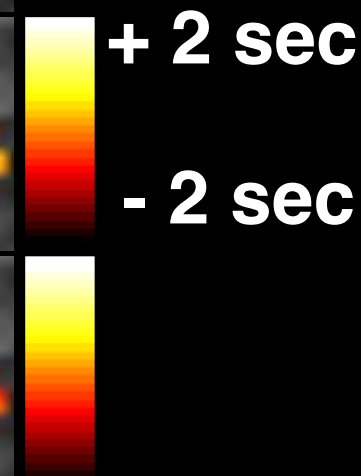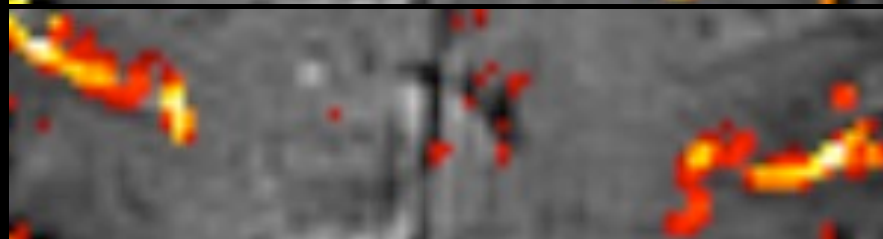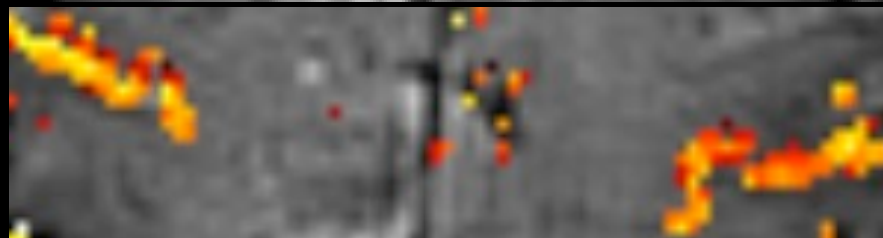


Time (sec)

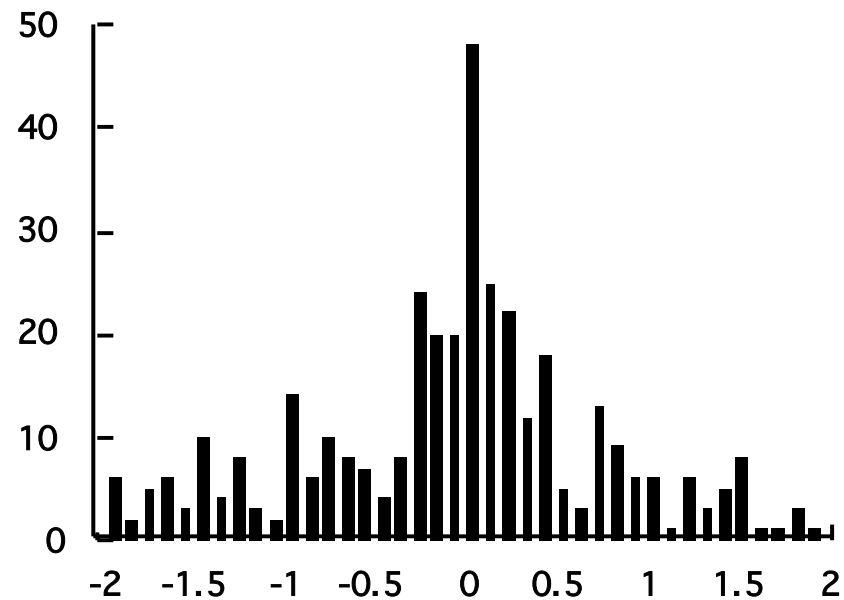


Delay (sec)

Regions of Interest Used for Hemi-Field Experiment

**Right
Hemisphere**

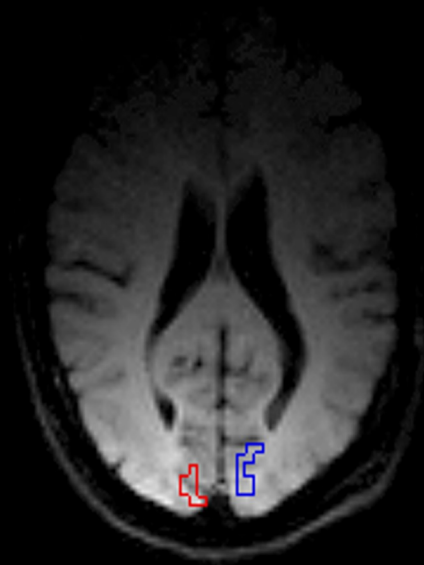



**Left
Hemisphere**

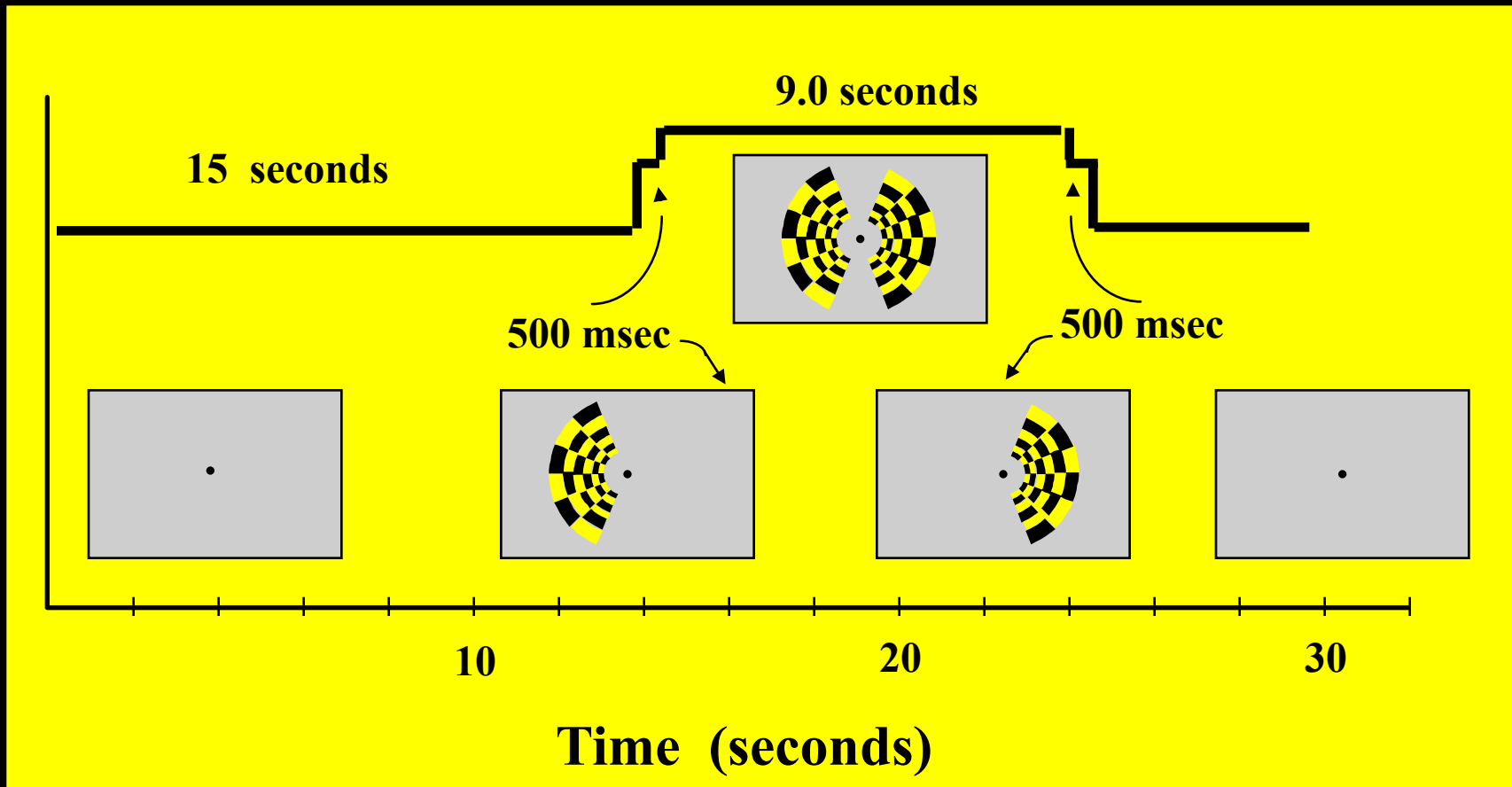

Hemi-field with 500 msec asynchrony

Average of 6 runs Standard Deviations Shown

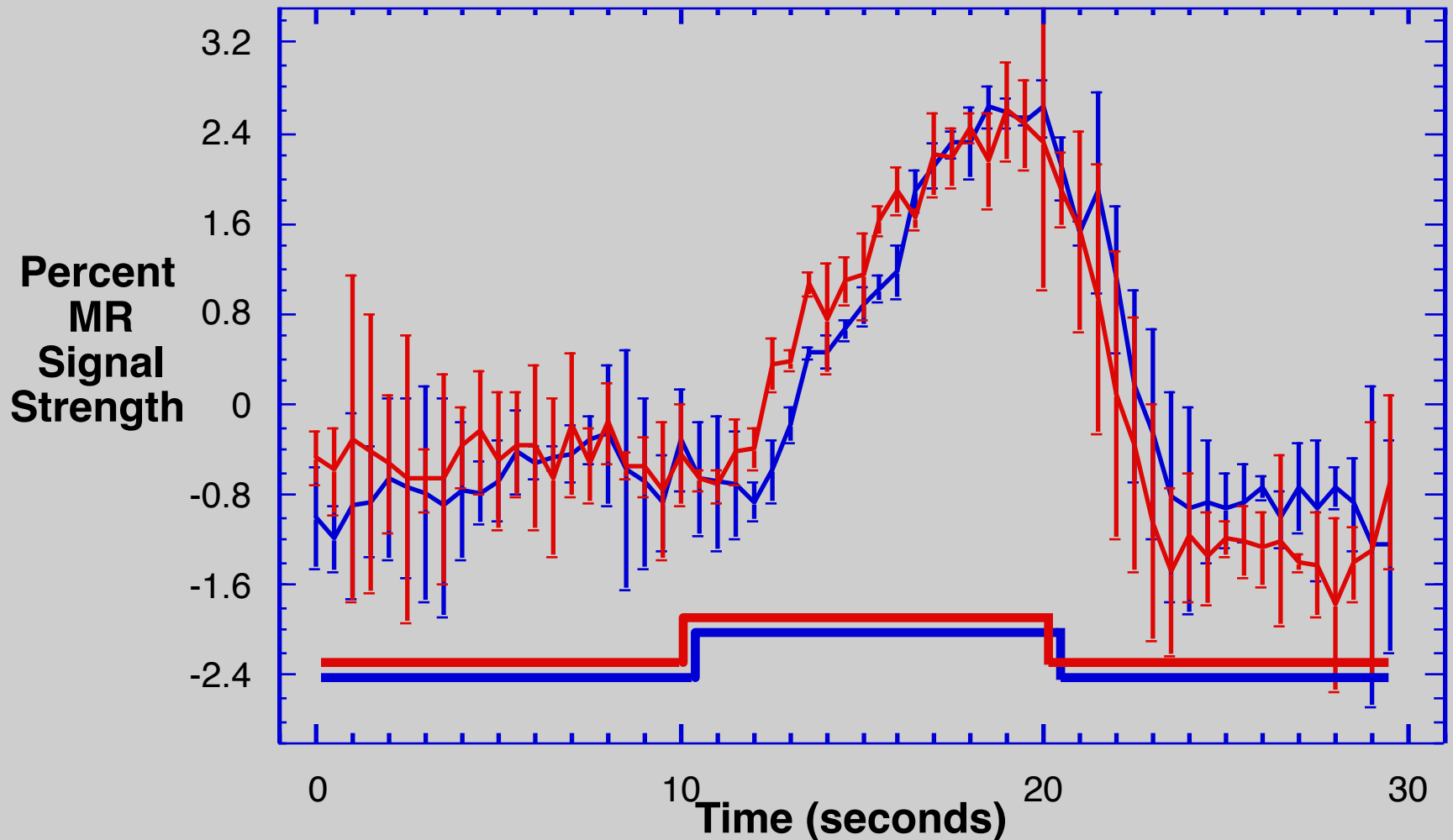

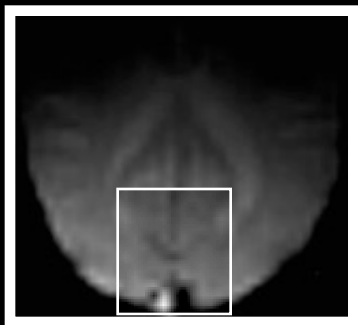

500 ms



500 ms




Right Hemifield

Left Hemifield


+ 2.5 s

0 s



- 2.5 s

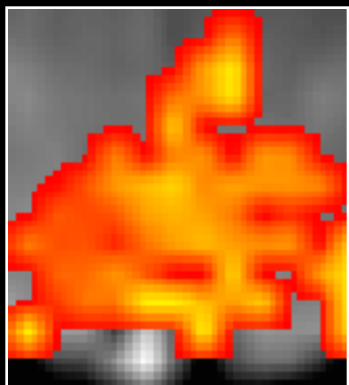


-

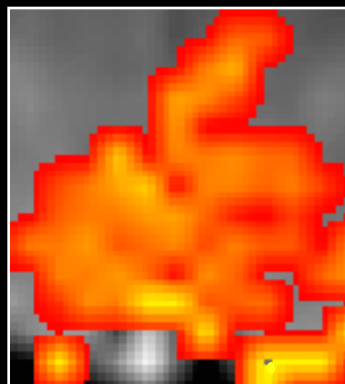

=





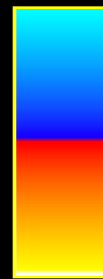
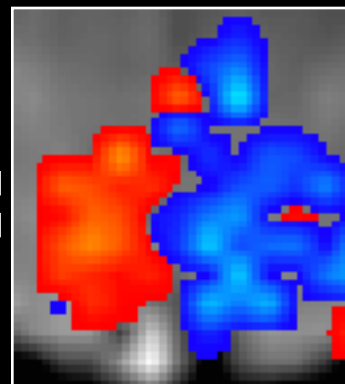


Different stimulus “ON” periods

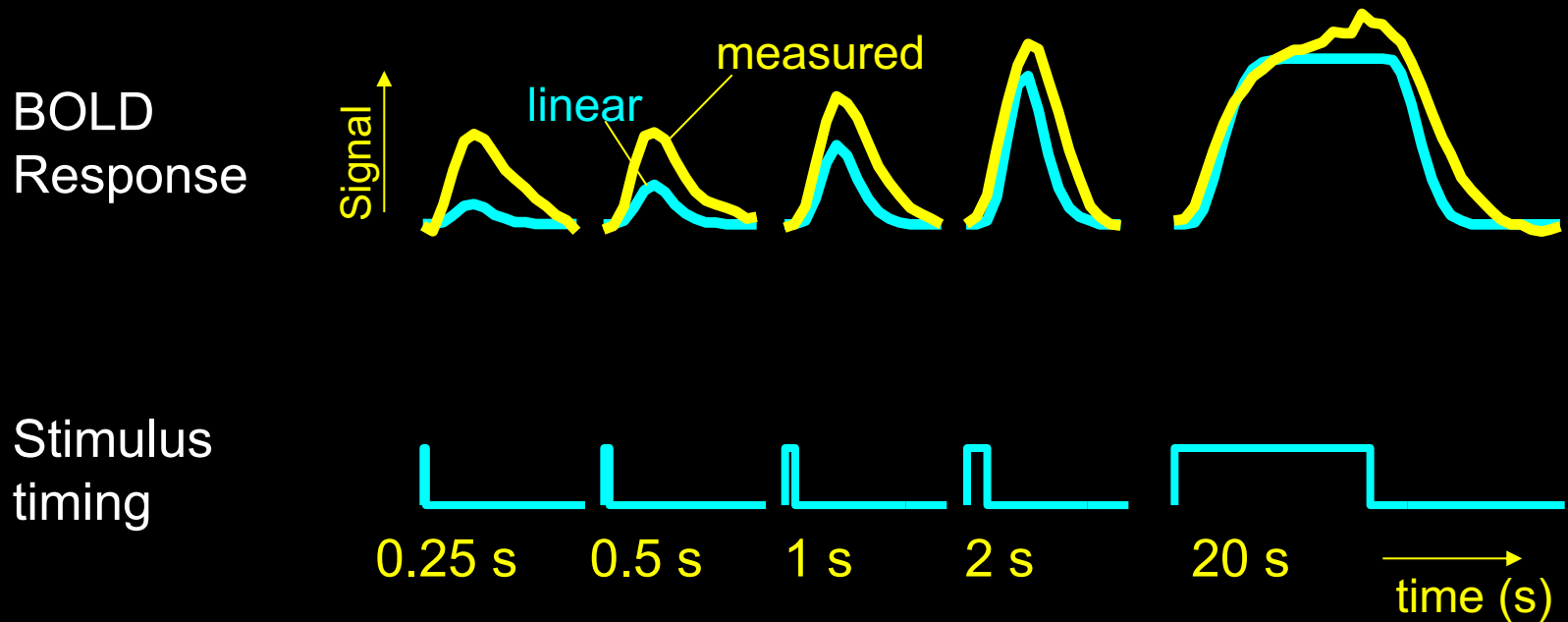


Brief stimuli produce larger responses than expected

Results — visual task

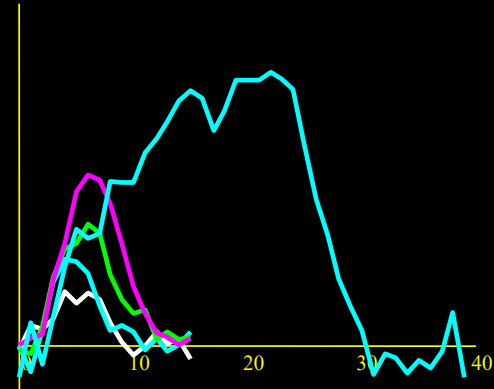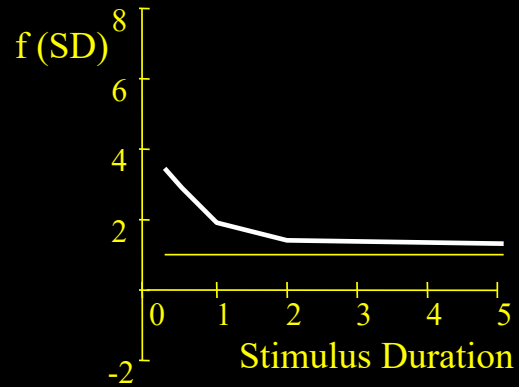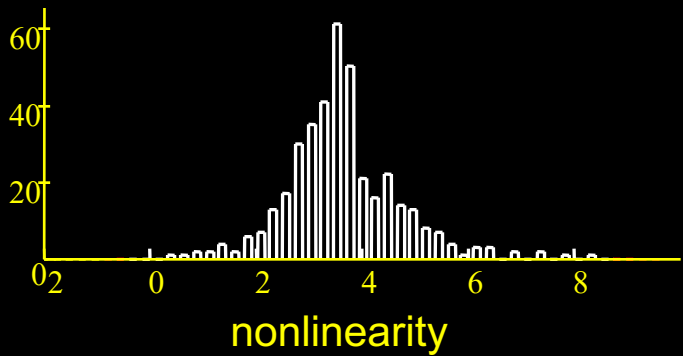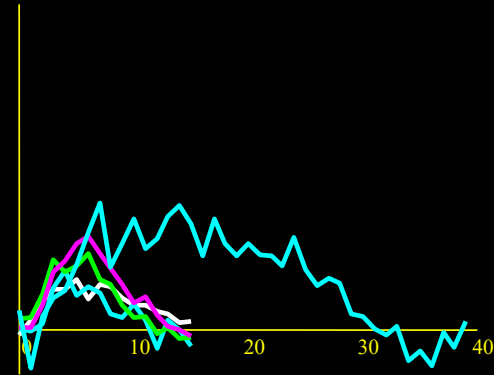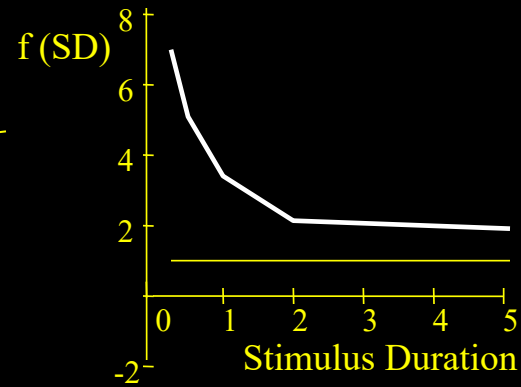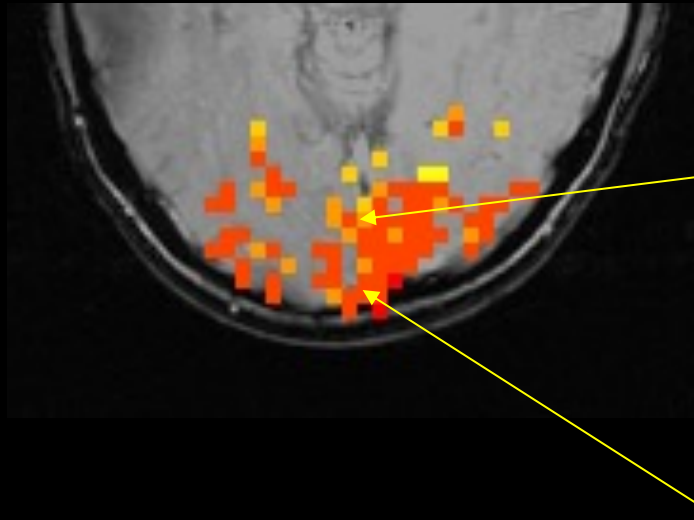


Results – visual task


Nonlinearity

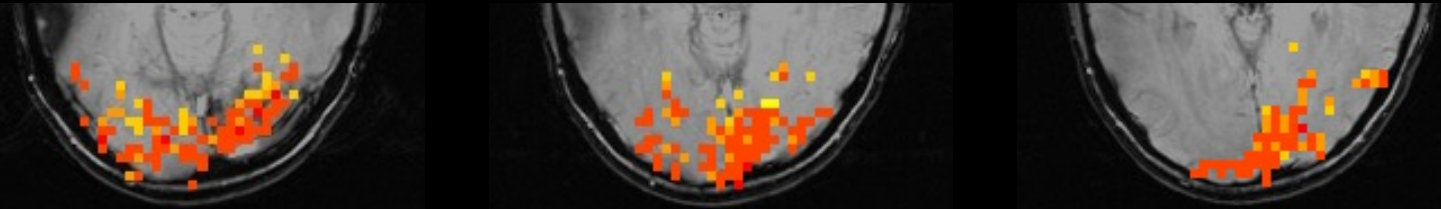







Magnitude

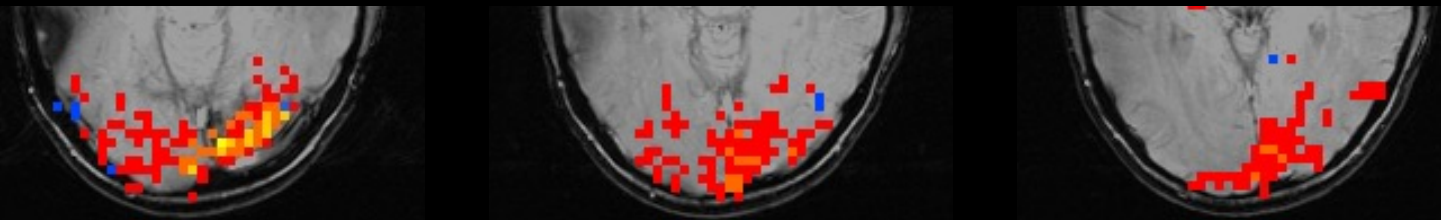

Latency

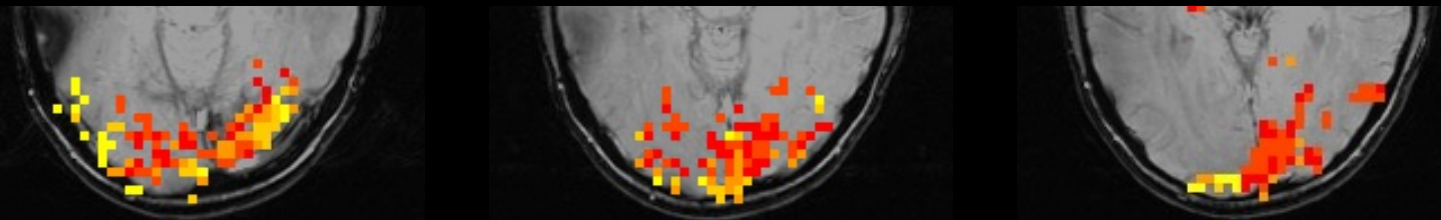


Results — motor task

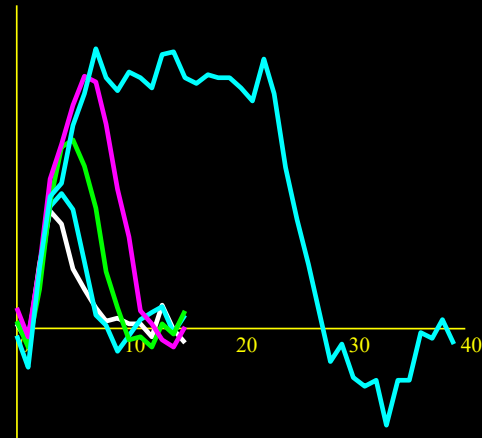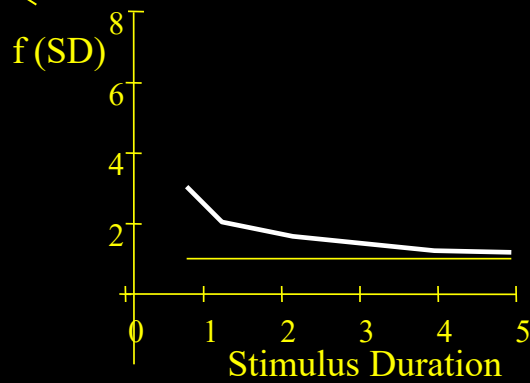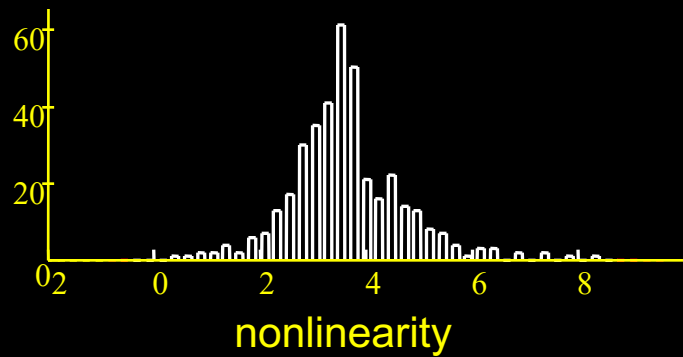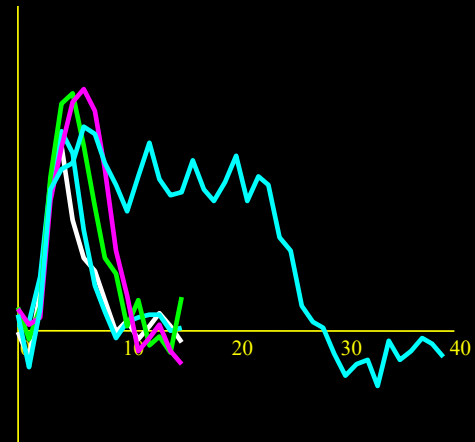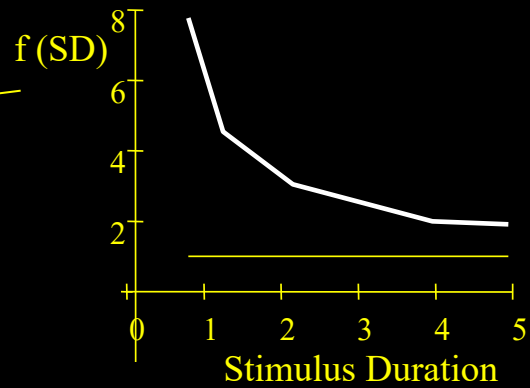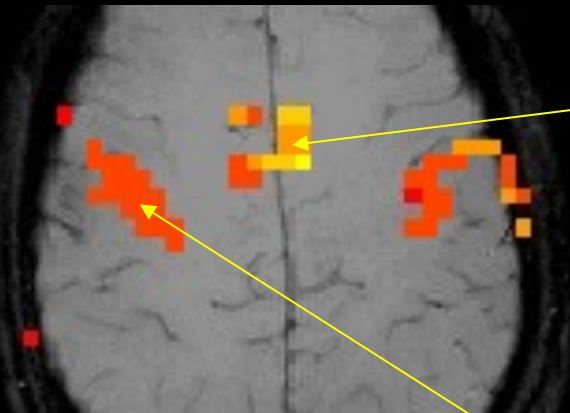


Results – motor task


Nonlinearity

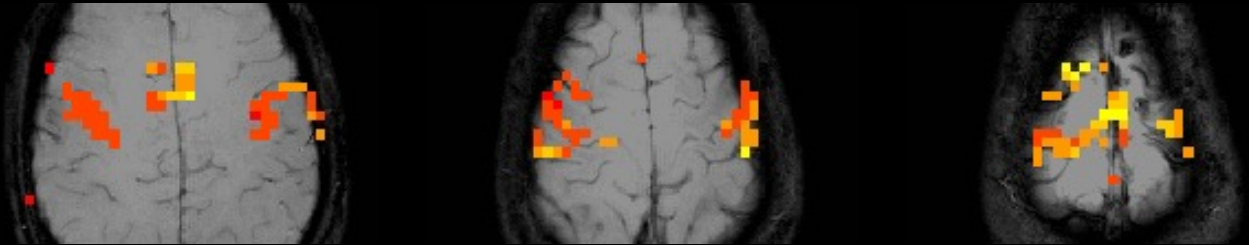


Magnitude

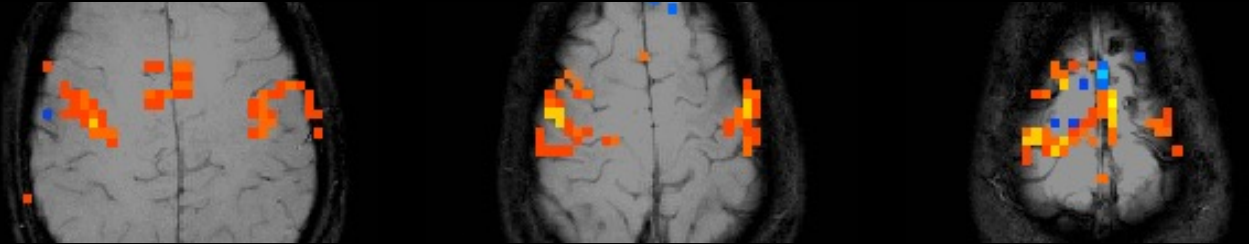

Latency

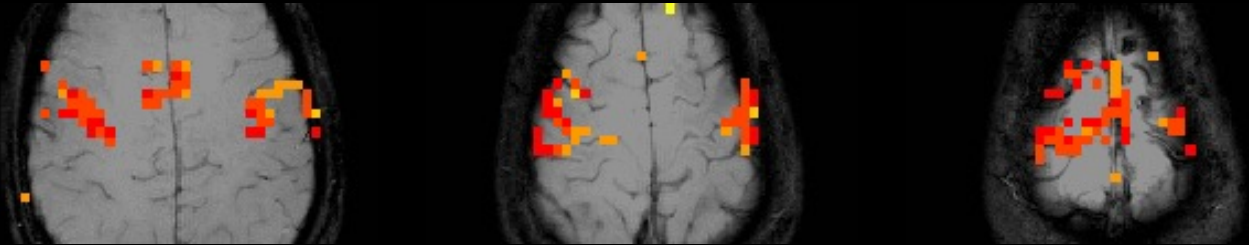



Different stimulus “ON” periods

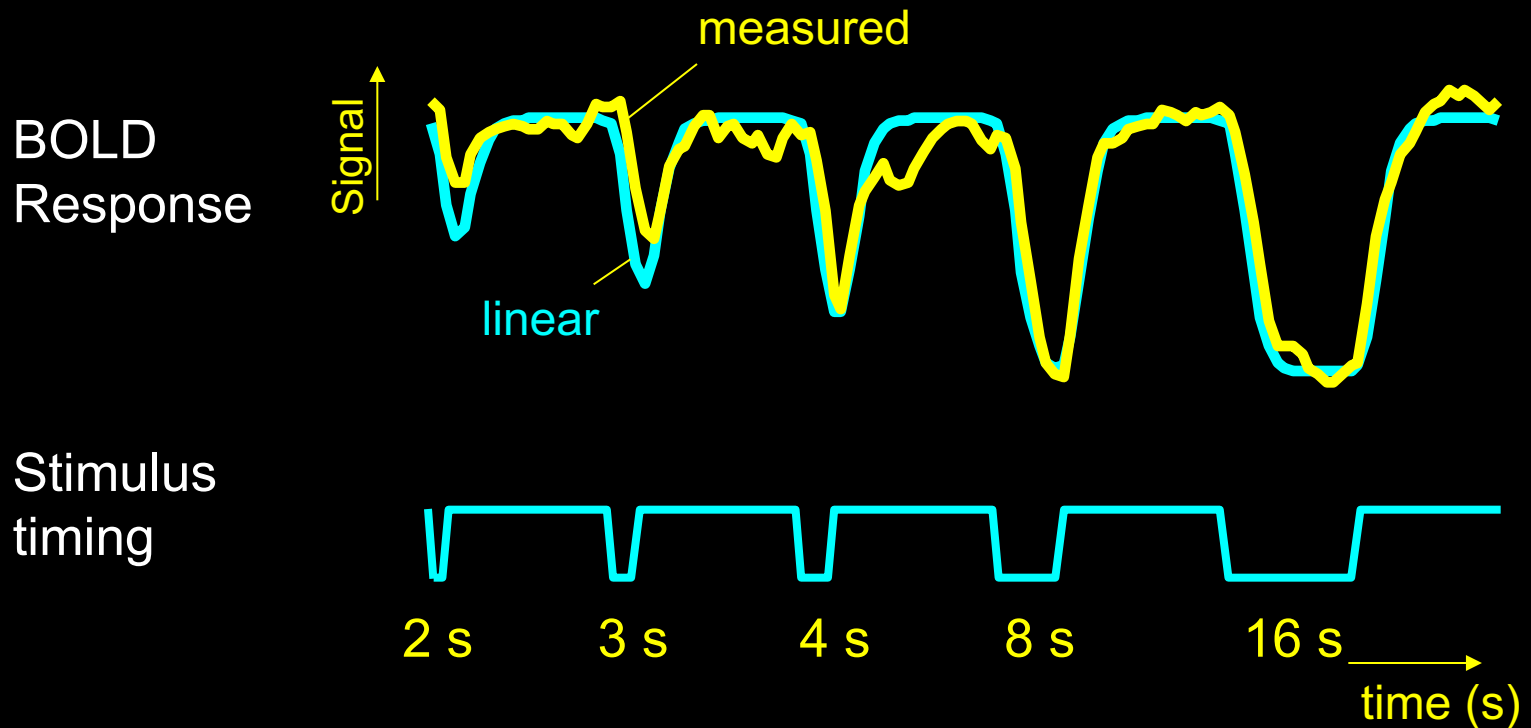

Brief stimulus OFF periods produce smaller decreases than expected


Sources of this Nonlinearity

- Neuronal

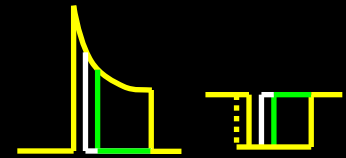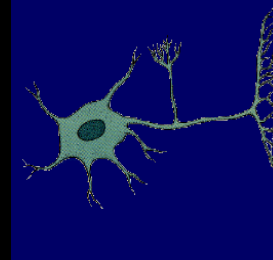


- Hemodynamic




- Oxygen extraction
- Blood volume dynamics

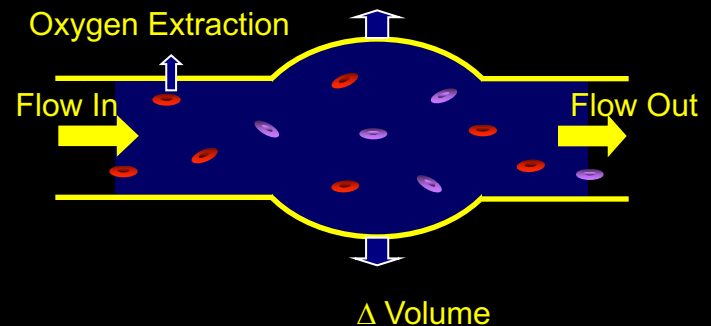


BOLD Correlation with Neuronal Activity

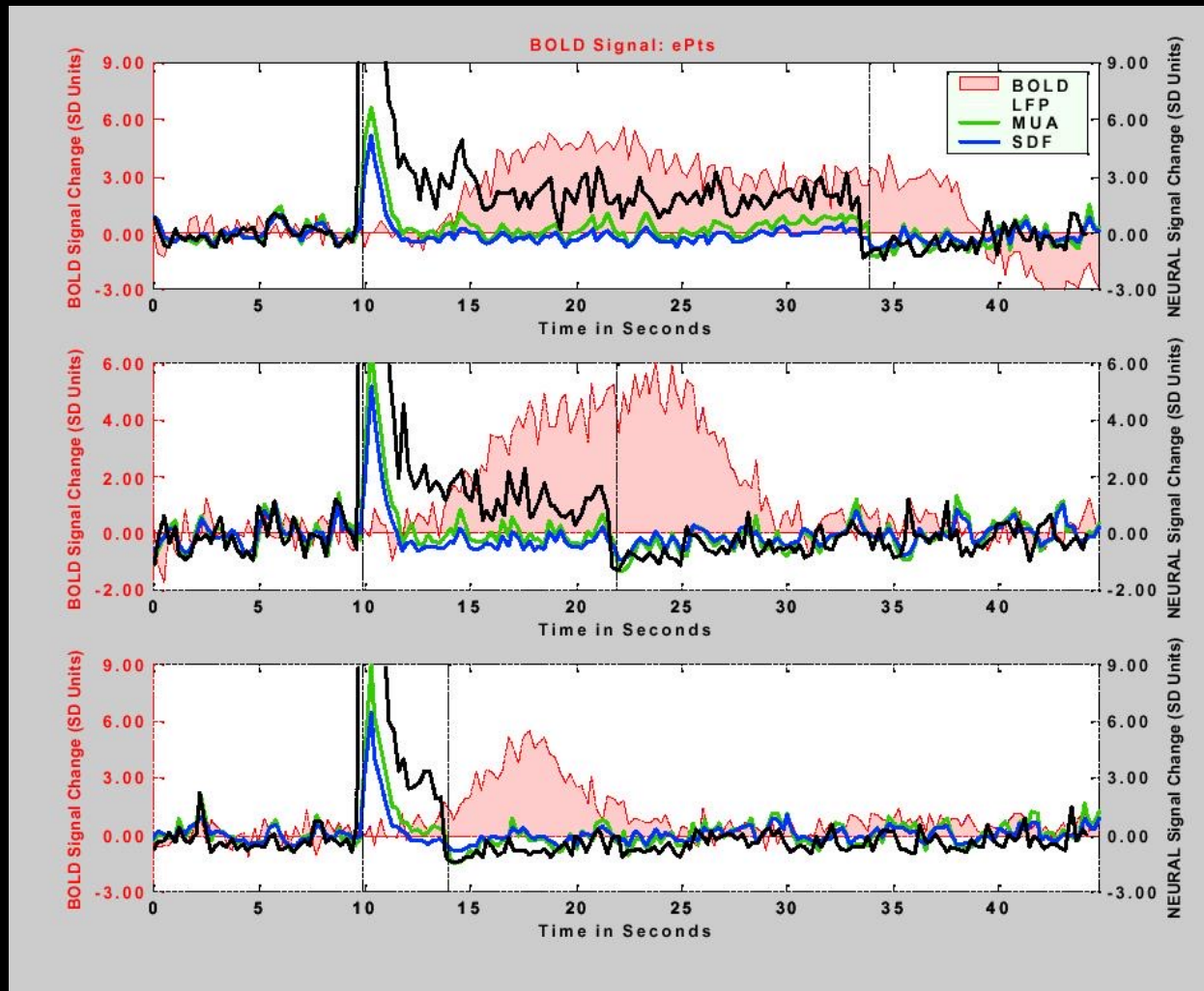

Logothetis et al. Nature, 412, 150-157

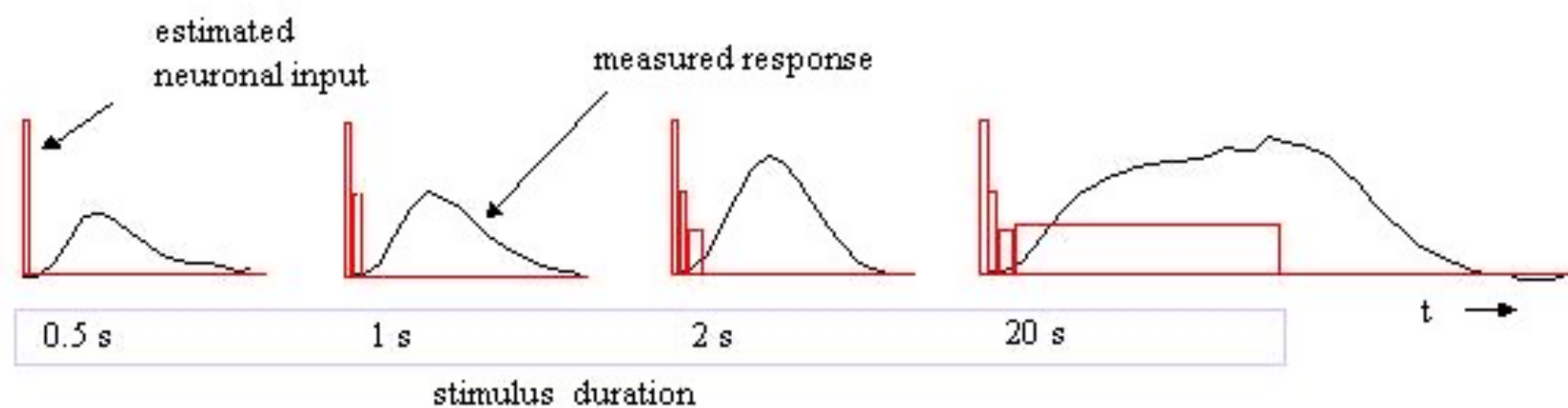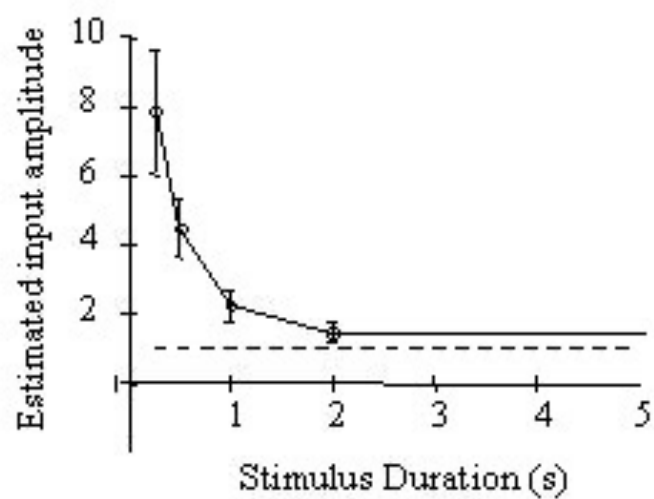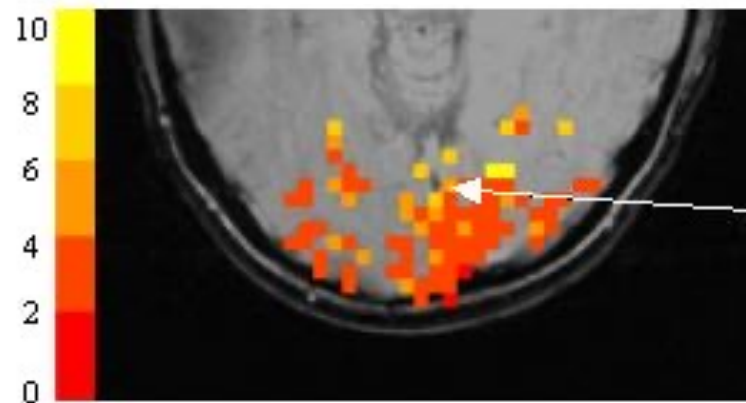

Functional Contrast

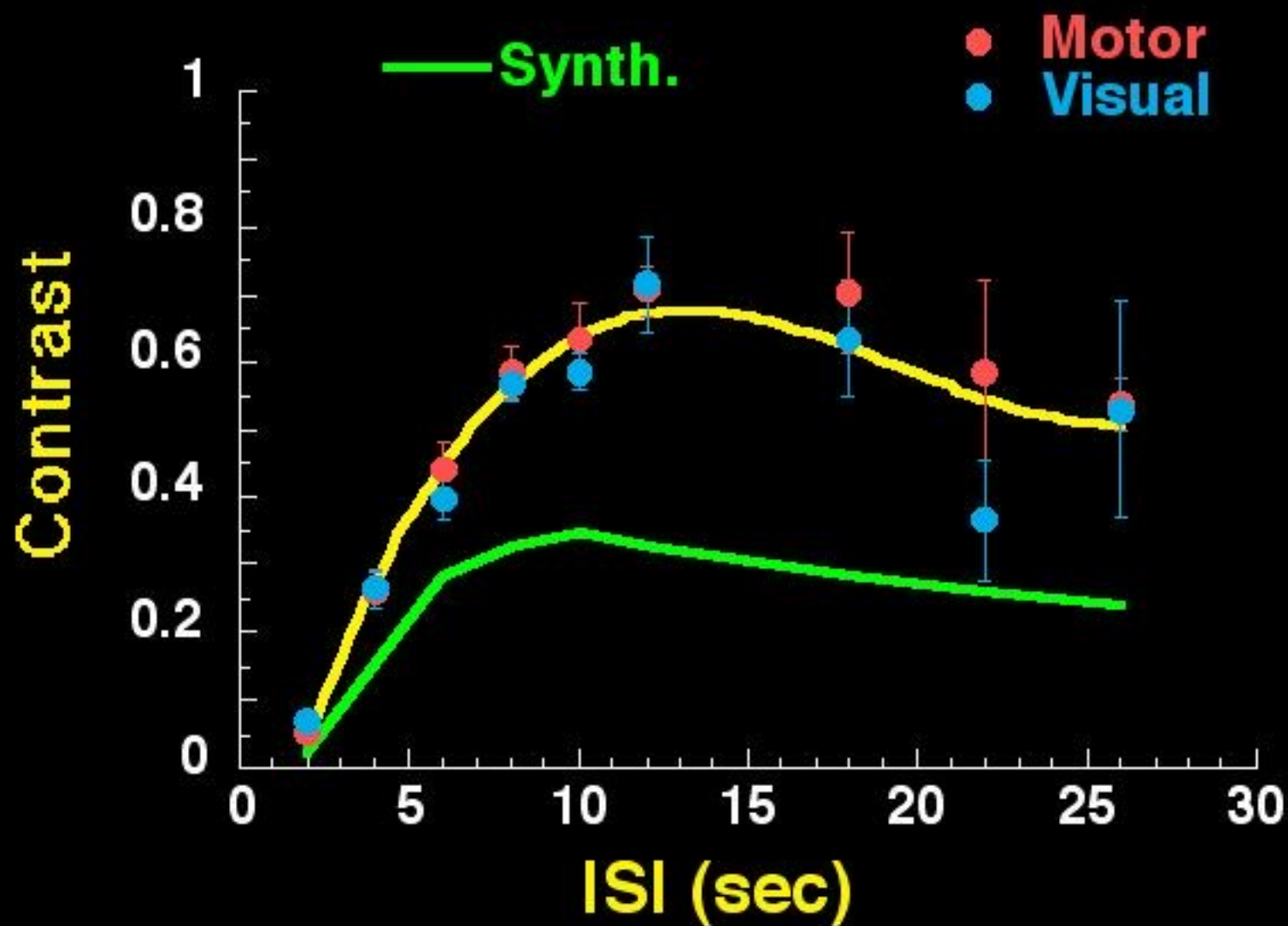

(Block design = 1)

Contrast to Noise Images

(ISI, SD)

20, 20

12, 2

10, 2


8, 2

6, 2


4, 2

2, 2

S1

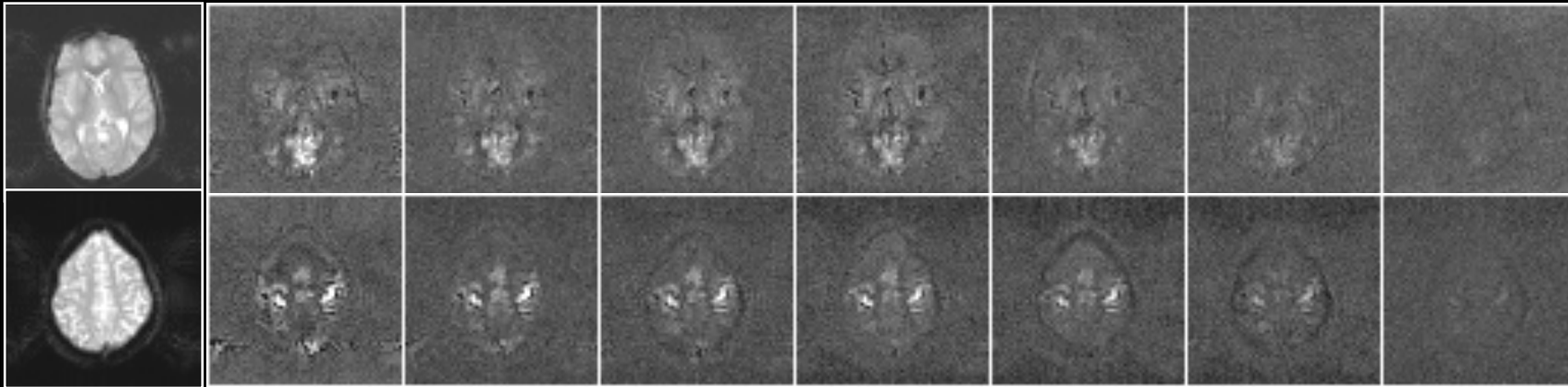



S2

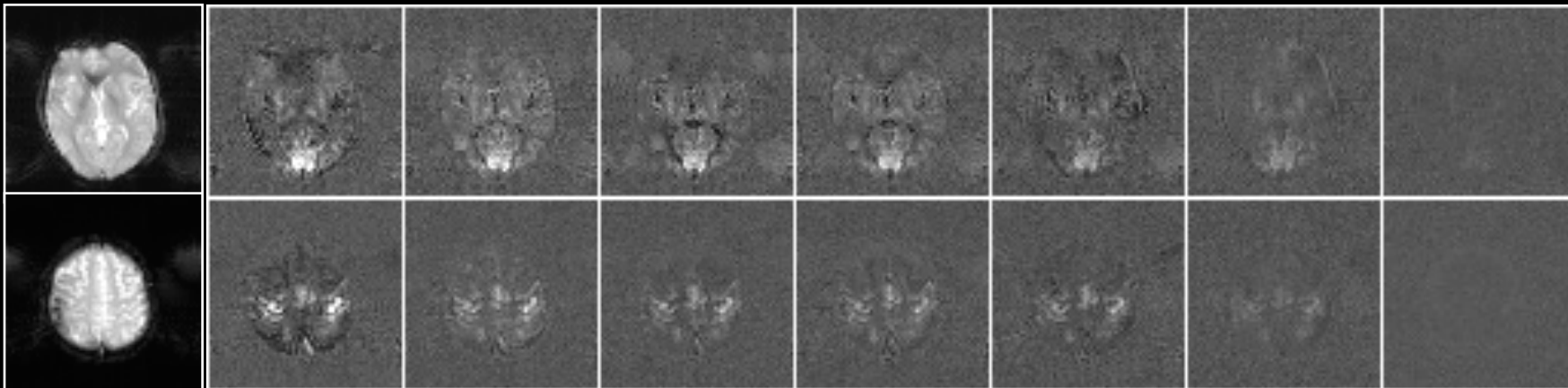


Detectability – constant ISI

SD – stimulus duration

ISI – inter-stimulus interval

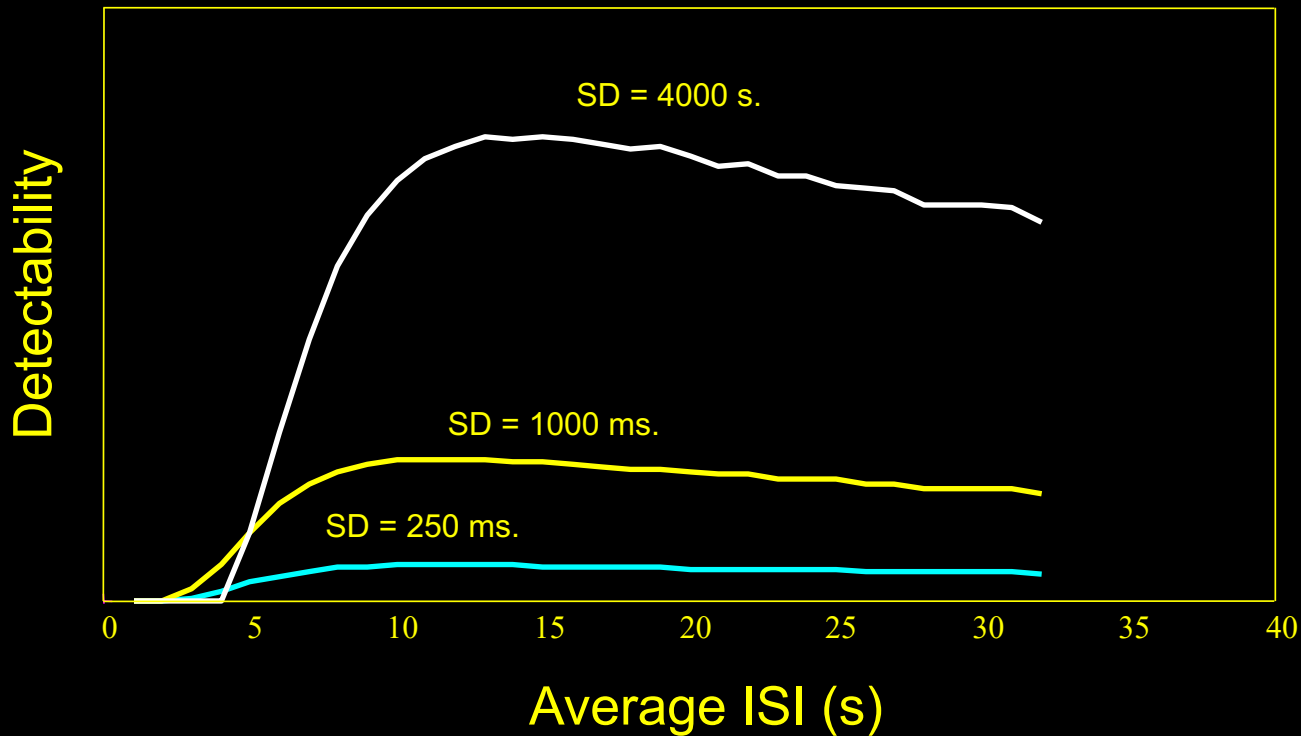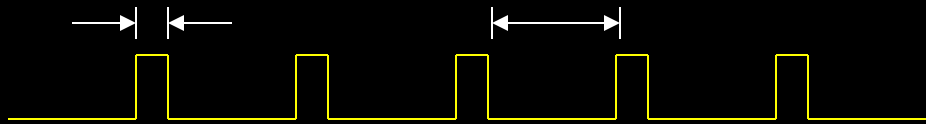


Visual Activation Paradigm: 1 , 2, & 3 Trials

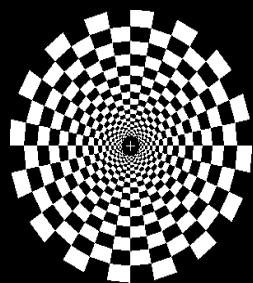


0 sec

20 sec

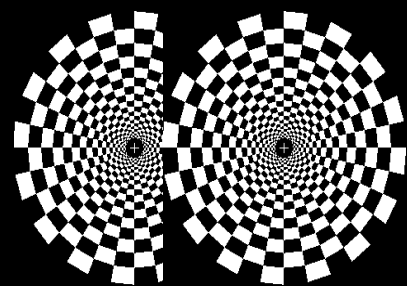

0 sec 2 sec

20 sec

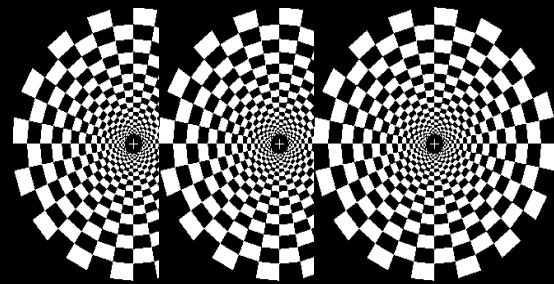


0 sec 2 sec 4 sec


20 sec

Response to Multiple Trials: Subject RW


RAW DATA

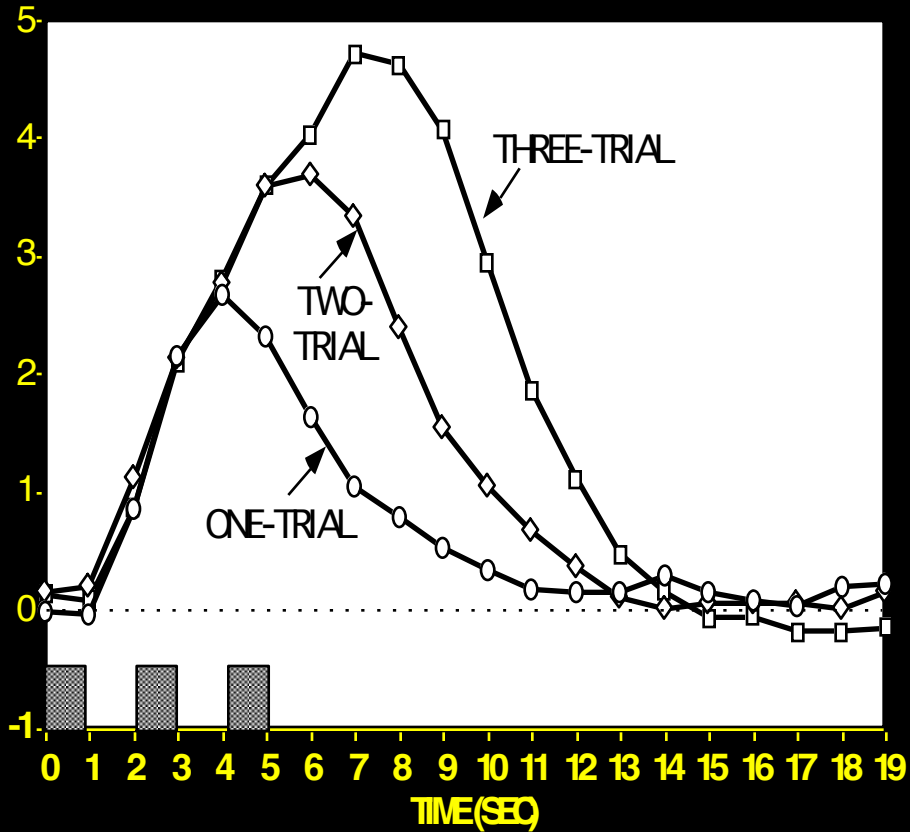


ESTIMATED RESPONSES

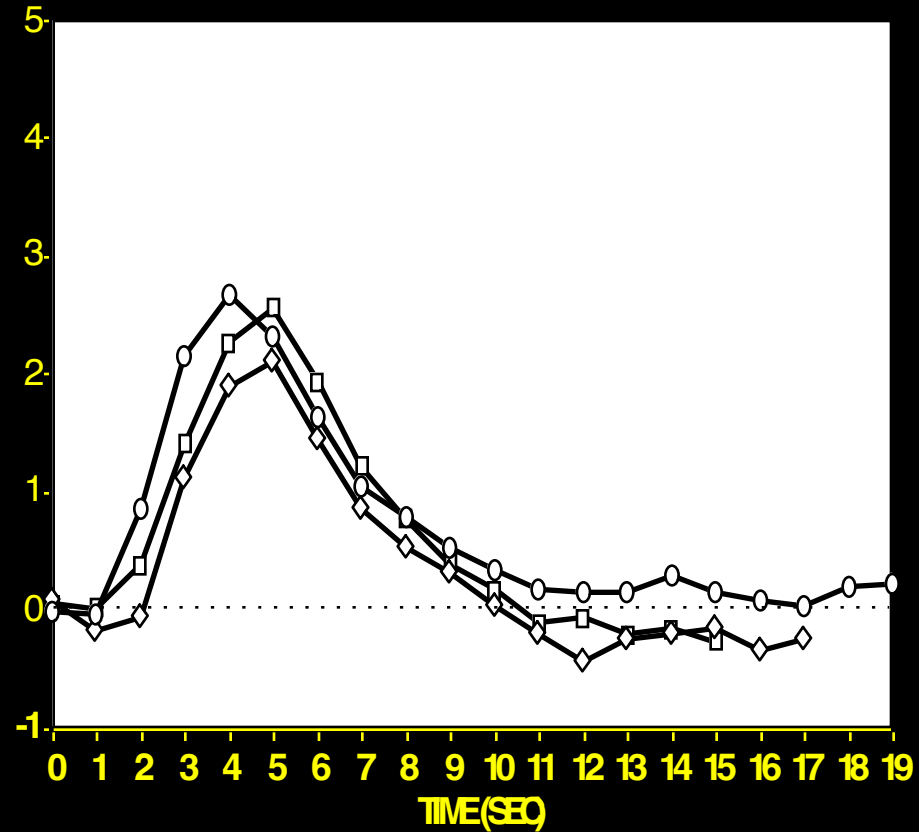

Detectability vs. Average ISI

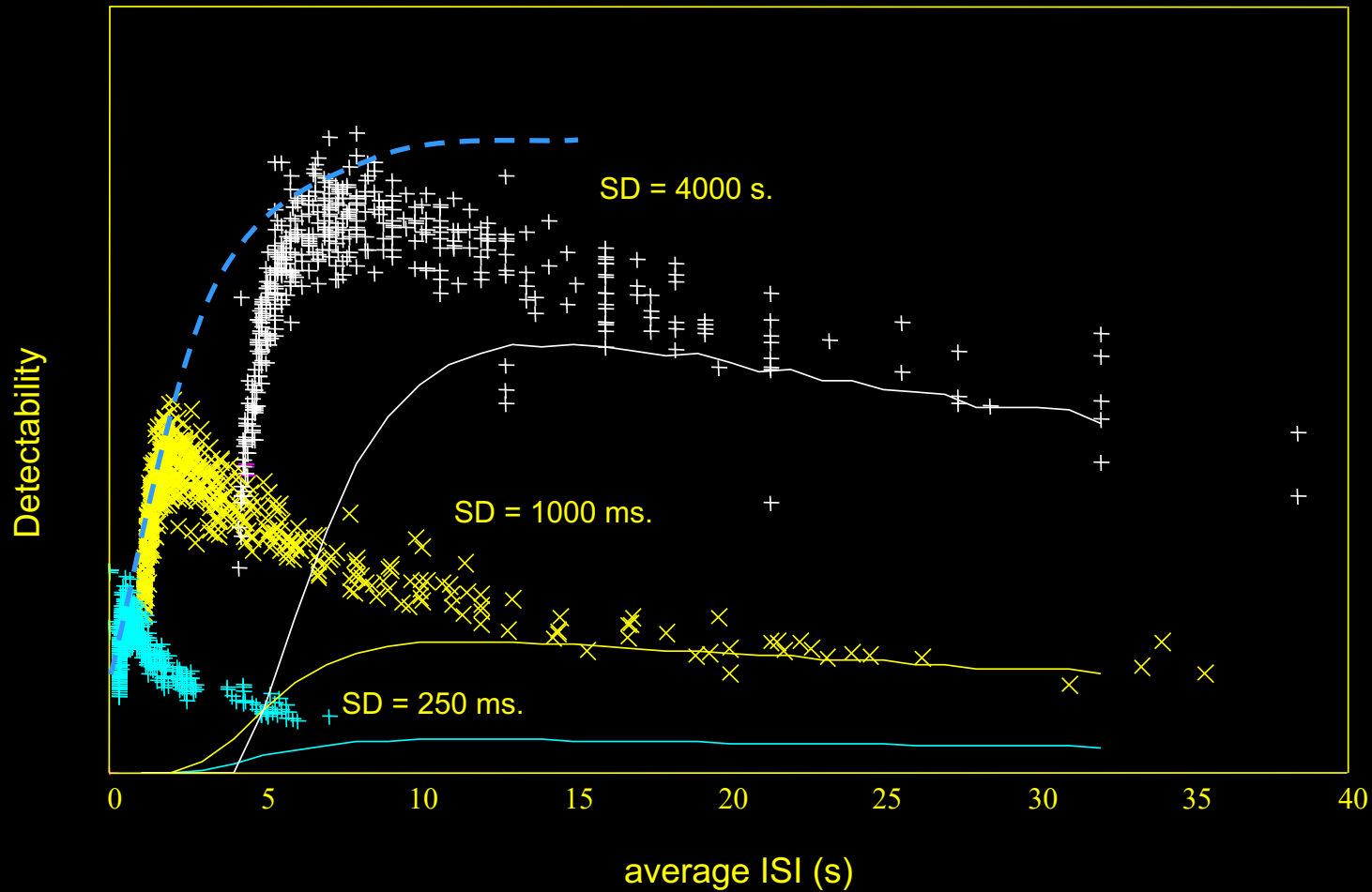

Estimation accuracy vs. average ISI

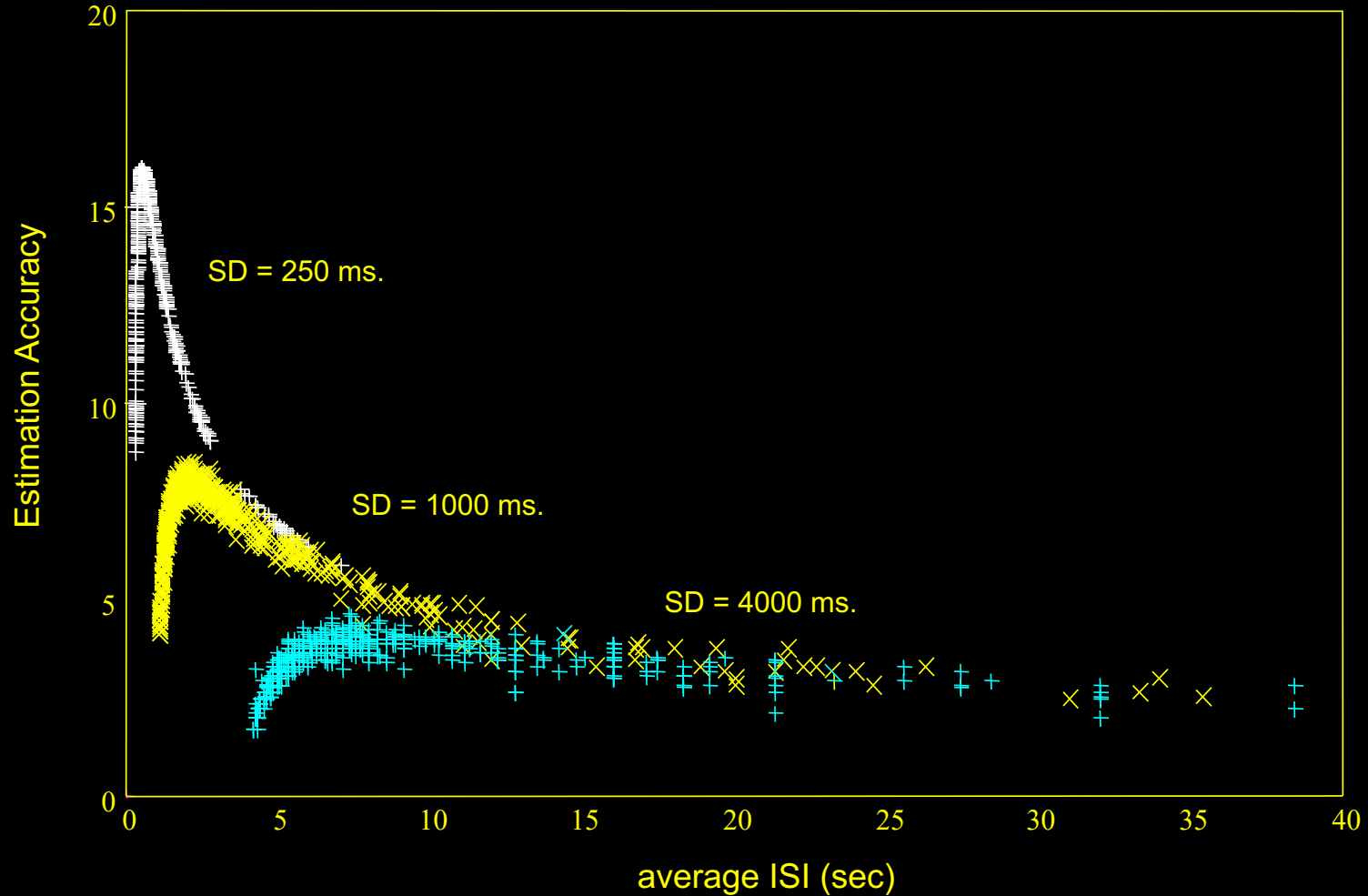

Motion

Recognize?

- Edge effects
- Shorter signal change latencies
- Unusually high signal changes
- External measuring devices

Correct?

- Image registration algorithms
- Orthogonalize to motion-related function (*cardiac, respiration, movement*)
- Navigator echo for k-space alignment
(*for multishot techniques*)
- Re-do scan


Bypass?

- Paradigm timing strategies..
- Gating (with T1-correction)

Suppress?


- Flatten image contrast
- Physical restraint
- Averaging, smoothing

Speaking - Blocked Trial

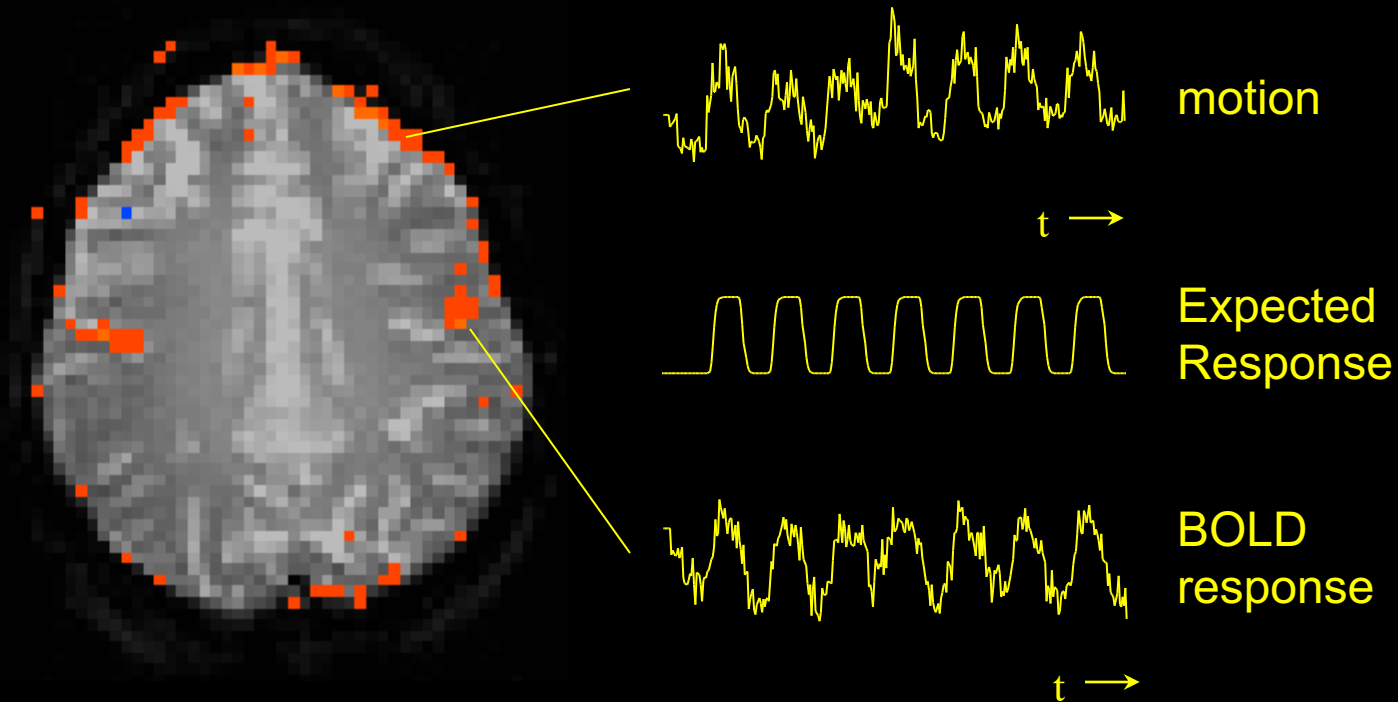


fMRI during tasks that involve brief motion


Blocked Design

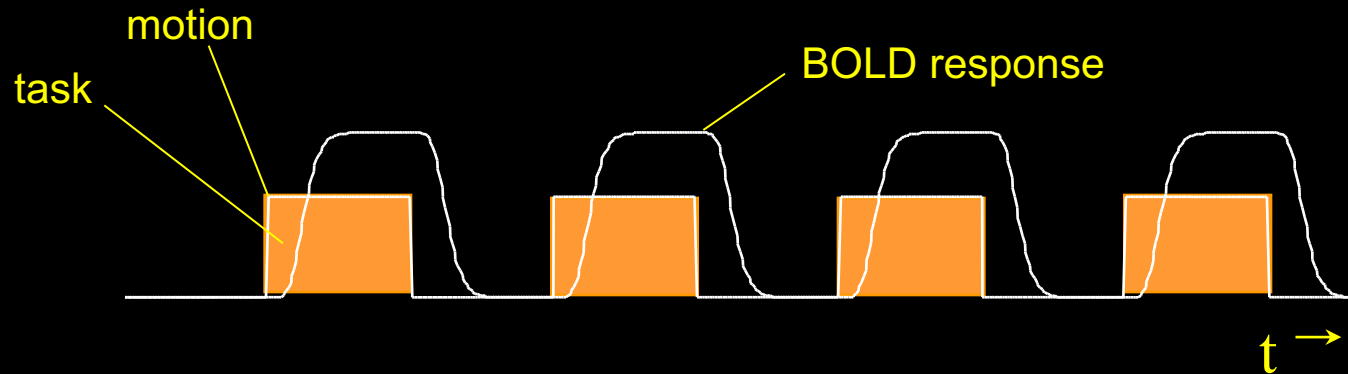

Event-Related Design

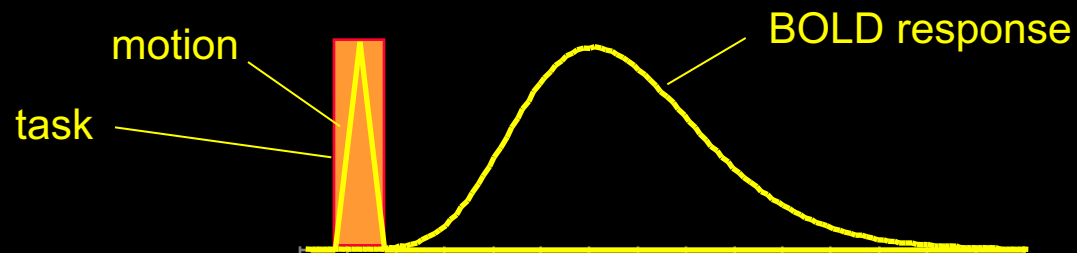

Overt Word Production

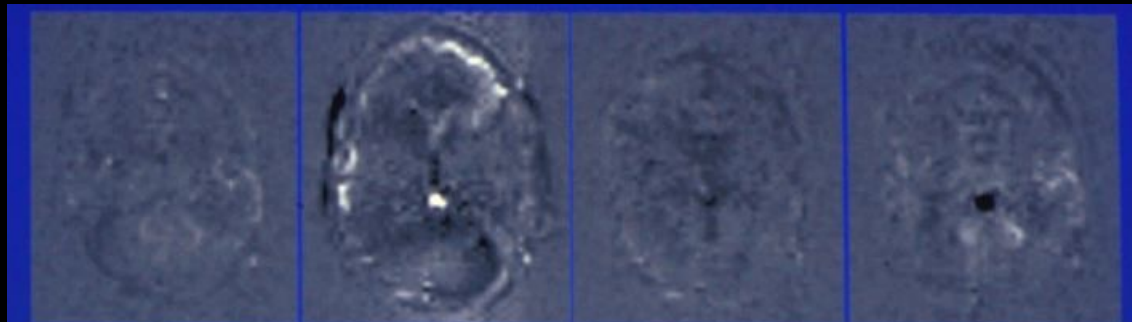

2

3

4

5



6

7

8

9




10

11

12


13

Motion-Decoupled fMRI: Functional MRI during of overt word production

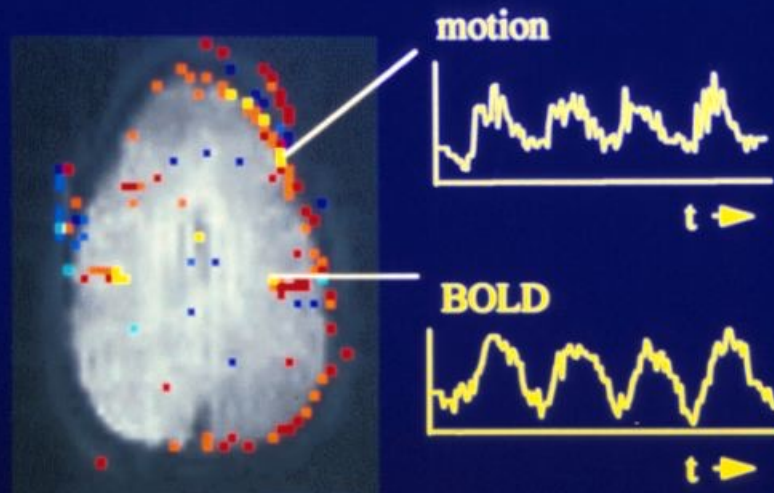

“block-trial” paradigm

Motion induced signal changes resemble functional (BOLD) signal changes

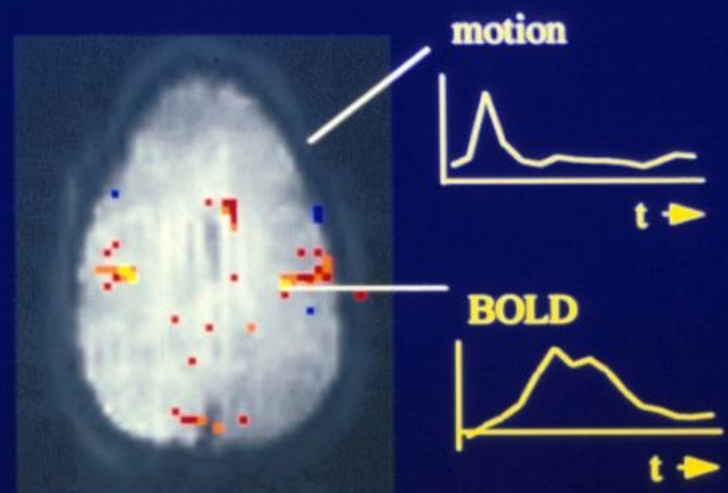


“single-trial” paradigm


Motion induced and BOLD signal changes are separated in time

R.M. Birn, et al.

Tongue Movement

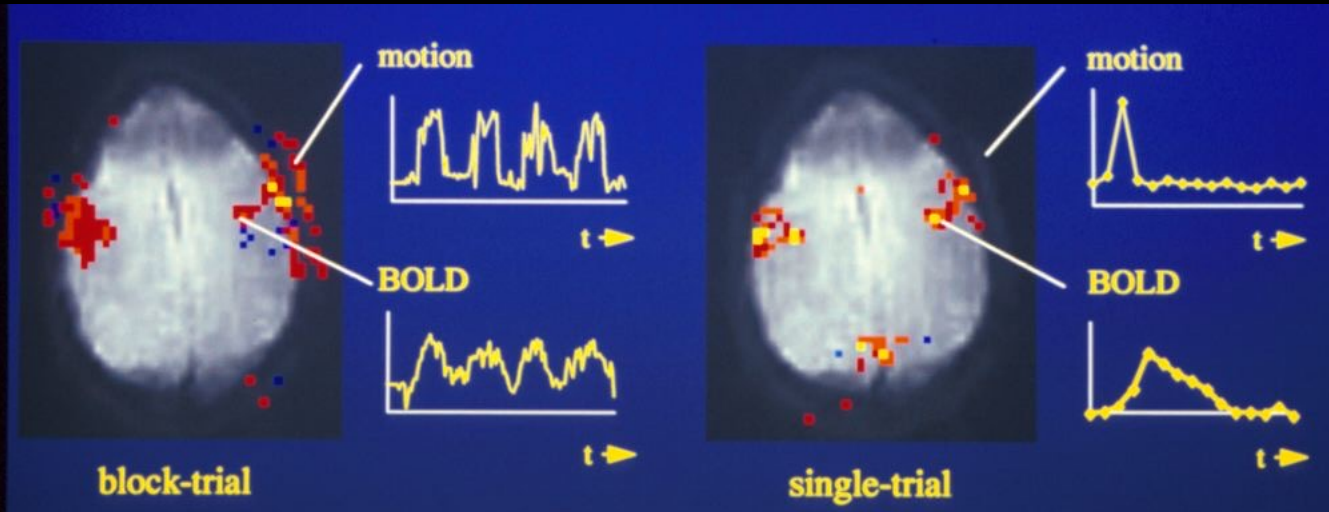


Jaw Clenching

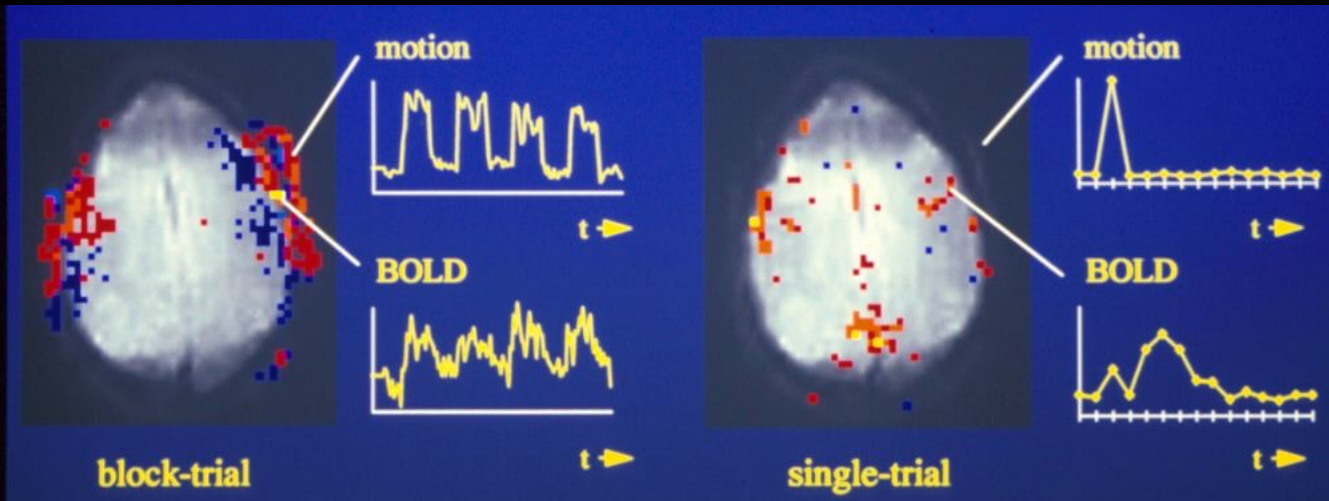


Constant ISI


Speaking - ER-fMRI

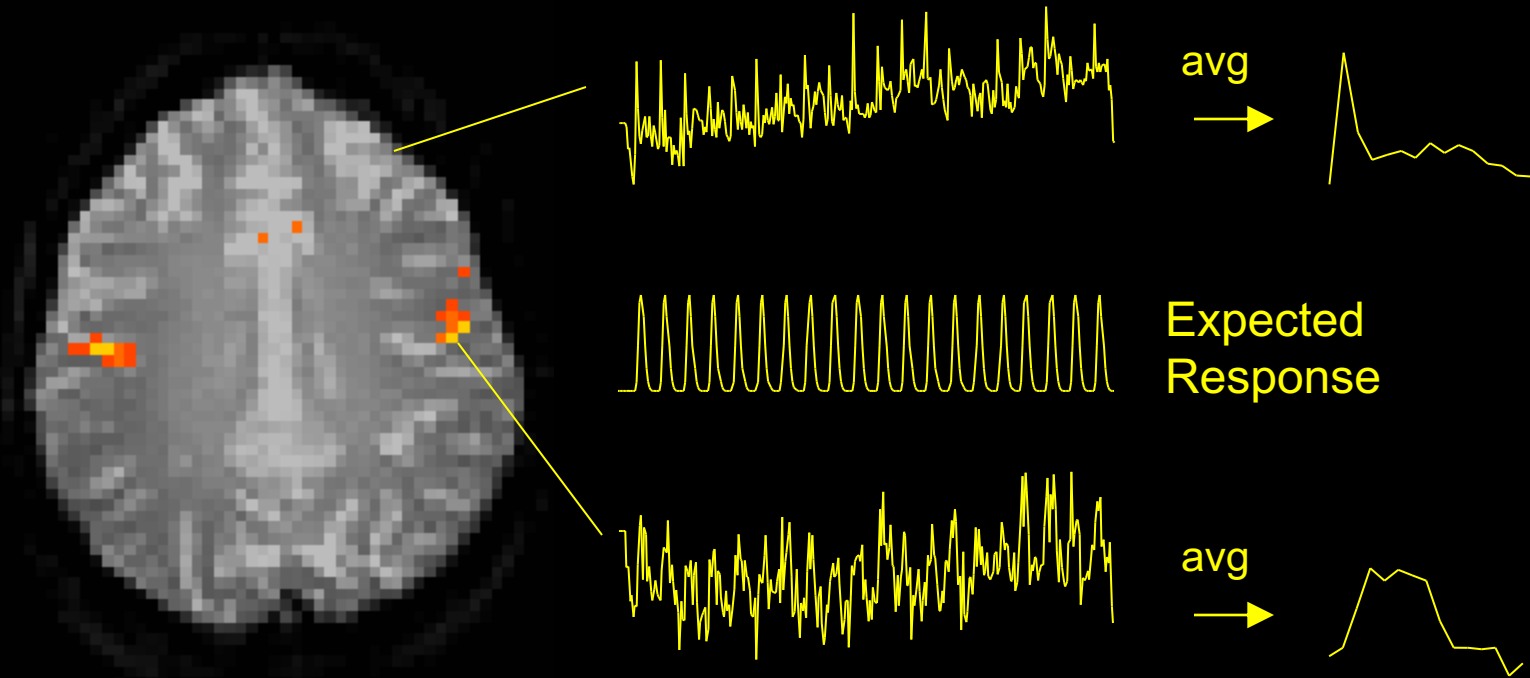

Swallowing - Event-Related

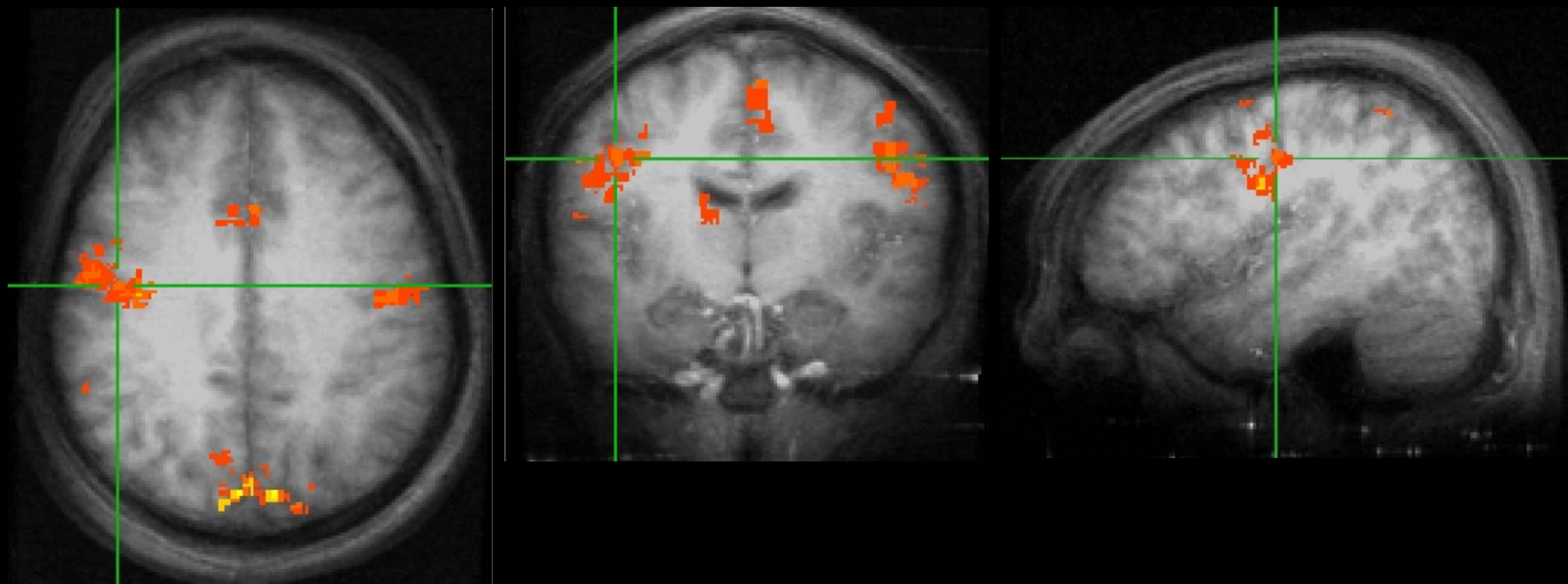

ODC Maps using fMRI

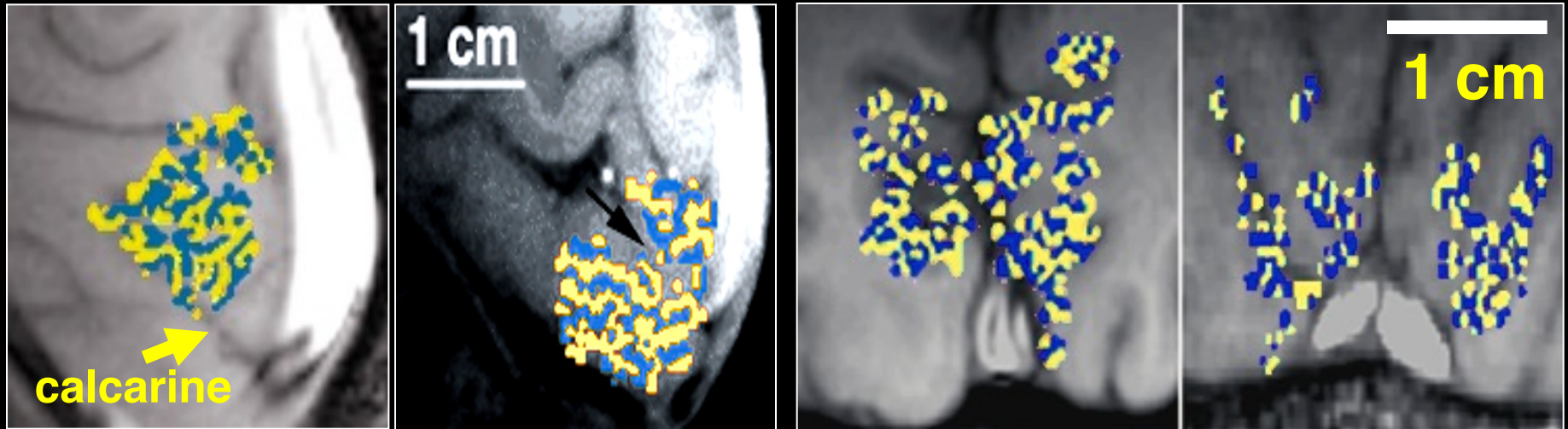

- Identical in size, orientation, and appearance to those obtained by optical imaging¹ and histology^{3,4}.


¹Malonek D, Grinvald A. *Science* 272, 551-4 (1996).

³Horton JC, Hocking DR. *J Neurosci* 16, 7228-39 (1996).

⁴Horton JC, et al. *Arch Ophthalmol* 108, 1025-31 (1990).

Why short is better than long


The vascular response to prolonged neural stimulation

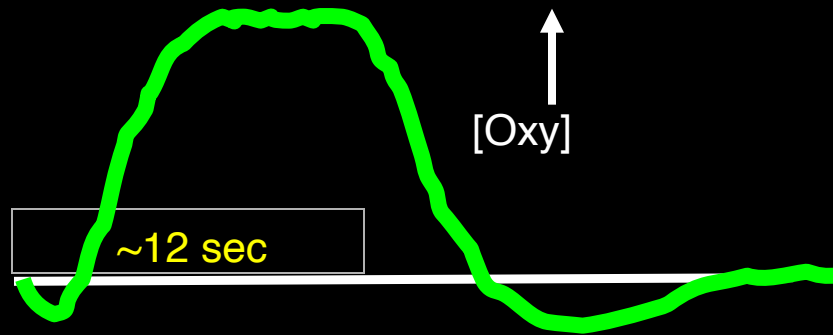

It is argued that fMRI cannot achieve submillimeter functional resolution because a saturated hyperoxic vascular response to neural activity spreads over many millimeters^{1,2}.

However, optical imaging has demonstrated that the hyperoxic response can yield well-localized maps when using short duration stimuli (<5 sec)¹.

The vascular response to brief neural stimulation

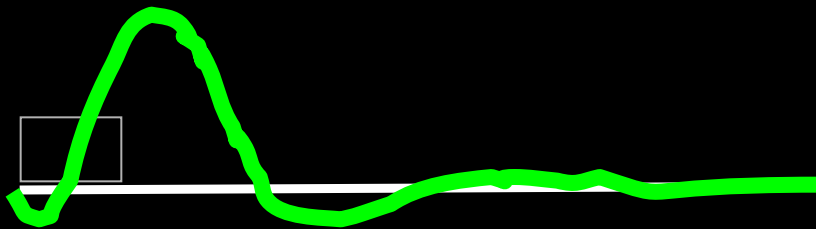


¹Malonek D, Grinvald A. Science 272, 551-4 (1996).

²Kim D-S, Duong T, Kim S-G. Nat Neurosci 3, 164-9 (2000).


Neuronal Activation Input Strategies

1. Block Design

2. Frequency Encoding

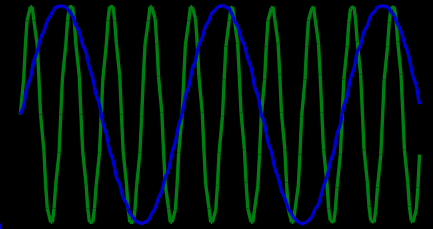


3. Phase Encoding

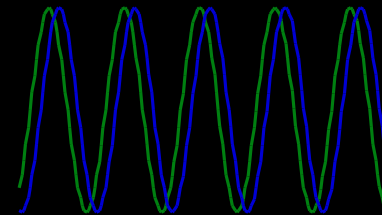


4. Single Event

5. Orthogonal Block Design

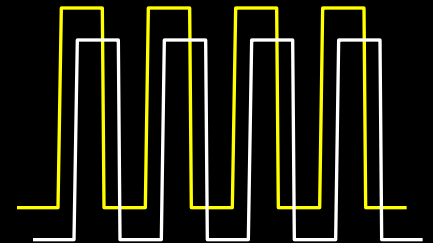

6. Free Behavior Design.


Example of a Set of Orthogonal Contrasts for Multiple Regression

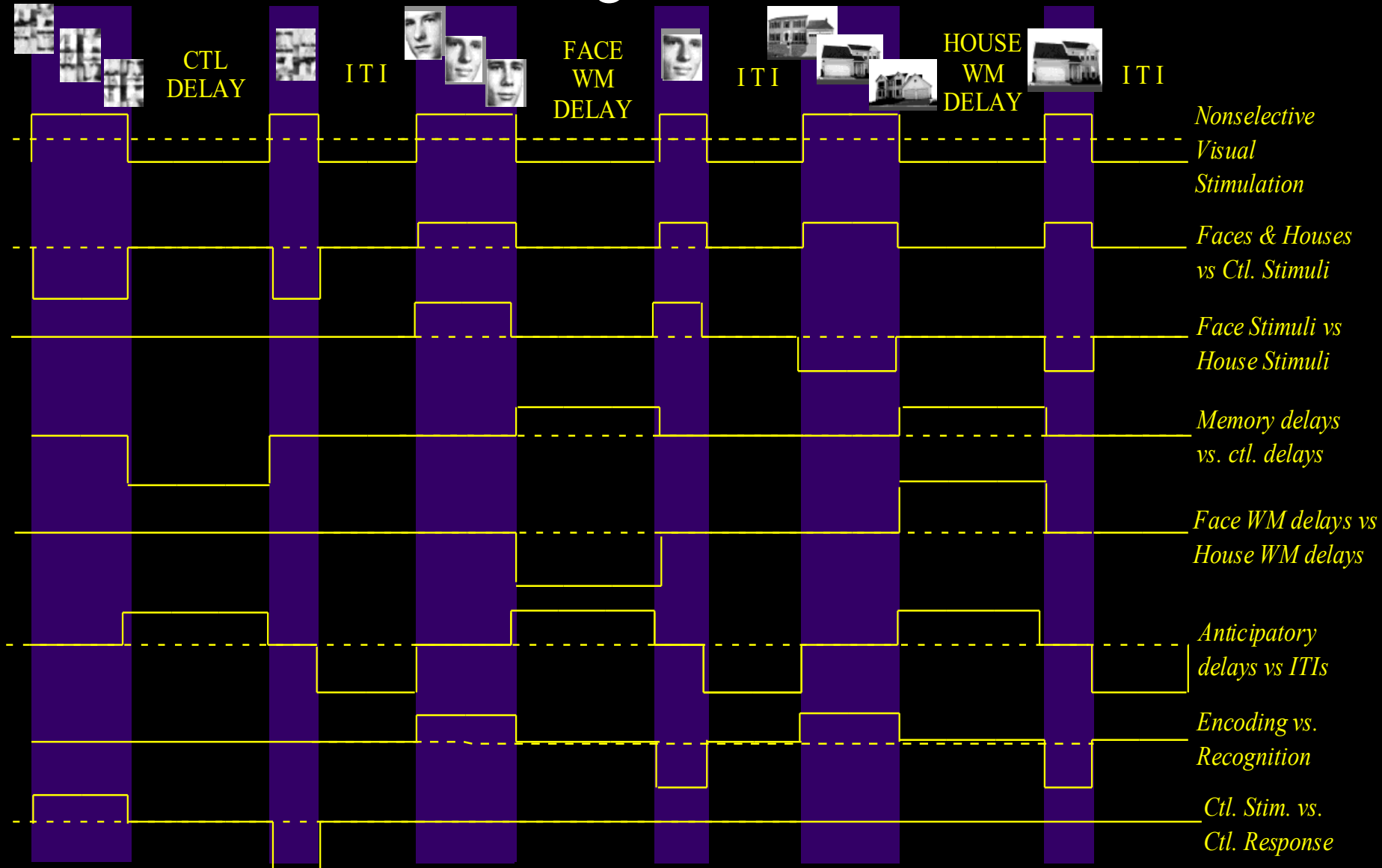


Neuronal Activation Input Strategies

1. Block Design

2. Frequency Encoding

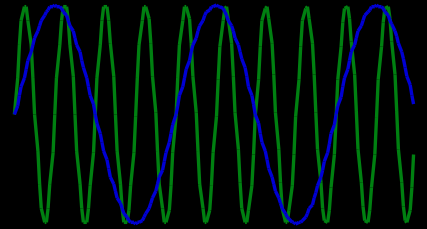


3. Phase Encoding

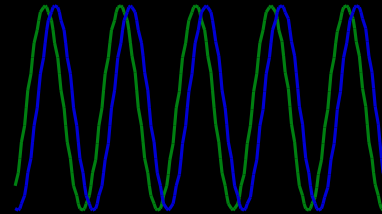

4. Single Event

5. Orthogonal Block Design

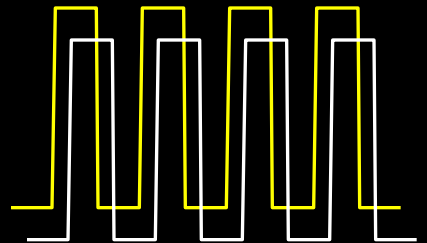



6. Free Behavior Design.

Free Behavior Design


Use a continuous measure as a reference function:

- Task performance
- Skin Conductance
- Heart, respiration rate..
- Eye position
- EEG


Resting Hemodynamic Autocorrelations

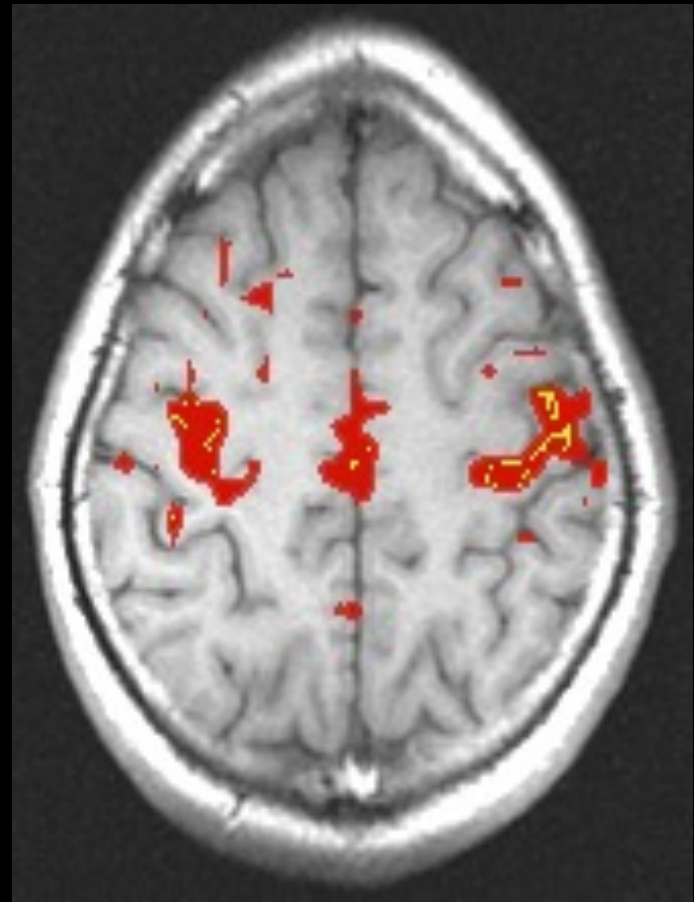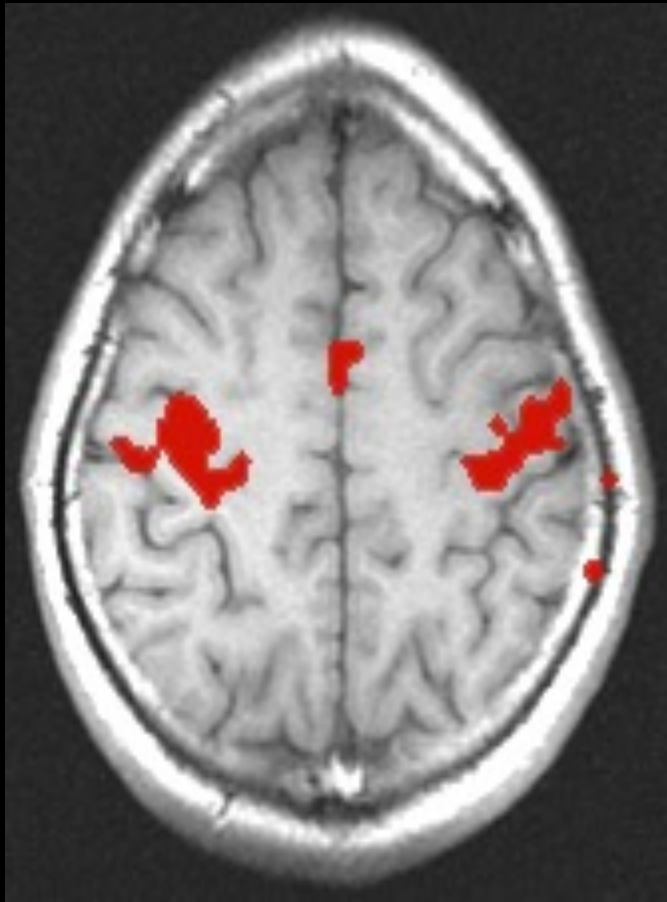


Brain activity correlated with SCR during “Rest”

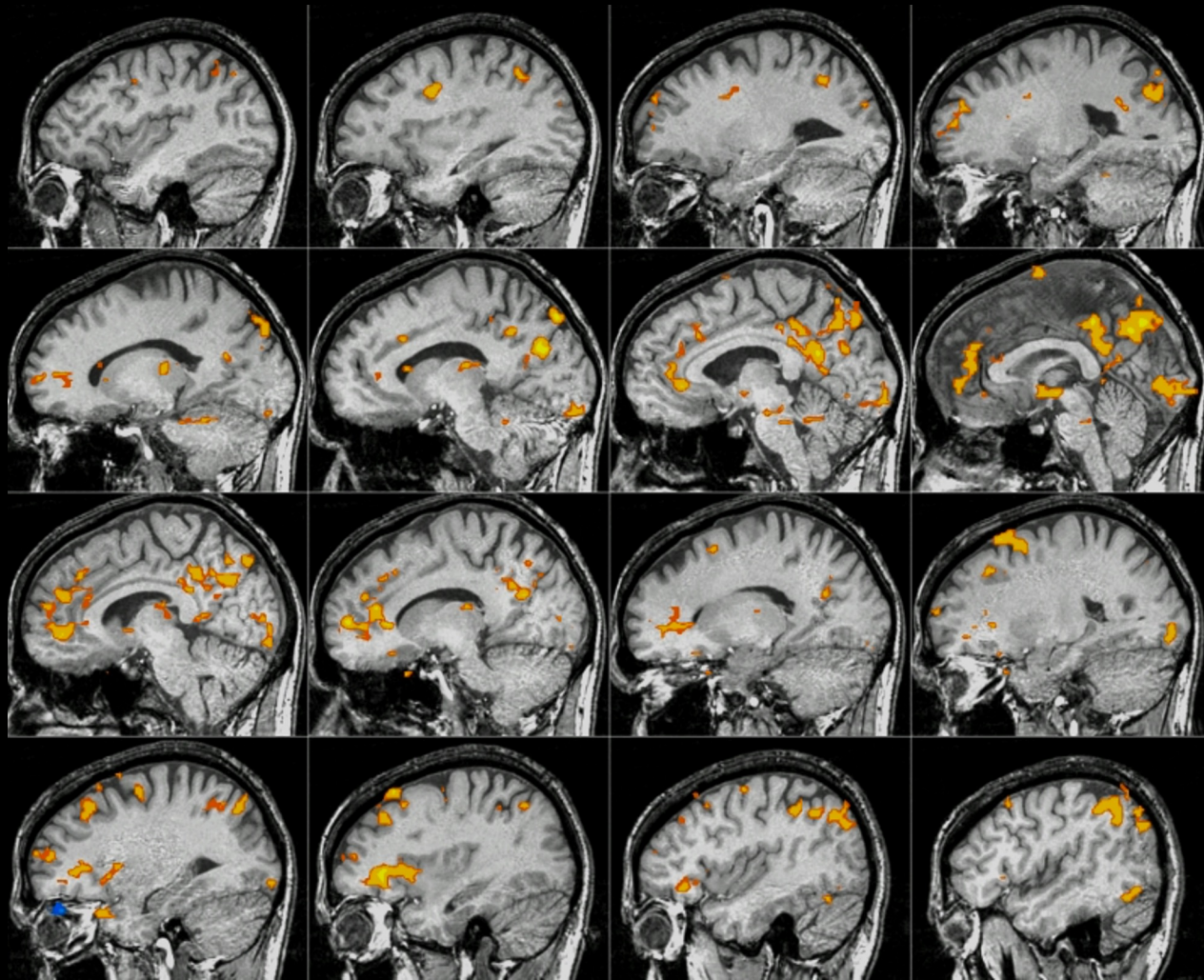

The Skin Conductance Response (SCR)

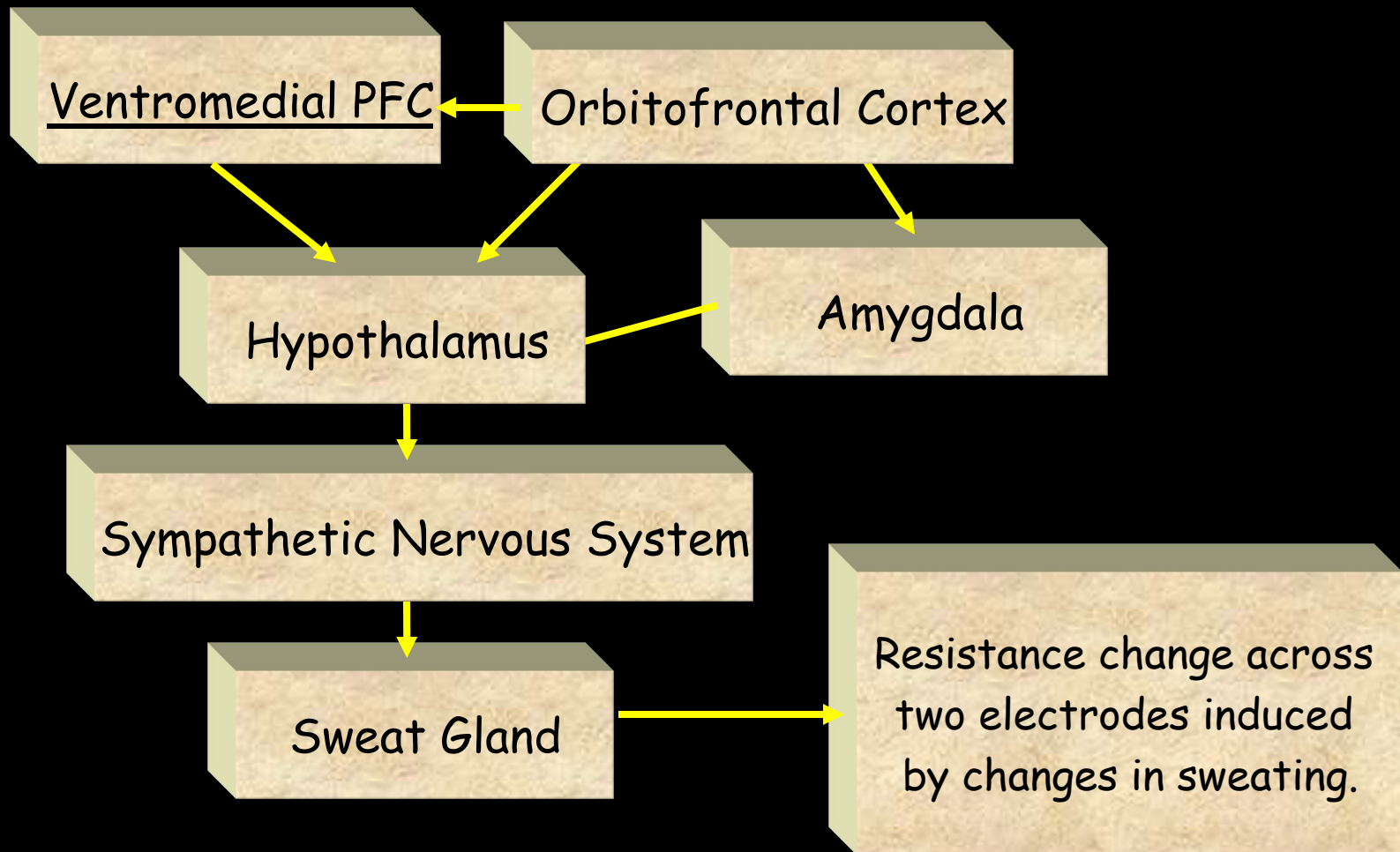

Skin Conductance Dynamics

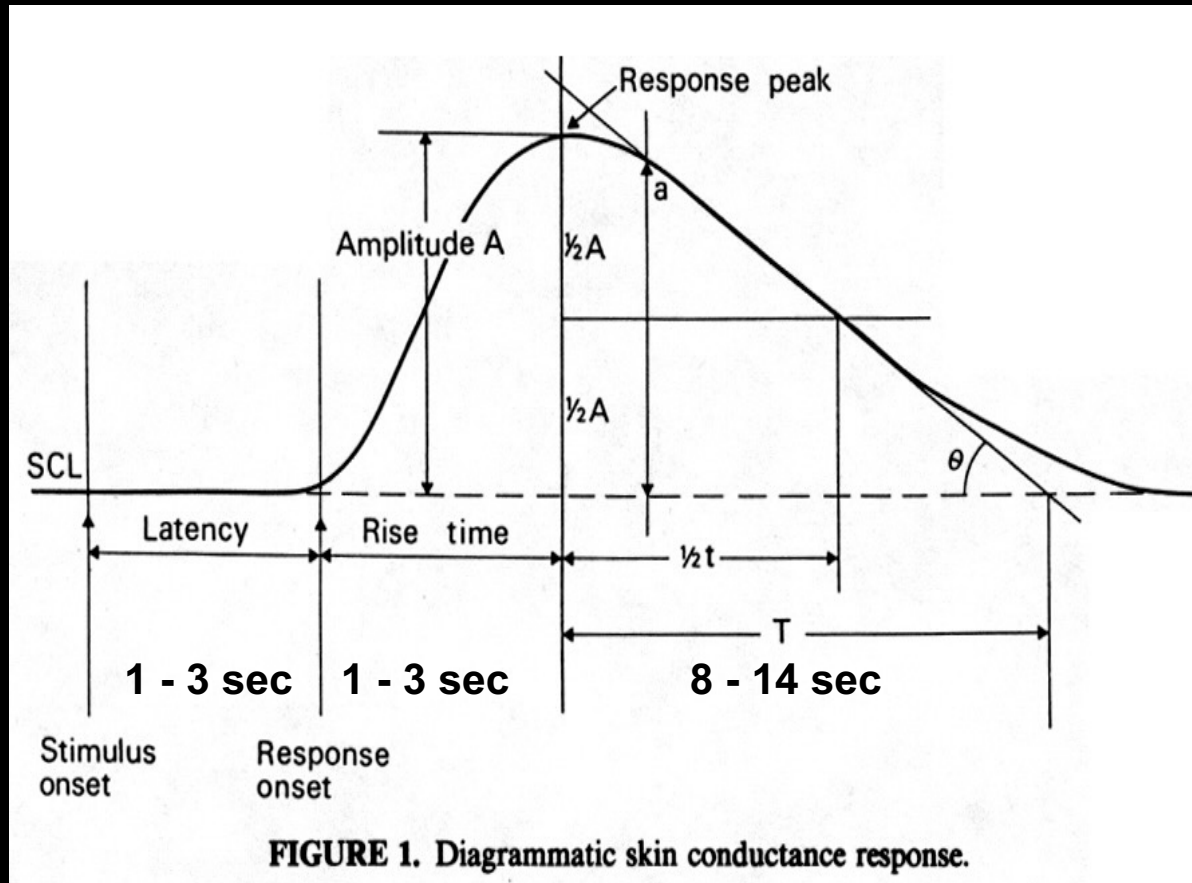

- Boucsein, Wolfram (1992). *Electrodermal Activity*. Plenum Press, NY
- Venables, Peter, (1991). *Autonomic Activity ANYAS 620:191-207.*

Functional Imaging Methods / 3T Group

Staff Scientists:

Sean Marrett

Jerzy Bodurka

Post Docs:

Rasmus Birn

Patrick Bellgowan

Ziad Saad

Graduate Student:


Natalia Petridou

Summer Student:

Dan Kelley

Program Assistant:

Kay Kuhns

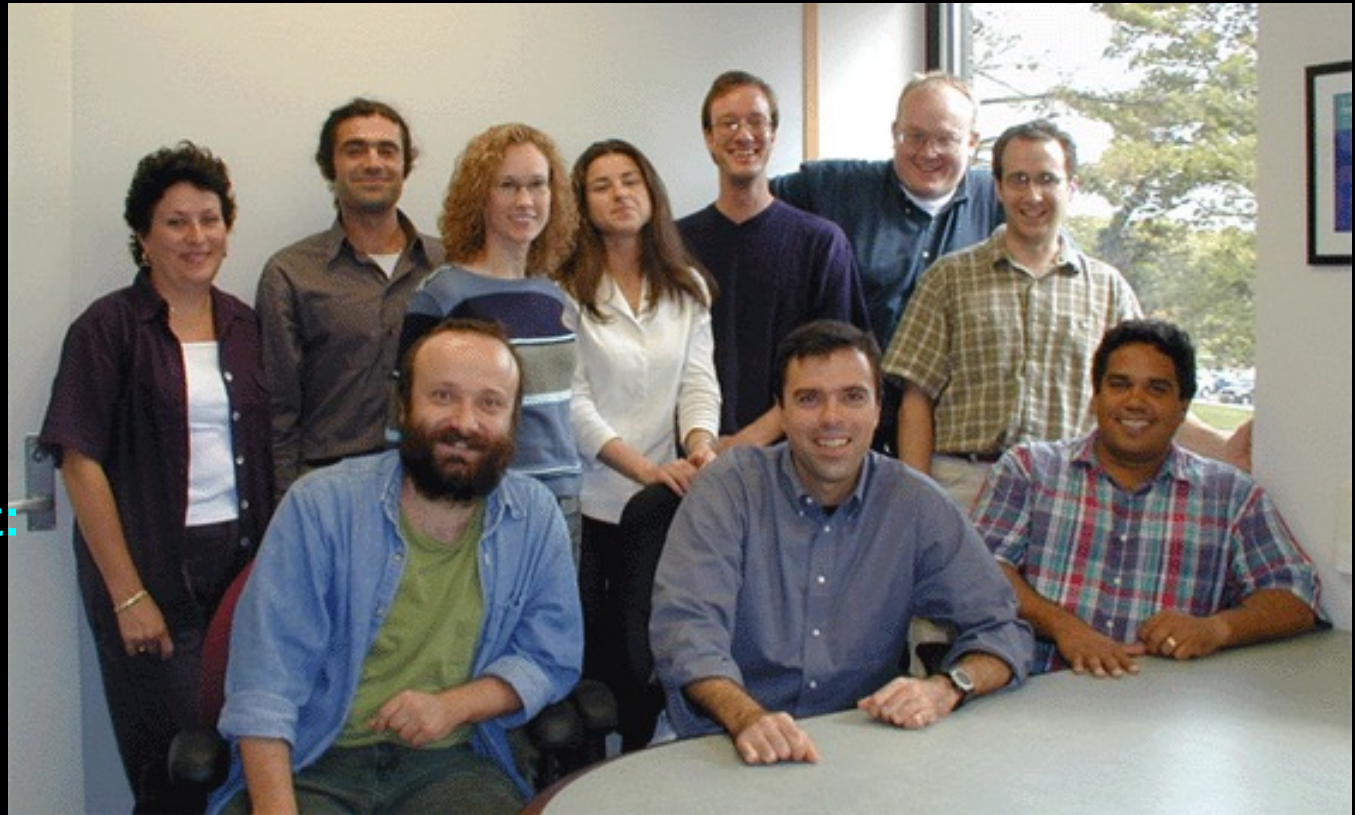

August, 2000

Additional Thanks To...

Eric Wong, UCSD

Robert Savoy, MGH

Peter Jezzard, Oxford

Robert Cox, NIH

Richard Hoge, MGH

Randy Buckner, Wash. U.

Ted DeYoe, MCW

Sue Courtney, Johns Hopkins

Mark Cohen, UCLA